



Veracode Detailed Report  
**Application Security Report**  
**As of 12 Aug 2026**

|                                |                            |
|--------------------------------|----------------------------|
| Prepared for:                  | Danaher Corporation        |
| Prepared on:                   | September 2, 2026          |
| Application:                   | ASG Tempo                  |
| Industry:                      | Not Specified              |
| Business Criticality:          | BC4 (High)                 |
| Required Analysis:             | Static                     |
| Type(s) of Analysis Conducted: | Static                     |
| Scope of Static Scan:          | 29 of 159 Modules Analyzed |

#### Inside This Report

|                              |    |
|------------------------------|----|
| Executive Summary            | 1  |
| Summary of Flaws by Severity | 1  |
| Action Items                 | 1  |
| Flaw Types by Category       | 10 |
| Policy Summary               | 11 |
| Findings & Recommendations   | 12 |
| Flaws in Common Modules      | 29 |
| Methodology                  | 42 |

*While every precaution has been taken in the preparation of this document, Veracode, Inc. assumes no responsibility for errors, omissions, or for damages resulting from the use of the information herein. The Veracode platform uses static and/or dynamic analysis techniques to discover potentially exploitable flaws. Due to the nature of software security testing, the lack of discoverable flaws does not mean the software is 100% secure.*

Veracode Detailed Report  
Application Security Report  
As of 12 Aug 2026

Mitigated Veracode Level: VL1  
Original Veracode Level: VL1  
Rated: Aug 12, 2026

Application:ASG Tempo

Target Level:VL4

Business Criticality:High

Adjusted/Published Rating:C\*/D

Scans Included in Report

| Static Scan                                                    | Dynamic Scan           | Manual Penetration Test |
|----------------------------------------------------------------|------------------------|-------------------------|
| Version:3.2 Build Number:22<br>Score: 78<br>Completed: 8/12/26 | Not Included in Report | Not Included in Report  |

Executive Summary

This report contains a summary of the security flaws identified in the application using manual penetration testing, automated static and/or automated dynamic security analysis techniques. This is useful for understanding the overall security quality of an individual application or for comparisons between applications.

Application Business Criticality: BC4 (High)


Impacts:Operational Risk (Medium), Financial Loss (Medium)

An application's business criticality is determined by business risk factors such as: reputation damage, financial loss, operational risk, sensitive information disclosure, personal safety, and legal violations. The Veracode Level and required assessment techniques are selected based on the policy assigned to the application.

Analyses Performed vs. Required

|            | Any | Static | Dynamic | Manual Penetration Test |
|------------|-----|--------|---------|-------------------------|
| Performed: |     | ●      | ○       | ○                       |
| Required:  | ○   | ●      | ○       | ○                       |

Summary of Flaws Found by Severity

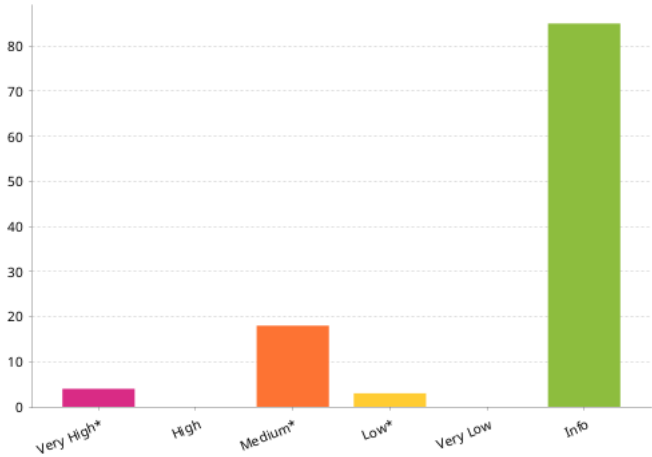

Action Items:

Veracode recommends the following approaches ranging from the most basic to the strong security measures that a vendor can undertake to increase the overall security level of the application.

Required Analysis

- ➔ Your policy requires periodic Static Scan. Your next analysis must be completed by 11/12/26. Please submit your application for Static Scan by the deadline and remediate the required detected flaws to conform to your assigned policy.

Flaws To Fix By Expires Date

A grace period is specified for any flaw that violates the rules contained in your policy. These include CWE, Rollup Category, Issue Severity, Industry Standards as well as any flaws that prevent an application from achieving a minimum Veracode Level and/or score. To maintain policy compliance you must fix these flaws and resubmit your application for scanning before the grace period expires. The detailed flaw listing will badge the flaws that must be fixed and show the fix by date as well.

- The grace period has expired [2/2/26] for 11 flaws that were found in your Static Scan.
- The grace period has expired [2/5/26] for 7 flaws that were found in your Static Scan.
- The grace period has expired [2/22/26] for 1 flaw that was found in your Static Scan.
- The grace period has expired [8/12/26] for 3 flaws that were found in your Static Scan.

## Flaws To Fix For Minimum Score

- A rule in your policy requires a minimum score of 80. Fix 1 Very High flaw to reach a score of 81.
- Your Static Scan was due on 8/12/26 for follow-up analysis to satisfy the grace period on your minimum score rule and your application is no longer compliant with your policy. Submit application for follow-up Static Scan once flaws have been remediated in order to regain compliance with your policy.


## Flaw Severities

- Medium severity flaws and above must be fixed for policy compliance.

## Longer Timeframe (6 - 12 months)

- Certify that software engineers have been trained on application security principles and practices.

Application Trend Data

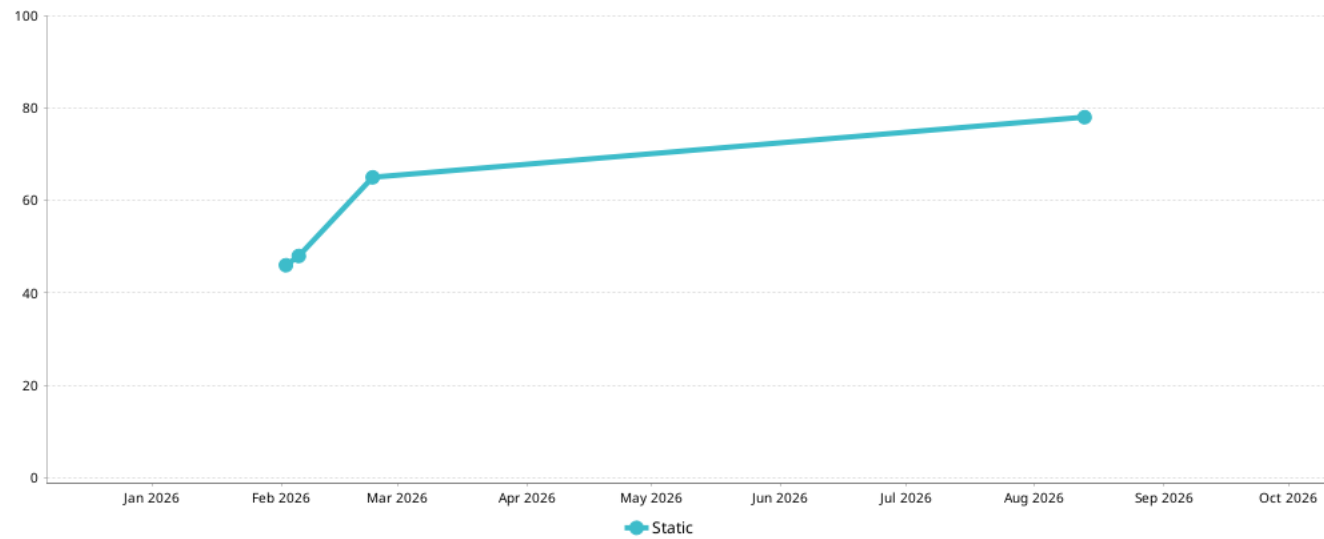

Scope of Static Scan

The following modules were included in the static scan because the scan submitter selected them as entry points, which are modules that accept external data.

Engine Version: 20260806120029

The following modules were included in the application scan:

| Module Name              | Compiler        | Operating Environment | Engine Version |
|--------------------------|-----------------|-----------------------|----------------|
| AmbiStoreClient.exe      | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| ASM.Client.exe           | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| BioShakeClient.exe       | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| CertusClient.exe         | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| CombiNLClient.exe        | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| CombiStdClient.exe       | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| CytomatClient.exe        | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| EchoMassSpec.Client.exe  | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| EnvisionReaderClient.exe | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| InhecoClient.exe         | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| IntellicapClient.exe     | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| KoyoClickPLC.exe         | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| KoyoPLC.exe              | MSIL_MSVC14_X86 | Windows               | 202608061      |

| Module Name              | Compiler        | Operating Environment | Engine Version     |
|--------------------------|-----------------|-----------------------|--------------------|
|                          |                 |                       | 20029              |
| Labcyte.POD.Service.exe  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| MicroServeClient.exe     | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| PlateLocClient.exe       | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| PlateStackerApp.exe      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| PressSealClient.exe      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| RobotTeachingUtility.exe | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| TeleshakeClient.exe      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| TempestClient.exe        | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Tempo.exe                | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| TempoAdmin.exe           | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| TempoScheduleViewer.exe  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| TempoTester.exe          | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| VCodeClient.exe          | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| VolumeCheck.Client.exe   | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| VSpinClient.exe          | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| XPeelClient.exe          | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |



**File Differences Between Scans**

The uploaded modules for this scan do not match the modules you uploaded for the previous scan. This disparity can affect the scan results even if Veracode did not find flaws in the files with differences. See appendix for more details.

The following modules were not selected for a full scan. Code paths in these modules that are not called from a scanned module are not included in this report.

| Module Name                        | Compiler        | Operating Environment | Engine Version     |
|------------------------------------|-----------------|-----------------------|--------------------|
| AGCSN.dll                          | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| ASM.Interface.dll                  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| AxInterop.CentrifugeLib.dll        | MSIL_MSVC11_X86 | Windows               | 202608061<br>20029 |
| AxInterop.MicroplateLabelerLib.dll | MSIL_MSVC11_X86 | Windows               | 202608061<br>20029 |

| Module Name                                 | Compiler        | Operating Environment | Engine Version     |
|---------------------------------------------|-----------------|-----------------------|--------------------|
| AxInterop.PlateLocLib.dll                   | MSIL_MSVC11_X86 | Windows               | 202608061<br>20029 |
| Beckman.Clients.StoreX.dll                  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Beckman.Clients.Wpf.dll                     | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BeckmanCoulter.CoreDelegateCommands.dll     | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BeckmanCoulter.Devices.SampleStore.dll      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BeckmanCoulter.Devices.VolumeCheck.dll      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BeckmanCoulter.Logging.dll                  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BeckmanCoulter.Silas.dll                    | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Device.Plc.dll                    | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Device.Silas.dll                  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Device.Utilities.dll              | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Devices.EchoMS.dll                | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Devices.MobileRobotDoor.dll       | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Devices.SystemStatusIndicator.dll | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Devices.VolumeCheck.dll           | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Driver.BannerTowerLight.dll       | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Driver.BarcodeReader.Silas.dll    | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Driver.Metronome.dll              | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BECLS.POD.Driver.SilasModule.dll            | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| BTILHCRRunner.dll                           | MSIL_MSVC11_X86 | Windows               | 202608061<br>20029 |
| DevAge.Core.dll                             | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| DevAge.Windows.Forms.dll                    | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| EnvisionReaderInterface.dll                 | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| InhcoMTCdll.dll                             | MSIL_MSVC8_X86  | Windows               | 202608061<br>20029 |
| Interop.CentrifugeLib.dll                   | MSIL_MSVC11_X86 | Windows               | 202608061          |

| Module Name                                | Compiler        | Operating Environment | Engine Version     |
|--------------------------------------------|-----------------|-----------------------|--------------------|
|                                            |                 |                       | 20029              |
| Interop.MicroplateLabelerLib.dll           | MSIL_MSVC11_X86 | Windows               | 202608061<br>20029 |
| Interop.PlateLocLib.dll                    | MSIL_MSVC11_X86 | Windows               | 202608061<br>20029 |
| Labcyte.Common.Mailer.dll                  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.Echo.API.dll                       | MSIL_MSVC11_X86 | Windows               | 202608061<br>20029 |
| Labcyte.Echo.Common.dll                    | MSIL_MSVC11_X86 | Windows               | 202608061<br>20029 |
| Labcyte.Echo.Labware.dll                   | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.Echo.Net.dll                       | MSIL_MSVC10_X86 | Windows               | 202608061<br>20029 |
| Labcyte.Echo.Net.Internal.dll              | MSIL_MSVC10_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.API.dll                        | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Base.dll                       | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Common.dll                     | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Database.dll                   | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.AmbiStore.dll          | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.ASM.dll                | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.BarcodeReader.dll      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.BioShake.dll           | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.BioTekWasher.dll       | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.BMG.dll                | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.BulkFillDeviceBase.dll | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.Centrifuge.dll         | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.CertusFlex.dll         | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.CombiNL.dll            | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.CombiStd.dll           | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.Cytomat.dll            | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |

| Module Name                            | Compiler        | Operating Environment | Engine Version |
|----------------------------------------|-----------------|-----------------------|----------------|
| Labcyte.POD.Devices.DeckStorage.dll    | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.DeckTransfer.dll   | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Delidder.dll       | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.DropChute.dll      | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.DynamicStorage.dll | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Echo.dll           | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.EnvisionReader.dll | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.GX.dll             | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.InhecoCPAC.dll     | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.InhecoODTC.dll     | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Insert.dll         | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.IntelliXcap.dll    | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Kaps500.dll        | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Latitude.dll       | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.MicroServe.dll     | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.PlateLoc.dll       | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.PLC.dll            | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.PressSeal.dll      | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Regrip.dll         | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Robot.dll          | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.SampleStore.dll    | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.StoreX.dll         | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Teleshake.dll      | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.Tempest.dll        | MSIL_MSVC14_X86 | Windows               | 20260806120029 |
| Labcyte.POD.Devices.TempStorage.dll    | MSIL_MSVC14_X86 | Windows               | 202608061      |

| Module Name                            | Compiler             | Operating Environment | Engine Version     |
|----------------------------------------|----------------------|-----------------------|--------------------|
|                                        |                      |                       | 20029              |
| Labcyte.POD.Devices.ThermalStation.dll | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.VCode.dll          | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.Virtual.dll        | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Devices.XPeel.dll          | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.AmbiStore.dll       | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.BigBear.dll         | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.BioShake.dll        | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.BioTekWasher.dll    | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.CapDetector.dll     | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.CertusFlex.dll      | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.CombiNL.dll         | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.CombiStd.dll        | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.CommonRobot.dll     | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Cytomat.dll         | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Dmx.dll             | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Envision.dll        | MSIL_CPP_MSVC143_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.GX.dll              | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.HiG.dll             | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.InhecoODTC.dll      | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.InhecoTEC.dll       | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.IntelliXcap.dll     | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Kaps500.dll         | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.KeyenceBarcode.dll  | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Koyo.dll            | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |

| Module Name                         | Compiler        | Operating Environment | Engine Version     |
|-------------------------------------|-----------------|-----------------------|--------------------|
| Labcyte.POD.Driver.KoyoClickPLC.dll | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.MicroServe.dll   | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Precise.dll      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.PressSealer.dll  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Rorze.dll        | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.StoreX.dll       | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Teleshake.dll    | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.Tempest.dll      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.VSpin.dll        | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Driver.XPeel.dll        | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Drivers.Interfaces.dll  | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Interface.BMG.dll       | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.PlateStorageControl.dll | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.Remote.dll              | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.ScheduleViewer.dll      | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.StorageProtocol.dll     | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Labcyte.POD.UserInterface.dll       | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Latitude.Driver.dll                 | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| Metronome.Api.dll                   | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| PlateStackerService.dll             | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| RobotTeachClient.dll                | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| RobotTeachingInterface.dll          | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| SciexDeviceClient.dll               | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| SciexDeviceProtos.dll               | MSIL_MSVC14_X86 | Windows               | 202608061<br>20029 |
| SourceGrid.dll                      | MSIL_MSVC14_X86 | Windows               | 202608061          |

| Module Name              | Compiler             | Operating Environment | Engine Version     |
|--------------------------|----------------------|-----------------------|--------------------|
|                          |                      |                       | 20029              |
| TempoAmrDoors.Models.dll | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |
| TI50UsbDotNetWin32.dll   | MSIL_CPP_MSVC142_X86 | Windows               | 202608061<br>20029 |
| WpfUtils.dll             | MSIL_MSVC14_X86      | Windows               | 202608061<br>20029 |

Flaw Types by Severity and Category

| Static Scan<br>Security Quality Score = 78<br>(+13) from prior scan |      |       |  |
|---------------------------------------------------------------------|------|-------|--|
| Very High                                                           | 4*   |       |  |
| Command or Argument Injection                                       | 4*   |       |  |
| High                                                                | 0    | (-2)  |  |
| SQL Injection                                                       |      | (-2)  |  |
| Medium                                                              | 18*  | (-84) |  |
| CRLF Injection                                                      |      | (-16) |  |
| Credentials Management                                              | 3    | (-1)  |  |
| Cryptographic Issues                                                | 3*   |       |  |
| Directory Traversal                                                 | 5*   | (-69) |  |
| Information Leakage                                                 | 5*   | (+2)  |  |
| Untrusted Search Path                                               | 2*   |       |  |
| Low                                                                 | 3*   | (-1)  |  |
| Information Leakage                                                 | 3*   | (-1)  |  |
| Very Low                                                            | 0    |       |  |
| Informational                                                       | 85   | (+8)  |  |
| Code Quality                                                        | 84   | (+7)  |  |
| Potential Backdoor                                                  | 1    | (+1)  |  |
| Total                                                               | 110* | (-79) |  |

### Policy Evaluation

Policy Name: Veracode Recommended High

Revision: 1

Policy Status: Did Not Pass

Description: Veracode provides default policies to make it easier for organizations to begin measuring their applications against policies. Veracode Recommended Policies are available for customers as an option when they are ready to move beyond the initial bar set by the Veracode Transitional Policies. The policies are based on the Veracode Level definitions.

Rules

| Rule type                | Requirement | Findings         | Status       |
|--------------------------|-------------|------------------|--------------|
| Minimum Veracode Level   | VL4         | VL1*             | Did not pass |
| (VL4) Min Analysis Score | 80          | 78*              | Did not pass |
| (VL4) Max Severity       | Medium      | Flaws found: 22* | Did not pass |
| Max Severity             | Medium      | Flaws found: 22* | Did not pass |
| Min Analysis Score       | 80          | 78*              | Did not pass |

\* - Reflects violated rules that have mitigated flaws

Scan Requirements

| Scan Type | Frequency | Last performed | Status |
|-----------|-----------|----------------|--------|
| Static    | Quarterly | 8/12/26        | Passed |

Remediation

| Flaw Severity | Grace Period | Flaws Exceeding | Status       |
|---------------|--------------|-----------------|--------------|
| Very High     | 0 days       | 4               | Did not pass |
| High          | 0 days       | 0               | Passed       |
| Medium        | 0 days       | 18              | Did not pass |
| Low           | 0 days       | 0               | Passed       |
| Very Low      | 0 days       | 0               | Passed       |
| Informational | 0 days       | 0               | Passed       |

| Type               | Grace Period | Exceeding | Status       |
|--------------------|--------------|-----------|--------------|
| Min Analysis Score | 0 days       | 1         | Did not pass |

## Findings & Recommendations

### Detailed Flaws by Severity

#### Very High (4 flaws\*)

 **Fix Required by Policy**

#### → Command or Argument Injection(4 flaws\*)

##### Description

Command or argument injection vulnerabilities occur when data enters an application from an untrusted source and is used to dynamically construct and execute a command. In the case of OS command injection, an attacker may be able to either alter the command executed by the application or append additional commands. In the case of argument injection, the attacker may influence the behavior of the program in other ways, for example, changing the destination of an outbound network request or injecting additional commands into an argument or parameter. The command is typically executed with the privileges of the executing process and gives an attacker a privilege or capability that he would not otherwise have.

##### Recommendations

Careful handling of all untrusted data is critical in preventing injection attacks. Using one or more of the following techniques provides defense-in-depth and minimizes the likelihood of an vulnerability.

- \* If possible, use library calls rather than external processes to recreate the desired functionality.
- \* Validate user-supplied input using positive filters (white lists) to ensure that it conforms to the expected format, using centralized data validation routines when possible.
- \* Select safe API routines. Some APIs that execute system commands take an array of strings as input rather than a single string, which protects against some forms of command injection by ensuring that a user-supplied argument cannot be interpreted as part of the command.

#### Associated Flaws by CWE ID:

#### → Improper Neutralization of Special Elements used in an OS Command ('OS Command Injection') (CWE ID 78)(4 flaws\*)

##### Description

This call contains a command injection flaw. The argument to the function is constructed using untrusted input. If an attacker is allowed to specify all or part of the command, it may be possible to execute commands on the server with the privileges of the executing process. The level of exposure depends on the effectiveness of input validation routines, if any.

*Effort to Fix:* 3 - Complex implementation error. Fix is approx. 51-500 lines of code. Up to 5 days to fix.

##### Recommendations

Validate all untrusted input to ensure that it conforms to the expected format, using centralized data validation routines when possible. When using black lists, be sure that the sanitizing routine performs a sufficient number of iterations to remove all instances of disallowed characters. Most APIs that execute system commands also have a "safe" version of the method that takes an array of strings as input rather than a single string, which protects against some forms of command injection.

#### Instances found via Static Scan

| Flaw Id                                                                                | Module # | Class # | Module                                                | Location                               | Fix By  |
|----------------------------------------------------------------------------------------|----------|---------|-------------------------------------------------------|----------------------------------------|---------|
|  809  | 14       | -       | tempoadmin.exe#3.2.0.1                                | .../deviceinfoviewmodel.cs 342         | 2/5/26  |
|  582  | 35       | -       | tempo.exe#3.2.0.1                                     | .../launchprogramcommand.cs 42         | 2/5/26  |
|  1005 | 53       | -       | koyopl.exe#3.2.0.1/labcyte.pod.common.dll             | buildagent/.../device/poddevice.cs 772 | 8/12/26 |
|  823  | 81       | -       | tempotester.exe#3.2.0.1/labcyte.pod.userinterface.dll | .../runreportsform.cs 750              | 2/22/26 |

High (0 flaws)

No flaws of this type were found

Medium (18 flaws\*)

 **Fix Required by Policy**

➔ Credentials Management(3 flaws)

Description

Improper management of credentials, such as usernames and passwords, may compromise system security. In particular, storing passwords in plaintext or hard-coding passwords directly into application code are design issues that cannot be easily remedied. Not only does embedding a password allow all of the project's developers to view the password, it also makes fixing the problem extremely difficult. Once the code is in production, the password cannot be changed without patching the software. If a hard-coded password is compromised in a commercial product, all deployed instances may be vulnerable to attack, putting customers at risk.

One variation on hard-coding plaintext passwords is to hard-code a constant string which is the result of a cryptographic one-way hash. For example, instead of storing the word "secret," the application stores an MD5 hash of the word. This is a common mechanism for obscuring hard-coded passwords from casual viewing but does not significantly reduce risk. However, using cryptographic hashes for data stored outside the application code can be an effective practice.

Recommendations

Avoid storing passwords in easily accessible locations, and never store any type of sensitive data in plaintext. Avoid using hard-coded usernames, passwords, or hash constants whenever possible, particularly in relation to security-critical components. Store passwords out-of-band from the application code. Follow best practices for protecting credentials stored in alternate locations such as configuration or properties files.

Associated Flaws by CWE ID:

➔ Use of Hard-coded Password (CWE ID 259)(3 flaws)

Description

A method uses a hard-coded password that may compromise system security in a way that cannot be easily remedied. The use of a hard-coded password significantly increases the possibility that the account being protected will be compromised. Moreover, the password cannot be changed without patching the software. If a hard-coded password is compromised in a commercial product, all deployed instances may be vulnerable to attack.

Effort to Fix: 4 - Simple design error. Requires redesign and up to 5 days to fix.

### Recommendations

Store passwords out-of-band from the application code. Follow best practices for protecting credentials stored in locations such as configuration or properties files.

#### Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                           | Location                             | Fix By |
|---------|----------|---------|----------------------------------|--------------------------------------|--------|
| 384     | 41       | -       | koyopl.exe#3.2.0.1               | buildagent/.../koyopl/mainform.cs 36 | 2/5/26 |
| 682     | 58       | -       | robotteachingutility.exe#3.2.0.1 | .../preciseteachingform.cs 27        | 2/5/26 |
| 681     | 72       | -       | robotteachingutility.exe#3.2.0.1 | .../rorzeteachingfrm.cs 18           | 2/5/26 |

### → Cryptographic Issues(3 flaws\*)

#### Description

Applications commonly use cryptography to implement authentication mechanisms and to ensure the confidentiality and integrity of sensitive data, both in transit and at rest. The proper and accurate implementation of cryptography is extremely critical to its efficacy. Configuration or coding mistakes as well as incorrect assumptions may negate a large degree of the protection it affords, leaving the crypto implementation vulnerable to attack.

Common cryptographic mistakes include, but are not limited to, selecting weak keys or weak cipher modes, unintentionally exposing sensitive cryptographic data, using predictable entropy sources, and mismanaging or hard-coding keys.

Developers often make the dangerous assumption that they can improve security by designing their own cryptographic algorithm; however, one of the basic tenets of cryptography is that any cipher whose effectiveness is reliant on the secrecy of the algorithm is fundamentally flawed.

#### Recommendations

Select the appropriate type of cryptography for the intended purpose. Avoid proprietary encryption algorithms as they typically rely on "security through obscurity" rather than sound mathematics. Select key sizes appropriate for the data being protected; for high assurance applications, 256-bit symmetric keys and 2048-bit asymmetric keys are sufficient. Follow best practices for key storage, and ensure that plaintext data and key material are not inadvertently exposed.

#### Associated Flaws by CWE ID:

### → Use of a Broken or Risky Cryptographic Algorithm (CWE ID 327)(3 flaws)

#### Description

The use of a broken or risky cryptographic algorithm is an unnecessary risk that may result in the disclosure of sensitive information.

*Effort to Fix:* 1 - Trivial implementation error. Fix is up to 5 lines of code. One hour or less to fix.

#### Instances found via Static Scan

| Flaw Id | Module # | Class # | Module              | Location                               | Fix By |
|---------|----------|---------|---------------------|----------------------------------------|--------|
| 625     | 89       | -       | tempo.exe#3.2.0.1/l | .../protocol/selectechoprotocol.cs 203 | 2/5/26 |

| Flaw Id | Module # | Class # | Module                                                        | Location                             | Fix By |
|---------|----------|---------|---------------------------------------------------------------|--------------------------------------|--------|
|         |          |         | abcyte.pod.api.dll                                            |                                      |        |
| 524     | 95       | -       | tempo.exe#3.2.0.1/<br>abcyte.pod.userinter<br>face.dll        | .../ucechoapplications.cs 870        | 2/5/26 |
| 160     | -        | 8       | tempo.exe#3.2.0.1/<br>abcyte.pod.api.dll/bt<br>ilhcrunner.dll | void a(ref string /*1*/, string) 45% | 2/2/26 |

## Directory Traversal(5 flaws\*)

### Description

Allowing user input to control paths used in filesystem operations may enable an attacker to access or modify otherwise protected system resources that would normally be inaccessible to end users. In some cases, the user-provided input may be passed directly to the filesystem operation, or it may be concatenated to one or more fixed strings to construct a fully-qualified path.

When an application improperly cleanses special character sequences in user-supplied filenames, a path traversal (or directory traversal) vulnerability may occur. For example, an attacker could specify a filename such as "../etc/passwd", which resolves to a file outside of the intended directory that the attacker would not normally be authorized to view.

### Recommendations

Assume all user-supplied input is malicious. Validate all user-supplied input to ensure that it conforms to the expected format, using centralized data validation routines when possible. When using black lists, be sure that the sanitizing routine performs a sufficient number of iterations to remove all instances of disallowed characters and ensure that the end result is not dangerous.

### Associated Flaws by CWE ID:

## External Control of File Name or Path (CWE ID 73)(5 flaws\*)

### Description

This call contains a path manipulation flaw. The argument to the function is a filename constructed using untrusted input. If an attacker is allowed to specify all or part of the filename, it may be possible to gain unauthorized access to files on the server, including those outside the webroot, that would be normally be inaccessible to end users. The level of exposure depends on the effectiveness of input validation routines, if any.

*Effort to Fix:* 2 - Implementation error. Fix is approx. 6-50 lines of code. 1 day to fix.

### Recommendations

Validate all untrusted input to ensure that it conforms to the expected format, using centralized data validation routines when possible. When using black lists, be sure that the sanitizing routine performs a sufficient number of iterations to remove all instances of disallowed characters.

### Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                      | Location           | Fix By |
|---------|----------|---------|---------------------------------------------|--------------------|--------|
| 145     | -        | 5       | tempo.exe#3.2.0.1/<br>abcyte.pod.api.dll/bt | void a(object) 99% | 2/2/26 |

| Flaw Id | Module # | Class # | Module                                                 | Location                                           | Fix By |
|---------|----------|---------|--------------------------------------------------------|----------------------------------------------------|--------|
|         |          |         | ilhcrunner.dll                                         |                                                    |        |
| 141     | -        | 5       | tempo.exe#3.2.0.1/abcyte.pod.api.dll/bt ilhcrunner.dll | void a(string, object, object) 98%                 | 2/2/26 |
| 143     | -        | 5       | tempo.exe#3.2.0.1/abcyte.pod.api.dll/bt ilhcrunner.dll | void a(string, object, object, object, object) 99% | 2/2/26 |
| 140     | -        | 5       | tempo.exe#3.2.0.1/abcyte.pod.api.dll/bt ilhcrunner.dll | void ao() 96%                                      | 2/2/26 |
| 142     | -        | 5       | tempo.exe#3.2.0.1/abcyte.pod.api.dll/bt ilhcrunner.dll | void q(string, object) 98%                         | 2/2/26 |

## Information Leakage(5 flaws\*)

### Description

An information leak is the intentional or unintentional disclosure of information that is either regarded as sensitive within the product's own functionality or provides information about the product or its environment that could be useful in an attack. Information leakage issues are commonly overlooked because they cannot be used to directly exploit the application. However, information leaks should be viewed as building blocks that an attacker uses to carry out other, more complicated attacks.

There are many different types of problems that involve information leaks, with severities that can range widely depending on the type of information leaked and the context of the information with respect to the application. Common sources of information leakage include, but are not limited to:

- \* Source code disclosure
- \* Browsable directories
- \* Log files or backup files in web-accessible directories
- \* Unfiltered backend error messages
- \* Exception stack traces
- \* Server version information
- \* Transmission of uninitialized memory containing sensitive data

### Recommendations

Configure applications and servers to return generic error messages and to suppress stack traces from being displayed to end users. Ensure that errors generated by the application do not provide insight into specific backend issues.

Remove all backup files, binary archives, alternate versions of files, and test files from web-accessible directories of production servers. The only files that should be present in the application's web document root are files required by the application. Ensure that deployment procedures include the removal of these file types by an administrator. Keep web and application servers fully patched to minimize exposure to publicly-disclosed information leakage vulnerabilities.

### Associated Flaws by CWE ID:

## → Improper Restriction of XML External Entity Reference (CWE ID 611)(2 flaws\*)

### Description

The product processes an XML document that can contain XML entities with URLs that resolve to documents outside of the intended sphere of control, causing the product to embed incorrect documents into its output. By default, the XML entity resolver will attempt to resolve and retrieve external references. If attacker-controlled XML can be submitted to one of these functions, then the attacker could gain access to information about an internal network, local filesystem, or other sensitive data. This is known as an XML eXternal Entity (XXE) attack.

*Effort to Fix:* 2 - Implementation error. Fix is approx. 6-50 lines of code. 1 day to fix.

### Recommendations

Configure the XML parser to disable external entity resolution.

#### Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                               | Location                       | Fix By  |
|---------|----------|---------|--------------------------------------------------------------------------------------|--------------------------------|---------|
| 832     | 7        | -       | tempo.exe#3.2.0.1/l<br>abcyte.pod.api.dll/la<br>bcyte.pod.devices.c<br>ertusflex.dll | .../certusflexdevice.cs 520    | 8/12/26 |
| 828     | 14       | -       | tempoadmin.exe#3.<br>2.0.1                                                           | .../deviceinfoviewmodel.cs 432 | 8/12/26 |

## → Server-Side Request Forgery (SSRF) (CWE ID 918)(3 flaws)

### Description

The web server receives a URL or similar request from an upstream component and retrieves the contents of this URL, but it does not sufficiently ensure that the request is being sent to the expected destination.

*Effort to Fix:* 1 - Trivial implementation error. Fix is up to 5 lines of code. One hour or less to fix.

#### Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                                                                               | Location            | Fix By |
|---------|----------|---------|--------------------------------------------------------------------------------------------------------------------------------------|---------------------|--------|
| 211     | -        | 3       | tempo.exe#3.2.0.1/l<br>abcyte.pod.api.dll/la<br>bcyte.pod.devices.sa<br>mplestore.dll/beckm<br>ancoulter.devices.sa<br>mplestore.dll | void MoveNext() 23% | 2/2/26 |
| 210     | -        | 2       | tempo.exe#3.2.0.1/l<br>abcyte.pod.api.dll/la<br>bcyte.pod.devices.sa<br>mplestore.dll/beckm<br>ancoulter.devices.sa<br>mplestore.dll | void MoveNext() 27% | 2/2/26 |
| 209     | -        | 1       | tempo.exe#3.2.0.1/l<br>abcyte.pod.api.dll/la<br>bcyte.pod.devices.sa<br>mplestore.dll/beckm<br>ancoulter.devices.sa                  | void MoveNext() 60% | 2/2/26 |

| Flaw Id       | Module # | Class # | Module | Location | Fix By |
|---------------|----------|---------|--------|----------|--------|
| mplestore.dll |          |         |        |          |        |

→ **Untrusted Search Path(2 flaws\*)**

### Description

Executing commands or loading libraries from an untrusted source or in an untrusted environment can cause an application to execute malicious commands (and payloads) on behalf of an attacker.

If an attacker is allowed to specify all or part of a filename to certain API function, it may be possible to load arbitrary libraries. In addition, certain functions perform automatic path searching, iterating over a list of paths to search for desired resources, such as executables, libraries or configuration files. If an attacker can modify the path, for example, by manipulating an environment variable, he may be able to trick the program into referencing an attacker-controlled resource. Similarly, the search path is static but an attacker is able to place a malicious copy of the resource higher in the search order than the file the application intends to load, then the application will load the malicious version.

### Recommendations

Avoid using user-supplied filenames or paths. When calling methods that load libraries or launch processes, ensure that full paths are provided specifying the resource to be loaded.

### Associated Flaws by CWE ID:

→ **Uncontrolled Search Path Element (CWE ID 427)(2 flaws\*)**

### Description

Windows will search the entire path for a program if the full path is not set explicitly. This could allow trojans or other malicious code to be executed instead. Also, be sure that a file extension is specified, since otherwise multiple extensions will be tried by the operating system. For example, .com programs will run before .exe programs, providing an opportunity for program substitution.

*Effort to Fix:* 1 - Trivial implementation error. Fix is up to 5 lines of code. One hour or less to fix.

### Recommendations

Ensure that a full path is provided specifying the program to be executed.

### Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                         | Location                             | Fix By |
|---------|----------|---------|----------------------------------------------------------------|--------------------------------------|--------|
| 139     | -        | 5       | tempo.exe#3.2.0.1/l<br>abcyte.pod.api.dll/bt<br>ilhcrunner.dll | System.Type a(string, string) 29%    | 2/2/26 |
| 158     | -        | 7       | tempo.exe#3.2.0.1/l<br>abcyte.pod.api.dll/bt<br>ilhcrunner.dll | void b(object, System.EventArgs) 94% | 2/2/26 |

Low (3 flaws\*)

➔ Information Leakage(3 flaws\*)

Description

An information leak is the intentional or unintentional disclosure of information that is either regarded as sensitive within the product's own functionality or provides information about the product or its environment that could be useful in an attack. Information leakage issues are commonly overlooked because they cannot be used to directly exploit the application. However, information leaks should be viewed as building blocks that an attacker uses to carry out other, more complicated attacks.

There are many different types of problems that involve information leaks, with severities that can range widely depending on the type of information leaked and the context of the information with respect to the application. Common sources of information leakage include, but are not limited to:

- \* Source code disclosure
- \* Browsable directories
- \* Log files or backup files in web-accessible directories
- \* Unfiltered backend error messages
- \* Exception stack traces
- \* Server version information
- \* Transmission of uninitialized memory containing sensitive data

Recommendations

Configure applications and servers to return generic error messages and to suppress stack traces from being displayed to end users. Ensure that errors generated by the application do not provide insight into specific backend issues.

Remove all backup files, binary archives, alternate versions of files, and test files from web-accessible directories of production servers. The only files that should be present in the application's web document root are files required by the application. Ensure that deployment procedures include the removal of these file types by an administrator. Keep web and application servers fully patched to minimize exposure to publicly-disclosed information leakage vulnerabilities.

Associated Flaws by CWE ID:

➔ Insertion of Sensitive Information Into Sent Data (CWE ID 201)(3 flaws\*)

Description

Sensitive information may be exposed as a result of outbound network connections made by the application.

*Effort to Fix:* 2 - Implementation error. Fix is approx. 6-50 lines of code. 1 day to fix.

Recommendations

Ensure that the transfer of sensitive data is intended and that it does not violate application security policy or user expectations.

Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                                                                 | Location            | Fix By |
|---------|----------|---------|------------------------------------------------------------------------------------------------------------------------|---------------------|--------|
| 499     | -        | 2       | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.samplestore.dll/beckmancoulter.devices.samplestore.dll | void MoveNext() 27% |        |

| Flaw Id | Module # | Class # | Module                                                                                                                                | Location            | Fix By |
|---------|----------|---------|---------------------------------------------------------------------------------------------------------------------------------------|---------------------|--------|
| 308     | -        | 1       | tempoadmin.exe#3.2.0.1/labcyte.pod.ap<br>i.dll/labcyte.pod.dev<br>ices.samplestore.dll/<br>beckmancoulter.dev<br>ices.samplestore.dll | void MoveNext() 60% |        |
| 309     | -        | 1       | tempoadmin.exe#3.2.0.1/labcyte.pod.ap<br>i.dll/labcyte.pod.dev<br>ices.samplestore.dll/<br>beckmancoulter.dev<br>ices.samplestore.dll | void MoveNext() 60% |        |

Very Low (0 flaws)

No flaws of this type were found

Info (85 flaws)

→ Code Quality(84 flaws)

Description

Code quality issues stem from failure to follow good coding practices and can lead to unpredictable behavior. These may include but are not limited to:

- \* Neglecting to remove debug code or dead code
- \* Improper resource management, such as using a pointer after it has been freed
- \* Using the incorrect operator to compare objects
- \* Failing to follow an API or framework specification
- \* Using a language feature or API in an unintended manner

While code quality flaws are generally less severe than other categories and usually are not directly exploitable, they may serve as indicators that developers are not following practices that increase the reliability and security of an application. For an attacker, code quality issues may provide an opportunity to stress the application in unexpected ways.

Recommendations

The wide variance of code quality issues makes it impractical to generalize how these issues should be addressed. Refer to individual categories for specific recommendations.

Associated Flaws by CWE ID:

→ Improper Resource Shutdown or Release (CWE ID 404)(84 flaws)

Description

The application fails to release (or incorrectly releases) a system resource before it is made available for re-use. This condition often occurs with resources such as database connections or file handles. Most unreleased resource issues result in general software reliability problems, but if an attacker can intentionally trigger a resource leak, it may be possible to launch a denial of service attack by depleting the resource pool.

Effort to Fix: 2 - Implementation error. Fix is approx. 6-50 lines of code. 1 day to fix.

## Recommendations

When a resource is created or allocated, the developer is responsible for properly releasing the resource as well as accounting for all potential paths of expiration or invalidation. Ensure that all code paths properly release resources.

### Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                                                             | Location                                                                                                                                          | Fix By |
|---------|----------|---------|--------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------|--------|
| 890     | -        | 22      | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.echo.common.dll | bool a(string, string) 96%                                                                                                                        |        |
| 891     | -        | 22      | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.echo.common.dll | bool b(string, string) 9%                                                                                                                         |        |
| 313     | -        | 6       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/bt<br>ilhcrunner.dll                                                         | bool d(string) 88%                                                                                                                                |        |
| 889     | -        | 12      | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.echo.api.dll    | bool DetectCSVFormat(string, string, int, System.Collections.Generic.List<ImportHeadInfoItem>, ref string /*1*/) 97%                              |        |
| 888     | -        | 12      | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.echo.api.dll    | bool DetectCSVStrFormat(string, string, int, System.Collections.Generic.List<ImportHeadInfoItem>, ref string /*1*/) 97%                           |        |
| 887     | -        | 12      | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.echo.api.dll    | bool OpenCSV(string, int, int, int, int, PickListHeaderInfos, ref System.Collections.Generic.Dictionary<int, string> /*1*/, ref string /*1*/) 44% |        |
| 892     | -        | 22      | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.echo.common.dll | bool SaveExceptionsXML(string) 79%                                                                                                                |        |
| 790     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                                                                          | .../clsmousekeyboardevents.cs 2848                                                                                                                |        |
| 791     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                                                                          | .../clsmousekeyboardevents.cs 2849                                                                                                                |        |
| 792     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                                                                          | .../clsmousekeyboardevents.cs 2857                                                                                                                |        |
| 793     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                                                                          | .../clsmousekeyboardevents.cs 2858                                                                                                                |        |
| 794     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                                                                          | .../clsmousekeyboardevents.cs 2866                                                                                                                |        |

| Flaw Id | Module # | Class # | Module                                                                  | Location                              | Fix By |
|---------|----------|---------|-------------------------------------------------------------------------|---------------------------------------|--------|
|         |          |         | n.dll                                                                   |                                       |        |
| 795     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2867    |        |
| 796     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2875    |        |
| 797     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2876    |        |
| 798     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2884    |        |
| 799     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2885    |        |
| 800     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2893    |        |
| 801     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2894    |        |
| 802     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2902    |        |
| 803     | 8        | -       | temposcheduleviewer.exe#3.2.0.1/agcsn.dll                               | .../clsmousekeyboardevents.cs 2903    |        |
| 808     | 11       | -       | tempoadmin.exe#3.2.0.1                                                  | .../autogen/configfileencryptor.cs 61 |        |
| 677     | 10       | -       | robotteachingutility.exe#3.2.0.1/labbyte.pod.driver.commonrobodt.dll    | .../config/configfileencryptor.cs 65  |        |
| 436     | 20       | -       | cytomatclient.exe#3.2.0.1                                               | buildagent/.../frmdevcaps.cs 90       |        |
| 437     | 20       | -       | cytomatclient.exe#3.2.0.1                                               | buildagent/.../frmdevcaps.cs 150      |        |
| 538     | 25       | -       | tempo.exe#3.2.0.1/labbyte.pod.api.dll/labbyte.pod.devices.gx.dll        | buildagent/.../gxdevice.cs 531        |        |
| 553     | 25       | -       | tempo.exe#3.2.0.1/labbyte.pod.api.dll/labbyte.pod.devices.gx.dll        | buildagent/.../gxdevice.cs 542        |        |
| 537     | 27       | -       | tempo.exe#3.2.0.1/labbyte.pod.api.dll/labbyte.pod.devices.inhecodtc.dll | .../inhecomethodfile.cs 64            |        |
| 554     | 40       | -       | tempo.exe#3.2.0.1/labbyte.common.mailer.dll                             | buildagent/.../mailservice.cs 40      |        |
| 783     | 40       | -       | tempo.exe#3.2.0.1/labbyte.common.mailer.dll                             | buildagent/.../mailservice.cs 82      |        |

| Flaw Id | Module # | Class # | Module                                                                                                                           | Location                        | Fix By |
|---------|----------|---------|----------------------------------------------------------------------------------------------------------------------------------|---------------------------------|--------|
|         |          |         | er.dll                                                                                                                           |                                 |        |
| 385     | 48       | -       | koyopl.exe#3.2.0.1<br>/labcyte.pod.driver.k<br>oyo.dll                                                                           | .../modbusserialtransport.cs 30 |        |
| 736     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 27298      |        |
| 737     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 30720      |        |
| 738     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 31871      |        |
| 739     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 32188      |        |
| 740     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 33874      |        |
| 741     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 34359      |        |
| 742     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 34642      |        |
| 743     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 34957      |        |
| 744     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 35510      |        |
| 745     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl                                                              | .../podb.designer.cs 36168      |        |

| Flaw Id | Module # | Class # | Module                                                                                                                           | Location                   | Fix By |
|---------|----------|---------|----------------------------------------------------------------------------------------------------------------------------------|----------------------------|--------|
|         |          |         | l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll                                                                        |                            |        |
| 746     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 37785 |        |
| 747     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 38175 |        |
| 748     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 38584 |        |
| 749     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 38958 |        |
| 750     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 39242 |        |
| 751     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 39557 |        |
| 752     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 40039 |        |
| 753     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 40492 |        |
| 754     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d<br>atabase.dll | .../podb.designer.cs 40999 |        |
| 755     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.po<br>d.driver.ambistore.dl<br>l/labcyte.pod.comm<br>on.dll/labcyte.pod.d                | .../podb.designer.cs 41568 |        |

| Flaw Id | Module # | Class # | Module                                                                                                          | Location                   | Fix By |
|---------|----------|---------|-----------------------------------------------------------------------------------------------------------------|----------------------------|--------|
|         |          |         | atabase.dll                                                                                                     |                            |        |
| 756     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 42183 |        |
| 757     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 42762 |        |
| 758     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 43176 |        |
| 759     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 43856 |        |
| 760     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 44300 |        |
| 761     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 44811 |        |
| 762     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 45107 |        |
| 763     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 45619 |        |
| 764     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 46013 |        |
| 765     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podb.designer.cs 46340 |        |

| Flaw Id | Module # | Class # | Module                                                                                                             | Location                     | Fix By |
|---------|----------|---------|--------------------------------------------------------------------------------------------------------------------|------------------------------|--------|
| 766     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 46962 |        |
| 767     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 47386 |        |
| 768     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 47825 |        |
| 769     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 48073 |        |
| 770     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 48648 |        |
| 771     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 49226 |        |
| 772     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 50650 |        |
| 773     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 51984 |        |
| 774     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 52395 |        |
| 775     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll<br>labcyte.pod.common.dll/labcyte.pod.atabase.dll | .../podddb.designer.cs 52919 |        |
| 776     | 52       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod                                                                        | .../podddb.designer.cs 53658 |        |

| Flaw Id | Module # | Class # | Module                                                                                                        | Location                                             | Fix By |
|---------|----------|---------|---------------------------------------------------------------------------------------------------------------|------------------------------------------------------|--------|
|         |          |         | d.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.database.dll                                        |                                                      |        |
| 777     | 52       | -       | ambistoreclient.exe #3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.database.dll | .../podb.designer.cs 53985                           |        |
| 778     | 52       | -       | ambistoreclient.exe #3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.database.dll | .../podb.designer.cs 54357                           |        |
| 779     | 52       | -       | ambistoreclient.exe #3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.database.dll | .../podb.designer.cs 54962                           |        |
| 780     | 52       | -       | ambistoreclient.exe #3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.database.dll | .../podb.designer.cs 56232                           |        |
| 781     | 52       | -       | ambistoreclient.exe #3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.database.dll | .../podb.designer.cs 56596                           |        |
| 782     | 52       | -       | ambistoreclient.exe #3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.pod.database.dll | .../podb.designer.cs 57175                           |        |
| 886     | -        | 12      | ambistoreclient.exe #3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.echo.api.dll     | System.IO.MemoryStream c(string) 86%                 |        |
| 893     | -        | 21      | ambistoreclient.exe #3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.echo.api.dll     | System.Xml.XmlDocument a(System.Xml.XmlDocument) 25% |        |
| 535     | 95       | -       | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll                                                               | .../ucechoapplications.cs 868                        |        |
| 314     | -        | 8       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/bt ilhcrunner.dll                                                       | void a(ref string /*1*/, string) 78%                 |        |
| 683     | 97       | -       | robotteachingutility.exe#3.2.0.1                                                                              | .../xcfgfileviewer.cs 165                            |        |

| Flaw Id | Module # | Class # | Module                           | Location                  | Fix By |
|---------|----------|---------|----------------------------------|---------------------------|--------|
| 684     | 97       | -       | robotteachingutility.exe#3.2.0.1 | .../xcfgfileviewer.cs 197 |        |

➔ Potential Backdoor(1 flaw)

Description

Application backdoors are modifications to programs that are designed to bypass security mechanisms or inject malicious functionality. Backdoors are often inserted by rogue developers with legitimate access to the source code or distribution binaries. Backdoors can take many forms, such as hard-coded credentials, "easter egg" style functionality, rootkits, or time bombs, among others.

Recommendations

Investigate all potential backdoors thoroughly to ensure there is no undesirable functionality. If the backdoors are real, eliminate them, and initiate a broader effort to inspect the entire codebase for malicious code. This may require a detailed review of all code, as it is possible to hide a serious attack in only one or two lines of code. These lines may be located almost anywhere in an application and may have been intentionally obfuscated by the attacker.

Associated Flaws by CWE ID:

➔ Reliance on Security Through Obscurity (CWE ID 656)(1 flaw)

Description

The strength of a security mechanism depends heavily on its obscurity, such that knowledge of its algorithms or key data is sufficient to allow the mechanism to be compromised.

*Effort to Fix:* 3 - Complex implementation error. Fix is approx. 51-500 lines of code. Up to 5 days to fix.

Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                                | Location                      | Fix By |
|---------|----------|---------|---------------------------------------------------------------------------------------|-------------------------------|--------|
| 869     | 37       | -       | ambistoreclient.exe<br>#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.base.dll | buildagent/.../lickey.cs 1184 |        |

## Flaws in Common Modules

This section highlights the score impact of flaws in common modules in this application.

Module: Labcyte.POD.API.dll Used by 5 executables; score impact: 11

| Location                       | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|--------------------------------|----------|-------------|---------------------|--------|----------------|
| runexecutorreports.cs 111      | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorreports.cs 222      | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorreports.cs 242      | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorreports.cs 426      | 3        | 2           | Directory Traversal | 73     | Neutral        |
| runexecutorreports.cs 881      | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorreports.cs 886      | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorreports.cs 954      | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorreports.cs 959      | 3        | 3           | Directory Traversal | 73     | Neutral        |
| staggerplatejob.cs 217         | 3        | 3           | Directory Traversal | 73     | Neutral        |
| staggerplatejob.cs 225         | 3        | 3           | Directory Traversal | 73     | Neutral        |
| staggerplatejob.cs 262         | 3        | 3           | Directory Traversal | 73     | Neutral        |
| staggerplatejob.cs 343         | 3        | 3           | Directory Traversal | 73     | Neutral        |
| staggerplatejob.cs 426         | 3        | 3           | Directory Traversal | 73     | Neutral        |
| staggerplatejob.cs 447         | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorruncreation.cs 621  | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorruncreation.cs 626  | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorruncreation.cs 725  | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorruncreation.cs 3404 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runexecutorruncreation.cs 3415 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| savetimepointjob.cs 205        | 3        | 3           | Directory Traversal | 73     | Neutral        |
| savetimepointjob.cs 250        | 3        | 3           | Directory Traversal | 73     | Neutral        |
| savetimepointjob.cs 258        | 3        | 3           | Directory Traversal | 73     | Neutral        |
| savetimepointjob.cs 288        | 3        | 3           | Directory Traversal | 73     | Neutral        |
| savetimepointjob.cs 308        | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runreporthandler.cs 26         | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runreporthandler.cs 94         | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runreporthandler.cs 135        | 3        | 2           | Directory Traversal | 73     | Neutral        |
| delayplatejob.cs 348           | 3        | 3           | Directory Traversal | 73     | Neutral        |
| delayplatejob.cs 428           | 3        | 2           | Directory Traversal | 73     | Neutral        |
| podprotocol.cs 2756            | 3        | 1           | Directory Traversal | 73     | Neutral        |
| podprotocol.cs 2889            | 3        | 1           | Directory Traversal | 73     | Neutral        |
| podsystem.cs 152               | 3        | 2           | Directory Traversal | 73     | Neutral        |
| podsystem.cs 161               | 3        | 2           | Directory Traversal | 73     | Neutral        |
| querytimepointjob.cs 243       | 3        | 3           | Directory Traversal | 73     | Neutral        |
| querytimepointjob.cs 326       | 3        | 3           | Directory Traversal | 73     | Neutral        |

| Location                         | Severity | # Instances | Flaw Category        | CWE ID | Exploitability |
|----------------------------------|----------|-------------|----------------------|--------|----------------|
| removetimepointjob.cs 157        | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| removetimepointjob.cs 222        | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| rundefinition.cs 523             | 3        | 1           | Directory Traversal  | 73     | Neutral        |
| rundefinition.cs 679             | 3        | 1           | Directory Traversal  | 73     | Neutral        |
| runscheduler.log.cs 114          | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| runscheduler.log.cs 122          | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| runscheduler.scheduler.cs 2846   | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| runscheduler.scheduler.cs 2882   | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| runscheduler.validationlog.cs 64 | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| runscheduler.validationlog.cs 78 | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| selectechoprotocol.cs 201        | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| selectechoprotocol.cs 203        | 3        | 3           | Cryptographic Issues | 327    | Likely         |
| runexecutormisc.cs 157           | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| runscheduler.cs 201              | 3        | 1           | Directory Traversal  | 73     | Neutral        |
| runsummarygenerator.cs 129       | 3        | 3           | Directory Traversal  | 73     | Neutral        |

Module: Labcyte.POD.UserInterface.dll Used by 2 executables; score impact: 10

| Location                  | Severity | # Instances | Flaw Category                 | CWE ID | Exploitability |
|---------------------------|----------|-------------|-------------------------------|--------|----------------|
| runreportsform.cs 750     | 5        | 1           | Command or Argument Injection | 78     | Neutral        |
| runreportsform.cs 607     | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| runreportsform.cs 620     | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| runreportsform.cs 657     | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| runreportsform.cs 662     | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| uccapcheckimages.cs 350   | 5        | 2           | Command or Argument Injection | 78     | Neutral        |
| devicetestform.cs 2397    | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| devicetestform.cs 2585    | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| devicetestform.cs 2590    | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| barcodeform.cs 256        | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| barcodeform.cs 261        | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| ucechoapplications.cs 868 | 0        | 2           | Code Quality                  | 404    | Neutral        |
| ucechoapplications.cs 868 | 3        | 2           | Directory Traversal           | 73     | Neutral        |
| ucechoapplications.cs 870 | 3        | 2           | Cryptographic Issues          | 327    | Likely         |
| runstatusform.cs 2381     | 3        | 2           | Directory Traversal           | 73     | Neutral        |

Module: Labcyte.POD.Common.dll Used by 96 executables; score impact: 8

| Location         | Severity | # Instances | Flaw Category                 | CWE ID | Exploitability |
|------------------|----------|-------------|-------------------------------|--------|----------------|
| poddevice.cs 772 | 5        | 5           | Command or Argument Injection | 78     | Neutral        |
| podlogger.cs 123 | 3        | 3           | Directory Traversal           | 73     | Neutral        |

| Location                           | Severity | # Instances | Flaw Category        | CWE ID | Exploitability |
|------------------------------------|----------|-------------|----------------------|--------|----------------|
| podlogger.cs 125                   | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| podlogger.cs 194                   | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| podlogger.cs 198                   | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| podlogger.cs 239                   | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| podlogger.cs 275                   | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| jobpropertymanager.cs 66           | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| jobpropertymanager.cs 70           | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| jobpropertymanager.cs 109          | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| jobpropertymanager.cs 122          | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| jobpropertymanager.cs 150          | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| devicereportmanager.cs 387         | 3        | 1           | Directory Traversal  | 73     | Neutral        |
| devicereportmanager.cs 584         | 3        | 1           | Directory Traversal  | 73     | Neutral        |
| devicereportmanager.cs 586         | 3        | 1           | Directory Traversal  | 73     | Neutral        |
| plateimportexport.cs 61            | 3        | 2           | Directory Traversal  | 73     | Neutral        |
| plateimportexport.cs 325           | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| plateimportexport.cs 329           | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| jobmanagerbase.cs 743              | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| jobmanagerbase.cs 746              | 3        | 3           | Directory Traversal  | 73     | Neutral        |
| filestreamlock.cs 74               | 3        | 4           | Directory Traversal  | 73     | Neutral        |
| jobdevicepoolingcalculations.cs 72 | 3        | 3           | Cryptographic Issues | 331    | Unlikely       |
| jobmanager.cs 1670                 | 3        | 3           | Directory Traversal  | 73     | Neutral        |

Module: Labcyte.Echo.API.dll Used by 99 executables; score impact: 7

| Location                                                 | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|----------------------------------------------------------|----------|-------------|---------------------|--------|----------------|
| labcyte_echo_api_dll.Labcyte.Echo.API.OutputGenerator    |          |             |                     |        |                |
| SaveAppData 80%                                          | 3        | 11          | Directory Traversal | 73     | Neutral        |
| SaveAppData 76%                                          | 3        | 11          | Directory Traversal | 73     | Neutral        |
| SaveAppData 21%                                          | 3        | 3           | Directory Traversal | 73     | Neutral        |
| SaveAppData 40%                                          | 3        | 8           | Directory Traversal | 73     | Neutral        |
| CreateOutputDirectories 39%                              | 3        | 1           | Directory Traversal | 73     | Neutral        |
| SaveAppData 57%                                          | 3        | 12          | Directory Traversal | 73     | Neutral        |
| a 62%                                                    | 3        | 20          | Directory Traversal | 73     | Neutral        |
| CreateOutputDirectories 89%                              | 3        | 1           | Directory Traversal | 73     | Neutral        |
| CreateOutputDirectories 57%                              | 3        | 1           | Directory Traversal | 73     | Neutral        |
| SaveAppData 78%                                          | 3        | 3           | Directory Traversal | 73     | Neutral        |
| a 90%                                                    | 3        | 20          | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.ImportPickListForm |          |             |                     |        |                |
| a 47%                                                    | 3        | 2           | Directory Traversal | 73     | Neutral        |

| Location                                                             | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|----------------------------------------------------------------------|----------|-------------|---------------------|--------|----------------|
| b 81%                                                                | 3        | 21          | Directory Traversal | 73     | Neutral        |
| a 40%                                                                | 3        | 21          | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.BulkFillOperation     |          |             |                     |        |                |
| SaveCSV 56%                                                          | 3        | 3           | Directory Traversal | 73     | Neutral        |
| LoadCSV 91%                                                          | 3        | 3           | Directory Traversal | 73     | Neutral        |
| SaveCSV 18%                                                          | 3        | 3           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.EchoBase              |          |             |                     |        |                |
| Save 42%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| Save 87%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| Load 71%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.EchoProtocol          |          |             |                     |        |                |
| Load 47%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| Save 64%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| Save 87%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.SurveyLog.SurveyData           |          |             |                     |        |                |
| Save 33%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| Save 66%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| b 5%                                                                 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.PexLogger                      |          |             |                     |        |                |
| Add 48%                                                              | 3        | 3           | Directory Traversal | 73     | Neutral        |
| a 95%                                                                | 3        | 3           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.PlateInfo                      |          |             |                     |        |                |
| Save 25%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| Save 49%                                                             | 3        | 1           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.EchoPlateSurvey       |          |             |                     |        |                |
| get_SurveyResults 22%                                                | 3        | 6           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.EchoTransferMap       |          |             |                     |        |                |
| get_TransferResults 26%                                              | 3        | 6           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.SurveyLog.SurveyHistoryManager |          |             |                     |        |                |
| SetSurveyForPlate 65%                                                | 3        | 3           | Directory Traversal | 73     | Neutral        |
| labcyte_echo_api_dll.Labcyte.Echo.API.PickList.EchoPickList          |          |             |                     |        |                |
| c 86%                                                                | 0        | 21          | Code Quality        | 404    | Neutral        |
| DetectCSVStrFormat 97%                                               | 0        | 21          | Code Quality        | 404    | Neutral        |
| OpenCSV 44%                                                          | 0        | 23          | Code Quality        | 404    | Neutral        |
| DetectCSVFormat 97%                                                  | 0        | 23          | Code Quality        | 404    | Neutral        |
| labcyte_echo_api_dll.k                                               |          |             |                     |        |                |
| a 25%                                                                | 0        | 21          | Code Quality        | 404    | Neutral        |

Module: BECLS.POD.Devices.MobileRobotDoor.dll Used by 6 executables; score impact: 3

| Location                     | Severity | # Instances | Flaw Category                 | CWE ID | Exploitability |
|------------------------------|----------|-------------|-------------------------------|--------|----------------|
| mobilerobotdoordevice.cs 104 | 5        | 3           | Command or Argument Injection | 78     | Neutral        |

Module: BTILHCRRunner.dll Used by 8 executables; score impact: 2

| Location           | Severity | # Instances | Flaw Category         | CWE ID | Exploitability |
|--------------------|----------|-------------|-----------------------|--------|----------------|
| btilhcrunner_dll.g |          |             |                       |        |                |
| a 29%              | 3        | 3           | Untrusted Search Path | 427    | Neutral        |
| q 98%              | 3        | 3           | Directory Traversal   | 73     | Neutral        |
| a 98%              | 3        | 3           | Directory Traversal   | 73     | Neutral        |
| a 99%              | 3        | 3           | Directory Traversal   | 73     | Neutral        |
| ao 96%             | 3        | 3           | Directory Traversal   | 73     | Neutral        |
| a 99%              | 3        | 3           | Directory Traversal   | 73     | Neutral        |
| btilhcrunner_dll.m |          |             |                       |        |                |
| b 94%              | 3        | 3           | Untrusted Search Path | 427    | Neutral        |
| btilhcrunner_dll.n |          |             |                       |        |                |
| a 78%              | 0        | 3           | Code Quality          | 404    | Neutral        |
| a 45%              | 3        | 3           | Cryptographic Issues  | 327    | Likely         |
| btilhcrunner_dll.j |          |             |                       |        |                |
| d 88%              | 0        | 3           | Code Quality          | 404    | Neutral        |

Module: Labcyte.POD.Devices.InhecoODTC.dll Used by 6 executables; score impact: 2

| Location                | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-------------------------|----------|-------------|---------------------|--------|----------------|
| inhecomethodfile.cs 31  | 3        | 2           | Directory Traversal | 73     | Likely         |
| inhecomethodfile.cs 64  | 0        | 3           | Code Quality        | 404    | Neutral        |
| inhecomethodfile.cs 64  | 3        | 2           | Directory Traversal | 73     | Neutral        |
| inhecomethodfile.cs 68  | 3        | 2           | Directory Traversal | 73     | Neutral        |
| inhecomethodfile.cs 127 | 3        | 2           | Directory Traversal | 73     | Neutral        |
| inhecomethodfile.cs 162 | 3        | 2           | Directory Traversal | 73     | Likely         |
| inhecomethodfile.cs 245 | 3        | 2           | Directory Traversal | 73     | Neutral        |
| inhecomethodfile.cs 317 | 3        | 2           | Directory Traversal | 73     | Neutral        |
| inhecoodtcdevice.cs 834 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.VCode.dll Used by 6 executables; score impact: 1

| Location                      | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-------------------------------|----------|-------------|---------------------|--------|----------------|
| labelplateusingfilejob.cs 124 | 3        | 1           | Directory Traversal | 73     | Likely         |
| labelplateusingfilejob.cs 212 | 3        | 1           | Directory Traversal | 73     | Likely         |
| labelplateusingfilejob.cs 239 | 3        | 1           | Directory Traversal | 73     | Likely         |

| Location                      | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-------------------------------|----------|-------------|---------------------|--------|----------------|
| labelplateusingfilejob.cs 267 | 3        | 1           | Directory Traversal | 73     | Neutral        |
| labelplateusingfilejob.cs 292 | 3        | 1           | Directory Traversal | 73     | Likely         |
| lookuptableform.cs 221        | 3        | 3           | Directory Traversal | 73     | Neutral        |
| lookuptableform.cs 280        | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.Cytomat.dll Used by 7 executables; score impact: 1

| Location              | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-----------------------|----------|-------------|---------------------|--------|----------------|
| cytomatdevice.cs 1732 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| cytomatdevice.cs 1746 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| cytomatdevice.cs 1782 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| cytomatdevice.cs 1883 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| cytomatdevice.cs 1892 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| cytomatdevice.cs 1905 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: SciexDeviceProtos.dll Used by 9 executables; score impact: 1

| Location                                                                   | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|----------------------------------------------------------------------------|----------|-------------|---------------------|--------|----------------|
| sciexdeviceprotos_dll.BeckmanCoulter.Devices.EchoMS.Protos.BatchFileParser |          |             |                     |        |                |
| Write 32%                                                                  | 3        | 1           | Directory Traversal | 73     | V.Likely       |
| Write 22%                                                                  | 3        | 1           | Directory Traversal | 73     | V.Likely       |
| Write 55%                                                                  | 3        | 1           | Directory Traversal | 73     | V.Likely       |
| Write 59%                                                                  | 3        | 1           | Directory Traversal | 73     | V.Likely       |
| !ctor 40%                                                                  | 3        | 1           | Directory Traversal | 73     | V.Likely       |

Module: Labcyte.POD.Devices.Echo.dll Used by 8 executables; score impact: 1

| Location           | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|--------------------|----------|-------------|---------------------|--------|----------------|
| echodevice.cs 1126 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| echodevice.cs 1188 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| surveyjob.cs 167   | 3        | 2           | Directory Traversal | 73     | Likely         |
| surveyjob.cs 249   | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.IntelliXcap.dll Used by 6 executables; score impact: 1

| Location                  | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|---------------------------|----------|-------------|---------------------|--------|----------------|
| intellixcapdevice.cs 379  | 3        | 3           | Directory Traversal | 73     | Neutral        |
| intellixcapdevice.cs 1313 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| intellixcapdevice.cs 1327 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| intellixcapdevice.cs 1366 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Driver.CommonRobot.dll Used by 12 executables; score impact: 1

| Location                  | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|---------------------------|----------|-------------|---------------------|--------|----------------|
| robotconfig.cs 443        | 3        | 2           | Directory Traversal | 73     | Neutral        |
| robotconfig.cs 594        | 3        | 4           | Directory Traversal | 73     | Neutral        |
| robotconfig.cs 637        | 3        | 4           | Directory Traversal | 73     | Neutral        |
| robotconfig.cs 649        | 3        | 4           | Directory Traversal | 73     | Neutral        |
| configfileencryptor.cs 65 | 0        | 4           | Code Quality        | 404    | Neutral        |

Module: Labcyte.Echo.Labware.dll Used by 3 executables; score impact: 1

| Location                                                       | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|----------------------------------------------------------------|----------|-------------|---------------------|--------|----------------|
| labcyte_echo_labware_dll.Labcyte.Echo.Labware.PlateInfoManager |          |             |                     |        |                |
| Save 50%                                                       | 3        | 1           | Directory Traversal | 73     | Neutral        |
| Save 30%                                                       | 3        | 1           | Directory Traversal | 73     | Neutral        |
| Load 42%                                                       | 3        | 2           | Directory Traversal | 73     | V.Likely       |

Module: Labcyte.POD.Devices.AmbiStore.dll Used by 8 executables; score impact: 1

| Location                | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-------------------------|----------|-------------|---------------------|--------|----------------|
| ambistoredevice.cs 2021 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| ambistoredevice.cs 2035 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| ambistoredevice.cs 2071 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.CertusFlex.dll Used by 6 executables; score impact: 1

| Location                | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-------------------------|----------|-------------|---------------------|--------|----------------|
| certusflexdevice.cs 520 | 3        | 3           | Information Leakage | 611    | Neutral        |
| certusflexdevice.cs 698 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| certusflexdevice.cs 838 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.CombiNL.dll Used by 6 executables; score impact: 1

| Location              | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-----------------------|----------|-------------|---------------------|--------|----------------|
| combinldevice.cs 1106 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| combinldevice.cs 1120 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| combinldevice.cs 1150 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.InhecoCPAC.dll Used by 6 executables; score impact: 1

| Location               | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|------------------------|----------|-------------|---------------------|--------|----------------|
| inhecopacdevice.cs 398 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| inhecopacdevice.cs 412 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| inhecopacdevice.cs 449 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.MicroServe.dll Used by 7 executables; score impact: 1

| Location                 | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|--------------------------|----------|-------------|---------------------|--------|----------------|
| microservedevice.cs 2881 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| microservedevice.cs 2895 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| microservedevice.cs 2931 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.Robot.dll Used by 6 executables; score impact: 1

| Location            | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|---------------------|----------|-------------|---------------------|--------|----------------|
| robotdevice.cs 1652 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| robotdevice.cs 1674 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| robotdevice.cs 1713 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.Tempest.dll Used by 6 executables; score impact: 1

| Location              | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-----------------------|----------|-------------|---------------------|--------|----------------|
| bulkfilljob.cs 603    | 3        | 3           | Directory Traversal | 73     | Neutral        |
| bulkfilljob.cs 658    | 3        | 3           | Directory Traversal | 73     | Neutral        |
| runprotocoljob.cs 361 | 3        | 2           | Directory Traversal | 73     | Neutral        |

Module: SciexDeviceClient.dll Used by 8 executables; score impact: 1

| Location                                                                                           | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|----------------------------------------------------------------------------------------------------|----------|-------------|---------------------|--------|----------------|
| sciexdeviceclient_dll.BeckmanCoulter.Devices.EchoMS.SciexDeviceClientWrapper._3CRunSamples_3Ed__19 |          |             |                     |        |                |
| MoveNext 34%                                                                                       | 3        | 3           | Directory Traversal | 73     | Neutral        |
| MoveNext 37%                                                                                       | 3        | 1           | Directory Traversal | 73     | Neutral        |
| MoveNext 36%                                                                                       | 3        | 1           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.Common.Mailer.dll Used by 2 executables; score impact: 1

| Location           | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|--------------------|----------|-------------|---------------------|--------|----------------|
| mailservice.cs 40  | 0        | 1           | Code Quality        | 404    | Neutral        |
| mailservice.cs 82  | 0        | 2           | Code Quality        | 404    | Neutral        |
| mailservice.cs 236 | 3        | 1           | Directory Traversal | 73     | Neutral        |
| mailservice.cs 505 | 3        | 2           | Directory Traversal | 73     | Neutral        |
| mailservice.cs 552 | 2        | 2           | Information Leakage | 209    | Neutral        |
| mailservice.cs 552 | 2        | 2           | Information Leakage | 201    | Unlikely       |

Module: BeckmanCoulter.Devices.SampleStore.dll Used by 7 executables; score impact: < 1

| Location                                                                                                               | Severity | # Instances | Flaw Category | CWE ID | Exploitability |
|------------------------------------------------------------------------------------------------------------------------|----------|-------------|---------------|--------|----------------|
| beckmancoulter_devices_samplestore_dll.BeckmanCoulter.Devices.SampleStore.OAuth.OAuthTokenGenerator._3CGetToken_3Ed__1 |          |             |               |        |                |

| Location                                                                                                                   | Severity | # Instances | Flaw Category       | CWE ID | Exploitability     |
|----------------------------------------------------------------------------------------------------------------------------|----------|-------------|---------------------|--------|--------------------|
| MoveNext 27%                                                                                                               | 3        | 4           | Information Leakage | 918    | Unlikely – Neutral |
| beckmancoulter_devices_samplestore_dll.BeckmanCoulter.Devices.SampleStore.OAuth.OAuthTokenGenerator._3CRefreshToken_3Ed__2 |          |             |                     |        |                    |
| MoveNext 23%                                                                                                               | 3        | 3           | Information Leakage | 918    | Neutral            |
| beckmancoulter_devices_samplestore_dll.BeckmanCoulter.Devices.SampleStore.OAuth.OAuthRequestHandler._3CSendRequest_3Ed_9   |          |             |                     |        |                    |
| MoveNext 60%                                                                                                               | 2        | 5           | Information Leakage | 201    | Unlikely – Neutral |

Module: BECLS.POD.Devices.VolumeCheck.dll Used by 6 executables; score impact: < 1

| Location                 | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|--------------------------|----------|-------------|---------------------|--------|----------------|
| readvolumesjob.cs 102    | 3        | 2           | Directory Traversal | 73     | Neutral        |
| volumecheckdevice.cs 163 | 3        | 2           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.Echo.Common.dll Used by 100 executables; score impact: < 1

| Location                                                      | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|---------------------------------------------------------------|----------|-------------|---------------------|--------|----------------|
| labcyte_echo_common_dll.Labcyte.Echo.Common.Output.OutputData |          |             |                     |        |                |
| b 9%                                                          | 0        | 41          | Code Quality        | 404    | Neutral        |
| SaveExceptionsXML 79%                                         | 0        | 41          | Code Quality        | 404    | Neutral        |
| a 96%                                                         | 3        | 41          | Directory Traversal | 73     | Neutral        |
| SaveXML 31%                                                   | 3        | 20          | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.ASM.dll Used by 8 executables; score impact: < 1

| Location               | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|------------------------|----------|-------------|---------------------|--------|----------------|
| remoteconnector.cs 315 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| remoteconnector.cs 355 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.BMG.dll Used by 6 executables; score impact: < 1

| Location               | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|------------------------|----------|-------------|---------------------|--------|----------------|
| remoteconnector.cs 331 | 3        | 2           | Directory Traversal | 73     | Neutral        |
| remoteconnector.cs 366 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.BioShake.dll Used by 6 executables; score impact: < 1

| Location              | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|-----------------------|----------|-------------|---------------------|--------|----------------|
| bioshakedevice.cs 660 | 3        | 3           | Directory Traversal | 73     | Neutral        |
| bioshakedevice.cs 727 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.EnvisionReader.dll Used by 7 executables; score impact: &lt; 1

| Location               | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|------------------------|----------|-------------|---------------------|--------|----------------|
| remoteconnector.cs 325 | 3        | 2           | Directory Traversal | 73     | Neutral        |
| remoteconnector.cs 371 | 3        | 4           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Devices.Teleshake.dll Used by 7 executables; score impact: &lt; 1

| Location               | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|------------------------|----------|-------------|---------------------|--------|----------------|
| teleshakedevice.cs 139 | 3        | 4           | Directory Traversal | 73     | Neutral        |
| teleshakedevice.cs 214 | 3        | 4           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Driver.AmbiStore.dll Used by 10 executables; score impact: &lt; 1

| Location                | Severity | # Instances | Flaw Category        | CWE ID | Exploitability |
|-------------------------|----------|-------------|----------------------|--------|----------------|
| ambistoredriver.cs 526  | 3        | 4           | Server Configuration | 441    | Neutral        |
| ambistoredriver2.cs 294 | 3        | 3           | Server Configuration | 441    | Neutral        |

Module: Labcyte.POD.StorageProtocol.dll Used by 3 executables; score impact: &lt; 1

| Location                 | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|--------------------------|----------|-------------|---------------------|--------|----------------|
| runoptionscontrol.cs 477 | 3        | 2           | Directory Traversal | 73     | Neutral        |
| runoptionscontrol.cs 482 | 3        | 2           | Directory Traversal | 73     | Neutral        |

Module: BeckmanCoulter.Devices.VolumeCheck.dll Used by 8 executables; score impact: &lt; 1

| Location                                                                                       | Severity | # Instances | Flaw Category        | CWE ID | Exploitability |
|------------------------------------------------------------------------------------------------|----------|-------------|----------------------|--------|----------------|
| beckmancoulter_devices_volumecheck_dll.BeckmanCoulter.Devices.VolumeCheck.SimVolumeCheckDevice |          |             |                      |        |                |
| GetPlateVolumes 81%                                                                            | 3        | 4           | Cryptographic Issues | 331    | Unlikely       |

Module: Labcyte.POD.Devices.SampleStore.dll Used by 6 executables; score impact: &lt; 1

| Location                   | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|----------------------------|----------|-------------|---------------------|--------|----------------|
| devicerundefinition.cs 300 | 3        | 3           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.Driver.MicroServe.dll Used by 9 executables; score impact: &lt; 1

| Location                | Severity | # Instances | Flaw Category        | CWE ID | Exploitability |
|-------------------------|----------|-------------|----------------------|--------|----------------|
| microservedriver.cs 442 | 3        | 4           | Server Configuration | 441    | Neutral        |

Module: Labcyte.POD.Driver.StoreX.dll Used by 9 executables; score impact: &lt; 1

| Location                  | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|---------------------------|----------|-------------|---------------------|--------|----------------|
| tcpclientsocketsim.cs 115 | 3        | 1           | Directory Traversal | 73     | Neutral        |

Module: Labcyte.POD.PlateStorageControl.dll Used by 3 executables; score impact: &lt; 1

| Location                 | Severity | # Instances | Flaw Category       | CWE ID | Exploitability |
|--------------------------|----------|-------------|---------------------|--------|----------------|
| platestoragegrid.cs 8933 | 3        | 2           | Directory Traversal | 73     | Neutral        |

Module: AGCSN.dll Used by 2 executables; score impact: &lt; 1

| Location                       | Severity | # Instances | Flaw Category | CWE ID | Exploitability |
|--------------------------------|----------|-------------|---------------|--------|----------------|
| clsmousekeyboardevents.cs 2848 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2849 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2857 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2858 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2866 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2867 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2875 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2876 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2884 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2885 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2893 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2894 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2902 | 0        | 1           | Code Quality  | 404    | Neutral        |
| clsmousekeyboardevents.cs 2903 | 0        | 1           | Code Quality  | 404    | Neutral        |

Module: Labcyte.POD.Base.dll Used by 113 executables; score impact: &lt; 1

| Location       | Severity | # Instances | Flaw Category      | CWE ID | Exploitability |
|----------------|----------|-------------|--------------------|--------|----------------|
| lickey.cs 1184 | 0        | 23          | Potential Backdoor | 656    | Neutral        |

Module: Labcyte.POD.Database.dll Used by 97 executables; score impact: &lt; 1

| Location                | Severity | # Instances | Flaw Category | CWE ID | Exploitability |
|-------------------------|----------|-------------|---------------|--------|----------------|
| poddb.designer.cs 27298 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 30720 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 31871 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 32188 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 33874 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 34359 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 34642 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 34957 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 35510 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 36168 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 37785 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 38175 | 0        | 21          | Code Quality  | 404    | Neutral        |

| Location                | Severity | # Instances | Flaw Category | CWE ID | Exploitability |
|-------------------------|----------|-------------|---------------|--------|----------------|
| poddb.designer.cs 38584 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 38958 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 39242 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 39557 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 40039 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 40492 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 40999 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 41568 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 42183 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 42762 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 43176 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 43856 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 44300 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 44811 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 45107 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 45619 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 46013 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 46340 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 46962 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 47386 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 47825 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 48073 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 48648 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 49226 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 50650 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 51984 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 52395 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 52919 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 53658 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 53985 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 54357 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 54962 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 56232 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 56596 | 0        | 21          | Code Quality  | 404    | Neutral        |
| poddb.designer.cs 57175 | 0        | 21          | Code Quality  | 404    | Neutral        |

Module: Labcyte.POD.Devices.GX.dll Used by 6 executables; score impact: < 1

| Location        | Severity | # Instances | Flaw Category | CWE ID | Exploitability |
|-----------------|----------|-------------|---------------|--------|----------------|
| gxdevice.cs 531 | 0        | 3           | Code Quality  | 404    | Neutral        |
| gxdevice.cs 542 | 0        | 3           | Code Quality  | 404    | Neutral        |

Module: Labcyte.POD.Driver.Koyo.dll Used by 8 executables; score impact: < 1

| Location                    | Severity | # Instances | Flaw Category | CWE ID | Exploitability |
|-----------------------------|----------|-------------|---------------|--------|----------------|
| modbusserialtransport.cs 30 | 0        | 4           | Code Quality  | 404    | Neutral        |

## About Veracode's Methodology

The Veracode platform uses static and dynamic analysis (for web applications) to identify software security flaws in your applications. Using both static and dynamic analysis helps reduce false negatives and detect a broader range of security flaws. Veracode static analysis models the application into an intermediate representation, which is then analyzed for security flaws using a set of automated security tests. Dynamic analysis uses an automated penetration testing technique to detect security flaws at runtime. Once the automated process is complete, a security technician verifies the output to ensure the lowest false positive rates in the industry. The end result is an accurate list of security flaws for the classes of automated scans applied to the application.

## Veracode Rating System Using Multiple Analysis Techniques

Higher assurance applications require more comprehensive analysis to accurately score their security quality. Because each analysis technique (automated static, automated dynamic, manual penetration testing or manual review) has differing false negative (FN) rates for different types of security flaws, any single analysis technique or even combination of techniques is bound to produce a certain level of false negatives. Some false negatives are acceptable for lower business critical applications, so a less expensive analysis using only one or two analysis techniques is acceptable. At higher business criticality the FN rate should be close to zero, so multiple analysis techniques are recommended.

## Application Security Policies

The Veracode platform allows an organization to define and enforce a uniform application security policy across all applications in its portfolio. The elements of an application security policy include the target Veracode Level for the application; types of flaws that should not be in the application (which may be defined by flaw severity, flaw category, CWE, or a common standard including OWASP, CWE/SANS Top 25, or PCI); minimum Veracode security score; required scan types and frequencies; and grace period within which any policy-relevant flaws should be fixed.

### Policy constraints

Policies have three main constraints that can be applied: rules, required scans, and remediation grace periods.

### Evaluating applications against a policy

When an application is evaluated against a policy, it can receive one of four assessments:

**Not assessed** The application has not yet had a scan published

**Passed** The application has passed all the aspects of the policy, including rules, required scans, and grace period.

**Did not pass** The application has not completed all required scans; has not achieved the target Veracode Level; or has one or more policy relevant flaws that have exceeded the grace period to fix.

**Conditional pass** The application has one or more policy relevant flaws that have not yet exceeded the grace period to fix.

## Understand Veracode Levels

The Veracode Level (VL) achieved by an application is determined by type of testing performed on the application, and the severity and types of flaws detected. A minimum security score (defined below) is also required for each level.

There are five Veracode Levels denoted as VL1, VL2, VL3, VL4, and VL5. VL1 is the lowest level and is achieved by demonstrating that security testing, automated static or dynamic, is utilized during the SDLC. VL5 is the highest level and is achieved by performing automated and manual testing and removing all significant flaws. The Veracode Levels VL2, VL3, and VL4 form a continuum of increasing software assurance between VL1 and VL5.

For IT staff operating applications, Veracode Levels can be used to set application security policies. For deployment scenarios of different business criticality, differing VLs should be made requirements. For example, the policy for applications that handle credit card transactions, and therefore have PCI compliance requirements, should be VL5. A medium business criticality internal application could have a policy requiring VL3.

Software developers can decide which VL they want to achieve based on the requirements of their customers. Developers of software that is mission critical to most of their customers will want to achieve VL5. Developers of general purpose business software may want

to achieve VL3 or VL4. Once the software has achieved a Veracode Level it can be communicated to customers through a Veracode Report or through the Veracode Directory on the Veracode web site.

### Criteria for achieving Veracode Levels

The following table defines the details to achieve each Veracode Level. The criteria for all columns: Flaw Severities Not Allowed, Flaw Categories not Allowed, Testing Required, and Minimum Score.

\*Dynamic is only an option for web applications.

| Veracode Level | Flaw Severities Not Allowed | Testing Required*           | Minimum Score |
|----------------|-----------------------------|-----------------------------|---------------|
| VL5            | V.High, High, Medium        | Static AND Manual           | 90            |
| VL4            | V.High, High, Medium        | Static                      | 80            |
| VL3            | V.High, High                | Static                      | 70            |
| VL2            | V.High                      | Static OR Dynamic OR Manual | 60            |
| VL1            |                             | Static OR Dynamic OR Manual |               |

When multiple testing techniques are used it is likely that not all testing will be performed on the exact same build. If that is the case the latest test results from a particular technique will be used to calculate the current Veracode Level. After 6 months test results will be deemed out of date and will no longer be used to calculate the current Veracode Level.

### Business Criticality

The foundation of the Veracode rating system is the concept that more critical applications require higher security quality scores to be acceptable risks. Less business critical applications can tolerate lower security quality. The business criticality is dictated by the typical deployed environment and the value of data used by the application. Factors that determine business criticality are: reputation damage, financial loss, operational risk, sensitive information disclosure, personal safety, and legal violations.

US. Govt. OMB Memorandum M-04-04; NIST FIPS Pub. 199

Business Criticality Description

|           |                                                                                               |
|-----------|-----------------------------------------------------------------------------------------------|
| Very High | Mission critical for business/safety of life and limb on the line                             |
| High      | Exploitation causes serious brand damage and financial loss with long term business impact    |
| Medium    | Applications connected to the internet that process financial or private customer information |
| Low       | Typically internal applications with non-critical business impact                             |
| Very Low  | Applications with no material business impact                                                 |

### Business Criticality Definitions

**Very High (BC5)** This is typically an application where the safety of life or limb is dependent on the system; it is mission critical the application maintain 100% availability for the long term viability of the project or business. Examples are control software for industrial, transportation or medical equipment or critical business systems such as financial trading systems.

**High (BC4)** This is typically an important multi-user business application reachable from the internet and is critical that the application maintain high availability to accomplish its mission. Exploitation of high criticality applications cause serious brand damage and business/financial loss and could lead to long term business impact.

**Medium (BC3)** This is typically a multi-user application connected to the internet or any system that processes financial or private customer information. Exploitation of medium criticality applications typically result in material business impact resulting

in some financial loss, brand damage or business liability. An example is a financial services company's internal 401K management system.

**Low (BC2)** This is typically an internal only application that requires low levels of application security such as authentication to protect access to non-critical business information and prevent IT disruptions. Exploitation of low criticality applications may lead to minor levels of inconvenience, distress or IT disruption. An example internal system is a conference room reservation or business card order system.

**Very Low (BC1)** Applications that have no material business impact should its confidentiality, data integrity and availability be affected. Code security analysis is not required for applications at this business criticality, and security spending should be directed to other higher criticality applications.

## Scoring Methodology

The Veracode scoring system, Security Quality Score, is built on the foundation of two industry standards, the Common Weakness Enumeration (CWE) and Common Vulnerability Scoring System (CVSS). CWE provides the dictionary of security flaws and CVSS provides the foundation for computing severity, based on the potential Confidentiality, Integrity and Availability impact of a flaw if exploited.

The Security Quality Score is a single score from 0 to 100, where 0 is the most insecure application and 100 is an application with no detectable security flaws. The score calculation includes non-linear factors so that, for instance, a single Severity 5 flaw is weighted more heavily than five Severity 1 flaws, and so that each additional flaw at a given severity contributes progressively less to the score.

Veracode assigns a severity level to each flaw type based on three foundational application security requirements — Confidentiality, Integrity and Availability. Each of the severity levels reflects the potential business impact if a security breach occurs across one or more of these security dimensions.

### Confidentiality Impact

According to CVSS, this metric measures the impact on confidentiality if an exploit should occur using the vulnerability on the target system. At the weakness level, the scope of the Confidentiality in this model is within an application and is measured at three levels of impact -None, Partial and Complete.

### Integrity Impact

This metric measures the potential impact on integrity of the application being analyzed. Integrity refers to the trustworthiness and guaranteed veracity of information within the application. Integrity measures are meant to protect data from unauthorized modification. When the integrity of a system is sound, it is fully proof from unauthorized modification of its contents.

### Availability Impact

This metric measures the potential impact on availability if a successful exploit of the vulnerability is carried out on a target application. Availability refers to the accessibility of information resources. Almost exclusive to this domain are denial-of-service vulnerabilities. Attacks that compromise authentication and authorization for application access, application memory, and administrative privileges are examples of impact on the availability of an application.

## Security Quality Score Calculation

The overall Security Quality Score is computed by aggregating impact levels of all weaknesses within an application and representing the score on a 100 point scale. This score does not predict vulnerability potential as much as it enumerates the security weaknesses and their impact levels within the application code.

The Raw Score formula puts weights on each flaw based on its impact level. These weights are exponential and determined by empirical analysis by Veracode's application security experts with validation from industry experts. The score is normalized to a scale of 0 to 100, where a score of 100 is an application with 0 detected flaws using the analysis technique for the application's business criticality.

## Understand Severity, Exploitability, and Remediation Effort

Severity and exploitability are two different measures of the seriousness of a flaw. Severity is defined in terms of the potential impact to confidentiality, integrity, and availability of the application as defined in the CVSS, and exploitability is defined in terms of the likelihood

or ease with which a flaw can be exploited. A high severity flaw with a high likelihood of being exploited by an attacker is potentially more dangerous than a high severity flaw with a low likelihood of being exploited.

Remediation effort, also called Complexity of Fix, is a measure of the likely effort required to fix a flaw. Together with severity, the remediation effort is used to give Fix First guidance to the developer.

### Veracode Flaw Severities

Veracode flaw severities are defined as follows:

| Severity      | Description                                                                                                                                                                                                          |
|---------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Very High     | The offending line or lines of code is a very serious weakness and is an easy target for an attacker. The code should be modified immediately to avoid potential attacks.                                            |
| High          | The offending line or lines of code have significant weakness, and the code should be modified immediately to avoid potential attacks.                                                                               |
| Medium        | A weakness of average severity. These should be fixed in high assurance software. A fix for this weakness should be considered after fixing the very high and high for medium assurance software.                    |
| Low           | This is a low priority weakness that will have a small impact on the security of the software. Fixing should be consideration for high assurance software. Medium and low assurance software can ignore these flaws. |
| Very Low      | Minor problems that some high assurance software may want to be aware of. These flaws can be safely ignored in medium and low assurance software.                                                                    |
| Informational | Issues that have no impact on the security quality of the application but which may be of interest to the reviewer.                                                                                                  |

#### Informational findings

Informational severity findings are items observed in the analysis of the application that have no impact on the security quality of the application but may be interesting to the reviewer for other reasons. These findings may include code quality issues, API usage, and other factors.

Informational severity findings have no impact on the security quality score of the application and are not included in the summary tables of flaws for the application.

### Exploitability

Each flaw instance in a static scan may receive an exploitability rating. The rating is an indication of the intrinsic likelihood that the flaw may be exploited by an attacker. Veracode recommends that the exploitability rating be used to prioritize flaw remediation within a particular group of flaws with the same severity and difficulty of fix classification.

The possible exploitability ratings include:

| Exploitability | Description                   |
|----------------|-------------------------------|
| V. Unlikely    | Very unlikely to be exploited |
| Unlikely       | Unlikely to be exploited      |

| Exploitability | Description                                  |
|----------------|----------------------------------------------|
| Neutral        | Neither likely nor unlikely to be exploited. |
| Likely         | Likely to be exploited                       |
| V. Likely      | Very likely to be exploited                  |

Note: All reported flaws found via dynamic scans are assumed to be exploitable, because the dynamic scan actually executes the attack in question and verifies that it is valid.

Effort/Complexity of Fix

Each flaw instance receives an effort/complexity of fix rating based on the classification of the flaw. The effort/complexity of fix rating is given on a scale of 1 to 5, as follows:

| Effort/Complexity of Fix | Description                                                                             |
|--------------------------|-----------------------------------------------------------------------------------------|
| 5                        | Complex design error. Requires significant redesign.                                    |
| 4                        | Simple design error. Requires redesign and up to 5 days to fix.                         |
| 3                        | Complex implementation error. Fix is approx. 51-500 lines of code. Up to 5 days to fix. |
| 2                        | Implementation error. Fix is approx. 6-50 lines of code. 1 day to fix.                  |
| 1                        | Trivial implementation error. Fix is up to 5 lines of code. One hour or less to fix.    |

Flaw Types by Severity Level

The flaw types by severity level table provides a summary of flaws found in the application by Severity and Category. The table puts the Security Quality Score into context by showing the specific breakout of flaws by severity, used to compute the score as described above. If multiple analysis techniques are used, the table includes a breakout of all flaws by category and severity for each analysis type performed.

Flaws by Severity

The flaws by severity chart shows the distribution of flaws by severity. An application can get a mediocre security rating by having a few high risk flaws or many medium risk flaws.

Flaws in Common Modules

The flaws in common modules listing shows a summary of flaws in shared dependency modules in this application. A shared dependency is a dependency that is used by more than one analyzed module. Each module is listed with the number of executables that consume it as a dependency and a summary of the impact on the application's security score of the flaws found in the dependency.

The score impact represents the amount that the application score would increase if all the flaws in the shared dependency module were fixed. This information can be used to focus remediation efforts on common modules with a higher impact on the application security score.

Only common modules that were uploaded with debug information are included in the Flaws in Common Modules listing.

## Action Items

The Action Items section of the report provides guidance on the steps required to bring the application to a state where it passes its assigned policy. These steps may include fixing or mitigating flaws or performing additional scans. The section also includes best practice recommendations to improve the security quality of the application.

## Common Weakness Enumeration (CWE)

The Common Weakness Enumeration (CWE) is an industry standard classification of types of software weaknesses, or flaws, that can lead to security problems. CWE is widely used to provide a standard taxonomy of software errors. Every flaw in a Veracode report is classified according to a standard CWE identifier.

More guidance and background about the CWE is available at <http://cwe.mitre.org/data/index.html>.

## About Manual Assessments

The Veracode platform can include the results from a manual assessment (usually a penetration test or code review) as part of a report. These results differ from the results of automated scans in several important ways, including objectives, attack vectors, and common attack patterns.

A manual penetration assessment is conducted to observe the application code in a run-time environment and to simulate real-world attack scenarios. Manual testing is able to identify design flaws, evaluate environmental conditions, compound multiple lower risk flaws into higher risk vulnerabilities, and determine if identified flaws affect the confidentiality, integrity, or availability of the application.

### Objectives

The stated objectives of a manual penetration assessment are:

- Perform testing, using proprietary and/or public tools, to determine whether it is possible for an attacker to:
- Circumvent authentication and authorization mechanisms
- Escalate application user privileges
- Hijack accounts belonging to other users
- Violate access controls placed by the site administrator
- Alter data or data presentation
- Corrupt application and data integrity, functionality and performance
- Circumvent application business logic
- Circumvent application session management
- Break or analyze use of cryptography within user accessible components
- Determine possible extent access or impact to the system by attempting to exploit vulnerabilities
- Score vulnerabilities using the Common Vulnerability Scoring System (CVSS)
- Provide tactical recommendations to address security issues of immediate consequence

Provide strategic recommendations to enhance security by leveraging industry best practices

### Attack vectors

In order to achieve the stated objectives, the following tests are performed as part of the manual penetration assessment, when applicable to the platforms and technologies in use:

- Cross Site Scripting (XSS)
- SQL Injection
- Command Injection
- Cross Site Request Forgery (CSRF)
- Authentication/Authorization Bypass
- Session Management testing, e.g. token analysis, session expiration, and logout effectiveness
- Account Management testing, e.g. password strength, password reset, account lockout, etc.
- Directory Traversal
- Response Splitting
- Stack/Heap Overflows
- Format String Attacks

- Cookie Analysis
- Server Side Includes Injection
- Remote File Inclusion
- LDAP Injection
- XPATH Injection
- Internationalization attacks
- Denial of Service testing at the application layer only
- AJAX Endpoint Analysis
- Web Services Endpoint Analysis
- HTTP Method Analysis
- SSL Certificate and Cipher Strength Analysis
- Forced Browsing

## CAPEC Attack Pattern Classification

The following attack pattern classifications are used to group similar application flaws discovered during manual penetration testing. Attack patterns describe the general methods employed to access and exploit the specific weaknesses that exist within an application. CAPEC (Common Attack Pattern Enumeration and Classification) is an effort led by Cigital, Inc. and is sponsored by the United States Department of Homeland Security's National Cyber Security Division.

## Abuse of Functionality

Exploitation of business logic errors or misappropriation of programmatic resources. Application functions are developed to specifications with particular intentions, and these types of attacks serve to undermine those intentions.

Examples:

- Exploiting password recovery mechanisms
- Accessing unpublished or test APIs
- Cache poisoning

## Spoofing

Impersonation of entities or trusted resources. A successful attack will present itself to a verifying entity with an acceptable level of authenticity.

Examples:

- Man in the middle attacks
- Checksum spoofing
- Phishing attacks

## Probabilistic Techniques

Using predictive capabilities or exhaustive search techniques in order to derive or manipulate sensitive information. Attacks capitalize on the availability of computing resources or the lack of entropy within targeted components.

Examples:

- Password brute forcing
- Cryptanalysis
- Manipulation of authentication tokens

## Exploitation of Authentication

Circumventing authentication requirements to access protected resources. Design or implementation flaws may allow authentication checks to be ignored, delegated, or bypassed.

Examples:

- Cross-site request forgery
- Reuse of session identifiers
- Flawed authentication protocol

## Resource Depletion

Affecting the availability of application components or resources through symmetric or asymmetric consumption. Unrestricted access to computationally expensive functions or implementation flaws that affect the stability of the application can be targeted by an attacker in order to cause denial of service conditions.

Examples:

- Flooding attacks
- Unlimited file upload size
- Memory leaks

## Exploitation of Privilege/Trust

Undermining the application's trust model in order to gain access to protected resources or gain additional levels of access as defined by the application. Applications that implicitly extend trust to resources or entities outside of their direct control are susceptible to attack.

Examples:

- Insufficient access control lists
- Circumvention of client side protections
- Manipulation of role identification information

## Injection

Inserting unexpected inputs to manipulate control flow or alter normal business processing. Applications must contain sufficient data validation checks in order to sanitize tainted data and prevent malicious, external control over internal processing.

Examples:

- SQL Injection
- Cross-site scripting
- XML Injection

## Data Structure Attacks

Supplying unexpected or excessive data that results in more data being written to a buffer than it is capable of holding. Successful attacks of this class can result in arbitrary command execution or denial of service conditions.

Examples:

- Buffer overflow
- Integer overflow
- Format string overflow

## Data Leakage Attacks

Recovering information exposed by the application that may itself be confidential or may be useful to an attacker in discovering or exploiting other weaknesses. A successful attack may be conducted passive observation or active interception methods. This attack pattern often manifests itself in the form of applications that expose sensitive information within error messages.

Examples:

- Sniffing clear-text communication protocols
- Stack traces returned to end users
- Sensitive information in HTML comments

## Resource Manipulation

Manipulating application dependencies or accessed resources in order to undermine security controls and gain unauthorized access to protected resources. Applications may use tainted data when constructing paths to local resources or when constructing processing environments.

Examples:

- Carriage Return Line Feed log file injection
- File retrieval via path manipulation
- User specification of configuration files

## Time and State Attacks

Undermining state condition assumptions made by the application or capitalizing on time delays between security checks and performed operations. An application that does not enforce a required processing sequence or does not handle concurrency adequately will be susceptible to these attack patterns.

Examples:

- Bypassing intermediate form processing steps
- Time-of-check and time-of-use race conditions
- Deadlock triggering to cause a denial of service

## Terms of Use

Use and distribution of this report are governed by the agreement between Veracode and its customer. In particular, this report and the results in the report cannot be used publicly in connection with Veracode's name without written permission.

Appendix A: Changes from Last Scan

| Latest Scan        |                             | Prior Scan |                              |
|--------------------|-----------------------------|------------|------------------------------|
| <b>Static Scan</b> |                             |            |                              |
| Scan Name:         | Version:3.2 Build Number:22 | Scan Name: | Tempo Develop Build_20260223 |
| Completed:         | 8/12/26                     | Completed: | 2/22/26                      |
| Score:             | 78                          | Score:     | 65                           |

Changes in scope of scan (static)

New Modules

| Module Name        | Compiler        | Operating Environment | Engine Version |
|--------------------|-----------------|-----------------------|----------------|
| PlateLocClient.exe | MSIL_MSVC14_X86 | Windows               | 20260806120029 |

Removed modules

| Module Name                    | Compiler        | Operating Environment | Engine Version |
|--------------------------------|-----------------|-----------------------|----------------|
| Labcyte.POD.Driver.BigBear.dll | MSIL_MSVC14_X86 | Windows               | 20260219223424 |
| LightCycler1536Client.exe      | MSIL_MSVC14_X86 | Windows               | 20260219223424 |
| LightCycler480Client.exe       | MSIL_MSVC14_X86 | Windows               | 20260219223424 |
| RobotTeachClient.dll           | MSIL_MSVC14_X86 | Windows               | 20260219223424 |

Flaws not detected in current scan

The following is a list of all flaws found in the prior scan of this application that were not detected in the current scan.

Very High (1 flaw)  **Fix Required by Policy**

➔ Command or Argument Injection(1 flaw)

Associated Flaws by CWE ID:

➔ Improper Neutralization of Special Elements used in an OS Command ('OS Command Injection') (CWE ID 78)(1 flaw)

Instances found via Static Scan

| Flaw Id                                                                                 | Module # | Class # | Module                                     | Location                           | Fix By |
|-----------------------------------------------------------------------------------------|----------|---------|--------------------------------------------|------------------------------------|--------|
|  383 | 53       | -       | koyopl.exe#9.9.9.0 /labcyte.pod.common.dll | github/.../device/poddevice.cs 771 |        |

High (2 flaws)

→ SQL Injection(2 flaws)

Associated Flaws by CWE ID:

→ Improper Neutralization of Special Elements used in an SQL Command ('SQL Injection') (CWE ID 89)(2 flaws)

Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                   | Location                            | Fix By |
|---------|----------|---------|------------------------------------------|-------------------------------------|--------|
| 626     | 51       | -       | tempo.exe#9.9.9.0/labcyte.pod.common.dll | github/.../platestoreengine.cs 8313 |        |
| 627     | 51       | -       | tempo.exe#9.9.9.0/labcyte.pod.common.dll | github/.../platestoreengine.cs 9887 |        |

Medium (30 flaws)

 Fix Required by Policy

→ CRLF Injection(16 flaws)

Associated Flaws by CWE ID:

→ Improper Output Neutralization for Logs (CWE ID 117)(16 flaws)

Instances found via Static Scan

| Flaw Id                                                                                | Module # | Class # | Module                                                                                         | Location                              | Fix By |
|----------------------------------------------------------------------------------------|----------|---------|------------------------------------------------------------------------------------------------|---------------------------------------|--------|
|  22 | -        | 27      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | bool CheckOrderFile(string) 76%       |        |
|  21 | -        | 27      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | bool CheckRunTemplateFile(string) 76% |        |
|  20 | -        | 27      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | bool CheckSampleList(string) 76%      |        |
|  7  | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | int Login(string, string) 43%         |        |
|  4  | -        | 25      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | int SetConnection(string, int) 22%    |        |

| Flaw Id | Module # | Class # | Module                                                                                         | Location                                            | Fix By |
|---------|----------|---------|------------------------------------------------------------------------------------------------|-----------------------------------------------------|--------|
| 5       | -        | 25      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | int SetConnection(string, int) 91%                  |        |
| 17      | -        | 27      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | string CreateOutputPath(string) 92%                 |        |
| 2       | -        | 24      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void LogMethod(string[]) 92%                        |        |
| 3       | -        | 24      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void LogResult(int, string[]) 93%                   |        |
| 9       | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 43% |        |
| 10      | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 45% |        |
| 12      | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 73% |        |
| 13      | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 76% |        |
| 15      | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 92% |        |
| 16      | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 95% |        |
| 24      | -        | 27      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void ValidateSchema(string, string) 87%             |        |

→ Credentials Management(1 flaw)

Associated Flaws by CWE ID:

→ Use of Hard-coded Password (CWE ID 259)(1 flaw)

Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                                         | Location                   | Fix By |
|---------|----------|---------|------------------------------------------------------------------------------------------------|----------------------------|--------|
| 6       | -        | 25      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | string Encrypt(string) 26% |        |

→ Directory Traversal(13 flaws\*)

Associated Flaws by CWE ID:

→ External Control of File Name or Path (CWE ID 73)(13 flaws\*)

Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                                         | Location                              | Fix By |
|---------|----------|---------|------------------------------------------------------------------------------------------------|---------------------------------------|--------|
| 19      | -        | 27      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | bool OrderFileHasPlateIDs(string) 19% |        |
| 617     | 29       | -       | tempo.exe#9.9.9.0/labcyte.pod.devices.intellixcap.dll                                          | .../intellixcapdevice.cs 448          |        |
| 611     | 38       | -       | tempo.exe#9.9.9.0/labcyte.pod.devices.lightcycler.dll                                          | .../lightcyclerdevice.cs 1165         |        |
| 612     | 38       | -       | tempo.exe#9.9.9.0/labcyte.pod.devices.lightcycler.dll                                          | .../lightcyclerdevice.cs 1168         |        |
| 613     | 38       | -       | tempo.exe#9.9.9.0/labcyte.pod.devices.lightcycler.dll                                          | .../lightcyclerdevice.cs 1213         |        |
| 615     | 38       | -       | tempo.exe#9.9.9.0/labcyte.pod.devices.lightcycler.dll                                          | .../lightcyclerdevice.cs 1266         |        |
| 592     | 73       | -       | tempo.exe#9.9.9.0/labcyte.pod.api.dll/labcyte.pod.devices.lightcycler1536.dll                  | .../run1536experimentjob.cs 676       |        |
| 593     | 73       | -       | tempo.exe#9.9.9.0/labcyte.pod.api.dll/labcyte.pod.devices.lightcycler1536.dll                  | .../run1536experimentjob.cs 676       |        |
| 18      | -        | 27      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | string TryCreateFolders(string[]) 45% |        |

| Flaw Id | Module # | Class # | Module                                                                                         | Location                                            | Fix By |
|---------|----------|---------|------------------------------------------------------------------------------------------------|-----------------------------------------------------|--------|
|         |          |         | e.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll                                         |                                                     |        |
| 8       | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 36% |        |
| 11      | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 68% |        |
| 14      | -        | 26      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void SaveExperiment(Simulator.Experiment, bool) 87% |        |
| 23      | -        | 27      | lightcycler1536client.exe#9.9.9.0/labcyte.pod.driver.lightcycler.dll/roche.htc1536.aim.api.dll | void ValidateSchema(string, string) 33%             |        |

### Low (1 flaw)

#### → Information Leakage(1 flaw)

Associated Flaws by CWE ID:

#### → Insertion of Sensitive Information Into Sent Data (CWE ID 201)(1 flaw)

#### Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                                                           | Location            | Fix By |
|---------|----------|---------|------------------------------------------------------------------------------------------------------------------|---------------------|--------|
| 498     | -        | 2       | tempo.exe#9.9.9.0/labcyte.pod.api.dll/labcyte.pod.devices.samplestore.dll/beckmancoulter.devices.samplestore.dll | void MoveNext() 27% |        |

### Info (1 flaw)

#### → Code Quality(1 flaw)

Associated Flaws by CWE ID:

→ Improper Resource Shutdown or Release (CWE ID 404)(1 flaw)

Instances found via Static Scan

| Flaw Id | Module # | Class # | Module                                                                                | Location                  | Fix By |
|---------|----------|---------|---------------------------------------------------------------------------------------|---------------------------|--------|
| 784     | 37       | -       | ambistoreclient.exe<br>#9.9.9.0/labcyte.pod.driver.ambistore.dll/labcyte.pod.base.dll | github/.../lickey.cs 1167 |        |

Appendix B: Approved Mitigated Flaws (by Danaher Corporation)

NOTE: Except in circumstances where the Customer has purchased Veracode's Mitigation Proposal Review Solution, Veracode does not review the mitigation strategy described below and is not responsible for its contents or the accuracy of any statements provided.

Very High (2 flaws)

Fix Required by Policy:  Flaw no longer impacts results.  
 Flaw continues to impact results.

→ Command or Argument Injection(2 flaws)

Associated Flaws by CWE ID:

→ Improper Neutralization of Special Elements used in an OS Command ('OS Command Injection') (CWE ID 78)(2 flaws)

Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                               | Exploitability | Module # | Class # | Module                                                          | Location                         |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-----------------------------------------------------------------|----------------------------------|
|  871                                                                                                                                                                                                                                                                                                                                                                                                                                                 | Neutral        | 47       | -       | tempo.exe#3.2.0.1/b<br>ecls.pod.devices.mobi<br>lerobotdoor.dll | .../mobilerobotdoordevice.cs 104 |
| <p>↳ Potential False Positive (Danaher Corporation): The call launches a URL via Process.Start(ProcessStartInfo), no cmd.exe, no argument concatenation, no executable path from input. The value is passed as a single URL target. Before launch Non-URI and non-http(s) payloads are rejected using the validation added. The taint source originates from DB-loaded device config, which maps into PODDeviceInfo.Port. This is operational configuration managed by Tempo admin, not an unauthenticated runtime input channel.</p> |                |          |         |                                                                 |                                  |
| <p>↳ Approve Mitigation (Danaher Corporation): Approve</p>                                                                                                                                                                                                                                                                                                                                                                                                                                                                            |                |          |         |                                                                 |                                  |
|  817                                                                                                                                                                                                                                                                                                                                                                                                                                               | Neutral        | 94       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.userinterfa<br>ce.dll         | .../uccapcheckimages.cs 350      |
| <p>↳ Potential False Positive (Danaher Corporation): Whitelisted allowed extensions. Only image extensions (".jpg", ".jpeg", ".bmp", ".gif", ".png") are allowed in the executable path, no command line arguments passed here. Multiple validations are added. This makes sure that no other harmful files or commands are executed.</p>                                                                                                                                                                                             |                |          |         |                                                                 |                                  |
| <p>↳ Approve Mitigation (Danaher Corporation): Approved. Technically Mitigate by Design.</p>                                                                                                                                                                                                                                                                                                                                                                                                                                          |                |          |         |                                                                 |                                  |

Medium (245 flaws)

Fix Required by Policy:  Flaw no longer impacts results.  
 Flaw continues to impact results.

→ Directory Traversal(233 flaws)

Associated Flaws by CWE ID:

## → External Control of File Name or Path (CWE ID 73)(233 flaws)

### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                         | Exploitability | Module # | Class # | Module                                                      | Location                          |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------------|-----------------------------------|
|  511                                                                                                                                                                           | Neutral        | 1        | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.am<br>bistore.dll | .../ambistoredevice.cs 2021       |
| <span>↳ Mitigate by OS Environment (Danaher Corporation):</span> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                             |                                   |
| <span>↳ Potential False Positive (Danaher Corporation):</span> _logPath is nota direct user input, instead the value is originated from DB                                                                                                                      |                |          |         |                                                             |                                   |
| <span>↳ Approve Mitigation (Danaher Corporation):</span> Approved as documented.                                                                                                                                                                                |                |          |         |                                                             |                                   |
|  512                                                                                                                                                                           | Neutral        | 1        | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.am<br>bistore.dll | .../ambistoredevice.cs 2035       |
| <span>↳ Mitigate by OS Environment (Danaher Corporation):</span> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                             |                                   |
| <span>↳ Potential False Positive (Danaher Corporation):</span> A part of the filename is hardcoded and _logpath is fetched from DB.. Its not a direct UI input                                                                                                  |                |          |         |                                                             |                                   |
| <span>↳ Approve Mitigation (Danaher Corporation):</span> Approved as documented.                                                                                                                                                                                |                |          |         |                                                             |                                   |
|  513                                                                                                                                                                         | Neutral        | 1        | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.am<br>bistore.dll | .../ambistoredevice.cs 2071       |
| <span>↳ Mitigate by OS Environment (Danaher Corporation):</span> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                             |                                   |
| <span>↳ Potential False Positive (Danaher Corporation):</span> A part of the filename is hardcoded and _logpath is fetched from DB.. Its not a direct UI input                                                                                                  |                |          |         |                                                             |                                   |
| <span>↳ Approve Mitigation (Danaher Corporation):</span> Approved as documented.                                                                                                                                                                                |                |          |         |                                                             |                                   |
|  956                                                                                                                                                                         | Neutral        | 4        | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.userinterfa<br>ce.dll     | buildagent/.../barcodeform.cs 256 |
| <span>↳ Potential False Positive (Danaher Corporation):</span> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                             |                                   |
| <span>↳ Approve Mitigation (Danaher Corporation):</span> Approved                                                                                                                                                                                               |                |          |         |                                                             |                                   |
|  955                                                                                                                                                                         | Neutral        | 4        | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.userinterfa<br>ce.dll     | buildagent/.../barcodeform.cs 261 |
| <span>↳ Potential False Positive (Danaher Corporation):</span> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server                                |                |          |         |                                                             |                                   |

| Flaw Id                                                                                                                                                                                                                                 | Exploitability | Module # | Class # | Module                                                                                                      | Location                                                                           |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|
| path involved here.                                                                                                                                                                                                                     |                |          |         |                                                                                                             |                                                                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                                             |                                                                                    |
| 930                                                                                                                                                                                                                                     | Neutral        | 5        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.bioshake.dll                                      | .../bioshakedevice.cs 660                                                          |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                             |                                                                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                                             |                                                                                    |
| 929                                                                                                                                                                                                                                     | Neutral        | 5        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.bioshake.dll                                      | .../bioshakedevice.cs 727                                                          |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                             |                                                                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                                             |                                                                                    |
| 905                                                                                                                                                                                                                                     | Neutral        | -        | 22      | ambistoreclient.exe#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.echo.common.dll | bool a(string, string) 96%                                                         |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                             |                                                                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                                             |                                                                                    |
| 942                                                                                                                                                                                                                                     | Neutral        | -        | 19      | tempo.exe#3.2.0.1/labcyte.echo.api.dll                                                                      | bool b(string) 5%                                                                  |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                             |                                                                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                                             |                                                                                    |
| 906                                                                                                                                                                                                                                     | Neutral        | -        | 22      | ambistoreclient.exe#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.echo.common.dll | bool b(string, string) 9%                                                          |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                             |                                                                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                                             |                                                                                    |
| 962                                                                                                                                                                                                                                     | Neutral        | -        | 12      | tempo.exe#3.2.0.1/labcyte.echo.api.dll                                                                      | bool DetectCSVFormat(string, string, int, System.Collections.Generic.List<ImportHe |

| Flaw Id | Exploitability | Module # | Class # | Module                                                                                                       | Location                                                                                                                                                                                                                                       |
|---------|----------------|----------|---------|--------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|         |                |          |         |                                                                                                              | adInfoltem>, ref string /*1*/) 97%                                                                                                                                                                                                             |
|         |                |          |         |                                                                                                              | ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |
|         |                |          |         |                                                                                                              | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |
| 961     | Neutral        | -        | 12      | tempo.exe#3.2.0.1/labcyte.echo.api.dll                                                                       | bool OpenCSV(string, int, int, int, int, PickListHeaderInfos, ref System.Collections.Generic.Dictionary<int, string> /*1*/, ref string /*1*/) 44%                                                                                              |
|         |                |          |         |                                                                                                              | ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |
|         |                |          |         |                                                                                                              | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |
| 941     | Neutral        | -        | 19      | tempoadmin.exe#3.2.0.1/labcyte.pod.com mon.dll/labcyte.echo.api.dll                                          | bool Save(string) 33%                                                                                                                                                                                                                          |
|         |                |          |         |                                                                                                              | ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |
|         |                |          |         |                                                                                                              | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |
| 940     | Neutral        | -        | 19      | tempoadmin.exe#3.2.0.1/labcyte.pod.com mon.dll/labcyte.echo.api.dll                                          | bool Save(string) 66%                                                                                                                                                                                                                          |
|         |                |          |         |                                                                                                              | ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |
|         |                |          |         |                                                                                                              | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |
| 909     | Neutral        | -        | 22      | ambistoreclient.exe#3.2.0.1/labcyte.pod.d river.ambistore.dll/labcyte.pod.common.dll/labcyte.echo.common.dll | bool SaveExceptionsXML(string) 79%                                                                                                                                                                                                             |
|         |                |          |         |                                                                                                              | ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |
|         |                |          |         |                                                                                                              | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |
| 907     | Neutral        | -        | 22      | ambistoreclient.exe#3.2.0.1/labcyte.pod.d river.ambistore.dll/labcyte.pod.common.dll/labcyte.echo.common.dll | bool SaveXML(string) 31%                                                                                                                                                                                                                       |
|         |                |          |         |                                                                                                              | ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable                                                                                                                              |

| Flaw Id                                                                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                                                    | Location                                       |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------------------------|------------------------------------------------|
| in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                       |                |          |         |                                                                           |                                                |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                                           |                                                |
| 939                                                                                                                                                                                                                                                | Neutral        | -        | 20      | tempo.exe#3.2.0.1/labcyte.echo.api.dll                                    | bool SetSurveyForPlate(string, SurveyData) 65% |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                           |                                                |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                                           |                                                |
| 571                                                                                                                                                                                                                                                | Neutral        | 6        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.temp est.dll    | buildagent/.../bulkfilljob.cs 603              |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                                |
| ↳ Potential False Positive (Danaher Corporation): A part of dispenseListFilePath variable is hardcoded and the devicename is originated from DB. Devicename is allowed to modify by the Admin user alone while setting up Tempo initially.         |                |          |         |                                                                           |                                                |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                           |                                                |
| 572                                                                                                                                                                                                                                                | Neutral        | 6        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.temp est.dll    | buildagent/.../bulkfilljob.cs 658              |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                                |
| ↳ Potential False Positive (Danaher Corporation): A part of filepath variable is hardcoded and the devicename is originated from DB. Devicename is allowed to modify by the Admin user alone while setting up Tempo initially.                     |                |          |         |                                                                           |                                                |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                           |                                                |
| 494                                                                                                                                                                                                                                                | Neutral        | 7        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.certu sflex.dll | .../certusflexdevice.cs 698                    |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                                |
| ↳ Potential False Positive (Danaher Corporation): LastExperimentResultPath is not provided by user, instead fetched programmatically                                                                                                               |                |          |         |                                                                           |                                                |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                           |                                                |
| 495                                                                                                                                                                                                                                                | Neutral        | 7        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.certu sflex.dll | .../certusflexdevice.cs 838                    |

| Flaw Id | Exploitability | Module # | Class # | Module                                                                                                                                                                                                                                                    | Location                  |
|---------|----------------|----------|---------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------|
|         |                |          |         | bctype.pod.api.dll/labcyte.pod.devices.certsflex.dll                                                                                                                                                                                                      |                           |
|         |                |          |         | ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                           |
|         |                |          |         | ↳ <i>Potential False Positive (Danaher Corporation):</i> LastExperimentResultPath is not provided by user, instead fetched programmatically                                                                                                               |                           |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                |                           |
| 491     | Neutral        | 9        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.combinl.dll                                                                                                                                                                                     | .../combinldevice.cs 1106 |
|         |                |          |         | ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                           |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                           |
| 492     | Neutral        | 9        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.combinl.dll                                                                                                                                                                                     | .../combinldevice.cs 1120 |
|         |                |          |         | ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                           |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                           |
| 493     | Neutral        | 9        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.combinl.dll                                                                                                                                                                                     | .../combinldevice.cs 1150 |
|         |                |          |         | ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                           |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                           |
| 663     | Neutral        | 12       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.cytomat.dll                                                                                                                                                                                     | .../cytomatdevice.cs 1732 |
|         |                |          |         | ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                           |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined.                                                                                                                                              |                           |

| Flaw Id                                                                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                                                | Location                  |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-----------------------------------------------------------------------|---------------------------|
| No user input.                                                                                                                                                                                                                                     |                |          |         |                                                                       |                           |
| 664                                                                                                                                                                                                                                                | Neutral        | 12       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.cytomat.dll | .../cytomatdevice.cs 1746 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                       |                           |
| ↳ Approve Mitigation (Danaher Corporation): Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                |          |         |                                                                       |                           |
| 487                                                                                                                                                                                                                                                | Neutral        | 12       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.cytomat.dll | .../cytomatdevice.cs 1782 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                       |                           |
| ↳ Approve Mitigation (Danaher Corporation): Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                |          |         |                                                                       |                           |
| 488                                                                                                                                                                                                                                                | Neutral        | 12       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.cytomat.dll | .../cytomatdevice.cs 1883 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                       |                           |
| ↳ Potential False Positive (Danaher Corporation): logFullPath is created programmatically... No direct user input used here                                                                                                                        |                |          |         |                                                                       |                           |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                       |                           |
| 489                                                                                                                                                                                                                                                | Neutral        | 12       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.cytomat.dll | .../cytomatdevice.cs 1892 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                       |                           |
| ↳ Potential False Positive (Danaher Corporation): oldLogFileInfos is not a direct user input. File names are returned using an user defined pattern in the code.                                                                                   |                |          |         |                                                                       |                           |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                       |                           |
| 490                                                                                                                                                                                                                                                | Neutral        | 12       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.cytomat.dll | .../cytomatdevice.cs 1905 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has                                  |                |          |         |                                                                       |                           |

| Flaw Id                                                                                                                                                                                                                                                                                                                                  | Exploitability | Module # | Class # | Module                                          | Location                              |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------|---------------------------------------|
| direct file access via Explorer.                                                                                                                                                                                                                                                                                                         |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Filename is constructed programmatically in above lines of code using timestamps appended. This is not user supplied.                                                                                                                                                           |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                               |                |          |         |                                                 |                                       |
| 639                                                                                                                                                                                                                                                                                                                                      | Neutral        | 13       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll       | .../multijobs/delayplatejob.cs 348    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> pointsDirectory- a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                               |                |          |         |                                                 |                                       |
| 640                                                                                                                                                                                                                                                                                                                                      | Neutral        | 13       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll       | .../multijobs/delayplatejob.cs 428    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.                  |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                               |                |          |         |                                                 |                                       |
| 563                                                                                                                                                                                                                                                                                                                                      | Neutral        | 15       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br> | .../device/devicereportmanager.cs 387 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                                                                                                              |                |          |         |                                                 |                                       |
| 564                                                                                                                                                                                                                                                                                                                                      | Neutral        | 15       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br> | .../device/devicereportmanager.cs 584 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fileName is built dynamically, and extension is hardcoded. fileName is not a complete user input.                                                                                                                                                                               |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                              |                |          |         |                                                 |                                       |
| 565                                                                                                                                                                                                                                                                                                                                      | Neutral        | 15       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl     | .../device/devicereportmanager.cs 586 |

| Flaw Id                                                                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                                | Location                              |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------|---------------------------------------|
|                                                                                                                                                                                                                                                    |                |          |         |                                                       |                                       |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                       |                                       |
| ↳ Potential False Positive (Danaher Corporation): fileName is built dynamically, and extension is hardcoded. fileName is not a complete user input.                                                                                                |                |          |         |                                                       |                                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                       |                                       |
| 573                                                                                                                                                                                                                                                | Neutral        | 16       | -       | tempo.exe#3.2.0.1/labcyte.pod.devices.samplestore.dll | .../rundef/devicerundefinition.cs 300 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                       |                                       |
| ↳ Potential False Positive (Danaher Corporation): fullPath is originated from the variable fileFullPath.. This is dynamically constructed at run time and a part of it is hardcoded. There is not direct user input is used here.                  |                |          |         |                                                       |                                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                       |                                       |
| 959                                                                                                                                                                                                                                                | Neutral        | 17       | -       | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll       | .../devicetestform.cs 2397            |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                       |                                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                       |                                       |
| 958                                                                                                                                                                                                                                                | Neutral        | 17       | -       | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll       | .../devicetestform.cs 2585            |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                       |                                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                       |                                       |
| 957                                                                                                                                                                                                                                                | Neutral        | 17       | -       | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll       | .../devicetestform.cs 2590            |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                       |                                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                       |                                       |
| 661                                                                                                                                                                                                                                                | Neutral        | 18       | -       | tempo.exe#3.2.0.1/labcyte.pod.devices.echo.dll        | buildagent/.../echodevice.cs 1126     |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                              | Exploitability | Module # | Class # | Module                                                 | Location                          |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------------------------------------|-----------------------------------|
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                        |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved for different reason. Paths are program defined. No user input.                                                                                                                                                                                                                                                 |                |          |         |                                                        |                                   |
| ✓ 662                                                                                                                                                                                                                                                                                                                                                                | Neutral        | 18       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.ec<br>ho.dll | buildagent/.../echodevice.cs 1188 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                        |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved for different reason. Paths are program defined. No user input.                                                                                                                                                                                                                                                 |                |          |         |                                                        |                                   |
| ✓ 552                                                                                                                                                                                                                                                                                                                                                                | Neutral        | 19       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br>        | .../filestreamlock.cs 74          |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                        |                                   |
| ↳ Potential False Positive (Danaher Corporation): filepath is originated from DB and is not a direct user input from UI                                                                                                                                                                                                                                              |                |          |         |                                                        |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                        |                                   |
| ✓ 926                                                                                                                                                                                                                                                                                                                                                                | Neutral        | 20       | -       | cytomatclient.exe#3.<br>2.0.1                          | buildagent/.../frmdevcaps.cs 101  |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on the server is not applicable here, the opened file exists on the same machine.                                                                                                                                                                                                     |                |          |         |                                                        |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                        |                                   |
| ✓ 924                                                                                                                                                                                                                                                                                                                                                                | Neutral        | 20       | -       | cytomatclient.exe#3.<br>2.0.1                          | buildagent/.../frmdevcaps.cs 219  |
| ↳ Potential False Positive (Danaher Corporation): The path comes from a local SaveFileDialog in a desktop UI flow, not from remote or server-side input. The selected filename is then used only to write the XML document via XmlWriter.Create in a user-initiated Save As action. Unauthorized access to files on server is not applicable here                    |                |          |         |                                                        |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                        |                                   |
| ✓ 913                                                                                                                                                                                                                                                                                                                                                                | Neutral        | 21       | -       | microserveclient.exe<br>#3.2.0.1                       | buildagent/.../frmmain.cs 516     |
| ↳ Potential False Positive (Danaher Corporation): The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path. |                |          |         |                                                        |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                        |                                   |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                                                                                                      | Exploitability | Module # | Class # | Module                                                          | Location                           |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-----------------------------------------------------------------|------------------------------------|
| 928                                                                                                                                                                                                                                                                                                                                                                                                                                          | Neutral        | 22       | -       | ambistoreclient.exe#3.2.0.1                                     | buildagent/.../frmmain.cs 559      |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path; it only writes the log content to the file selected by the user. |                |          |         |                                                                 |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                                                                                                   |                |          |         |                                                                 |                                    |
| 927                                                                                                                                                                                                                                                                                                                                                                                                                                          | Neutral        | 23       | -       | cytomatclient.exe#3.2.0.1                                       | buildagent/.../frmmain.cs 841      |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path.                                                                  |                |          |         |                                                                 |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                                                                                                   |                |          |         |                                                                 |                                    |
| 434                                                                                                                                                                                                                                                                                                                                                                                                                                          | Neutral        | 23       | -       | cytomatclient.exe#3.2.0.1                                       | buildagent/.../frmmain.cs 1198     |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                                                                                    |                |          |         |                                                                 |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> lblFileName.Text is hardcoded in line number 142 and not user supplied.                                                                                                                                                                                                                                                                                                             |                |          |         |                                                                 |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                                                                                                                                   |                |          |         |                                                                 |                                    |
| 813                                                                                                                                                                                                                                                                                                                                                                                                                                          | Neutral        | 24       | -       | tempoadmin.exe#3.2.0.1                                          | buildagent/.../genrobotxcfg.cs 102 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                                                                                    |                |          |         |                                                                 |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> backup is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                                                                                                                                                                                                    |                |          |         |                                                                 |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                 |                                    |
| 812                                                                                                                                                                                                                                                                                                                                                                                                                                          | Neutral        | 24       | -       | tempoadmin.exe#3.2.0.1                                          | buildagent/.../genrobotxcfg.cs 102 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                                                                                    |                |          |         |                                                                 |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> filePath is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                 |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                 |                                    |
| 647                                                                                                                                                                                                                                                                                                                                                                                                                                          | Neutral        | 26       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhec | .../inhecopacdevice.cs 398         |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                                       | Location                    |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------------------------------------------------------------|-----------------------------|
| ocpac.dll                                                                                                                                                                                                                                                 |                |          |         |                                                                              |                             |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                              |                             |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                |          |         |                                                                              |                             |
| 648                                                                                                                                                                                                                                                       | Neutral        | 26       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhec<br>ocpac.dll | .../inhecocpacdevice.cs 412 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                              |                             |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                |          |         |                                                                              |                             |
| 649                                                                                                                                                                                                                                                       | Neutral        | 26       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhec<br>ocpac.dll | .../inhecocpacdevice.cs 449 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                              |                             |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                |          |         |                                                                              |                             |
| 632                                                                                                                                                                                                                                                       | Likely         | 27       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhec<br>oodtc.dll | .../inhecomethodfile.cs 31  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                              |                             |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> _methodFilePath is derived from the value fetched by built in file browser.                                                                                                                      |                |          |         |                                                                              |                             |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented                                                                                                                                                                                 |                |          |         |                                                                              |                             |
| 633                                                                                                                                                                                                                                                       | Neutral        | 27       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhec<br>oodtc.dll | .../inhecomethodfile.cs 64  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                              |                             |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> _methodFilePath is derived from the value fetched by built in file browser.                                                                                                                      |                |          |         |                                                                              |                             |

| Flaw Id                                                                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                                                   | Location                    |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------------------------------------------------------|-----------------------------|
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented                                                                                                                                                                                 |                |          |         |                                                                          |                             |
| 634                                                                                                                                                                                                                                                | Neutral        | 27       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhecoodtc.dll | .../inhecomethodfile.cs 68  |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                          |                             |
| ↳ Potential False Positive (Danaher Corporation): _methodFilePath is derived from the value fetched by built in file browser.                                                                                                                      |                |          |         |                                                                          |                             |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented                                                                                                                                                                                 |                |          |         |                                                                          |                             |
| 643                                                                                                                                                                                                                                                | Neutral        | 27       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhecoodtc.dll | .../inhecomethodfile.cs 127 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                          |                             |
| ↳ Potential False Positive (Danaher Corporation): methodfilePath uses using the built-in file browser.                                                                                                                                             |                |          |         |                                                                          |                             |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented                                                                                                                                                                                 |                |          |         |                                                                          |                             |
| 644                                                                                                                                                                                                                                                | Likely         | 27       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhecoodtc.dll | .../inhecomethodfile.cs 162 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                          |                             |
| ↳ Potential False Positive (Danaher Corporation): _methodFilePath uses the built-in file browser.                                                                                                                                                  |                |          |         |                                                                          |                             |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented                                                                                                                                                                                 |                |          |         |                                                                          |                             |
| 645                                                                                                                                                                                                                                                | Neutral        | 27       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhecoodtc.dll | .../inhecomethodfile.cs 245 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                          |                             |
| ↳ Potential False Positive (Danaher Corporation): _methodFilePath uses the built-in file browser.                                                                                                                                                  |                |          |         |                                                                          |                             |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented                                                                                                                                                                                 |                |          |         |                                                                          |                             |
| 646                                                                                                                                                                                                                                                | Neutral        | 27       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.inhecoodtc.dll | .../inhecomethodfile.cs 317 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has                                  |                |          |         |                                                                          |                             |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                                               | Location                      |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------------------------------------------------------------------|-------------------------------|
| direct file access via Explorer.                                                                                                                                                                                                                          |                |          |         |                                                                                      |                               |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> _methodFilePath uses the built-in file browser.                                                                                                                                                  |                |          |         |                                                                                      |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented                                                                                                                                                                                 |                |          |         |                                                                                      |                               |
| 631                                                                                                                                                                                                                                                       | Neutral        | 28       | -       | tempo.exe#3.2.0.1/lab<br>byte.pod.api.dll/lab<br>byte.pod.devices.inhec<br>oodtc.dll | .../inhecoodtcdevice.cs 834   |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                                      |                               |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> summaryFilePath is derived dynamically at runtime and is not a direct user input from UI.                                                                                                        |                |          |         |                                                                                      |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                      |                               |
| 616                                                                                                                                                                                                                                                       | Neutral        | 29       | -       | tempo.exe#3.2.0.1/lab<br>byte.pod.devices.int<br>ellixcap.dll                        | .../intellixcapdevice.cs 379  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                                      |                               |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> dirPath is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                |                |          |         |                                                                                      |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                      |                               |
| 628                                                                                                                                                                                                                                                       | Neutral        | 29       | -       | tempo.exe#3.2.0.1/lab<br>byte.pod.devices.int<br>ellixcap.dll                        | .../intellixcapdevice.cs 1313 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                                      |                               |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> _deviceLogPath is not a direct user input instaed its created dynamically at runtime.                                                                                                            |                |          |         |                                                                                      |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                      |                               |
| 629                                                                                                                                                                                                                                                       | Neutral        | 29       | -       | tempo.exe#3.2.0.1/lab<br>byte.pod.devices.int<br>ellixcap.dll                        | .../intellixcapdevice.cs 1327 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                                      |                               |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> _deviceLogFilename is derived dynamically at runtime and doesn't use direct user input value here                                                                                                |                |          |         |                                                                                      |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                      |                               |
| 630                                                                                                                                                                                                                                                       | Neutral        | 29       | -       | tempo.exe#3.2.0.1/lab<br>byte.pod.devices.int                                        | .../intellixcapdevice.cs 1366 |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                          | Location                              |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------|---------------------------------------|
| ellixcap.dll                                                                                                                                                                                                                                              |                |          |         |                                                 |                                       |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> _deviceLogFilename is derived dynamically at runtime and doesn't use direct user input value here                                                                                                |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                 |                                       |
| 550                                                                                                                                                                                                                                                       | Neutral        | 31       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br> | buildagent/.../job/jobmanager.cs 1670 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> filename is program generated and no direct user input used here.                                                                                                                                |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                 |                                       |
| 548                                                                                                                                                                                                                                                       | Neutral        | 32       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br> | .../job/jobmanagerbase.cs 743         |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> directory name is hardcoded here and the filename is program generated. No direct user input from UI is used here                                                                                |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                 |                                       |
| 549                                                                                                                                                                                                                                                       | Neutral        | 32       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br> | .../job/jobmanagerbase.cs 746         |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> directory name is hardcoded here and the filename is originated from DB. No direct user input from UI is used here                                                                               |                |          |         |                                                 |                                       |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                 |                                       |
| 531                                                                                                                                                                                                                                                       | Neutral        | 33       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br> | .../job/jobpropertymanager.cs 66      |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                 |                                       |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> dir is programmatically created in code, and no user                                                                                                                                             |                |          |         |                                                 |                                       |

| Flaw Id                                                                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                                                    | Location                          |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------------------------|-----------------------------------|
| input                                                                                                                                                                                                                                              |                |          |         |                                                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                           |                                   |
| 532                                                                                                                                                                                                                                                | Neutral        | 33       | -       | tempo.exe#3.2.0.1/labcyte.pod.common.dll                                  | .../job/jobpropertymanager.cs 70  |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                   |
| ↳ Potential False Positive (Danaher Corporation): fileName is programmatically created in code, and no user input                                                                                                                                  |                |          |         |                                                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                           |                                   |
| 533                                                                                                                                                                                                                                                | Neutral        | 33       | -       | tempo.exe#3.2.0.1/labcyte.pod.common.dll                                  | .../job/jobpropertymanager.cs 109 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                   |
| ↳ Potential False Positive (Danaher Corporation): dir is programmatically created in code, and no user input                                                                                                                                       |                |          |         |                                                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                           |                                   |
| 546                                                                                                                                                                                                                                                | Neutral        | 33       | -       | tempo.exe#3.2.0.1/labcyte.pod.common.dll                                  | .../job/jobpropertymanager.cs 122 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                   |
| ↳ Potential False Positive (Danaher Corporation): filename is not a direct user input, its constructed dynamically at runtime using time stamps                                                                                                    |                |          |         |                                                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                                           |                                   |
| 547                                                                                                                                                                                                                                                | Neutral        | 33       | -       | tempo.exe#3.2.0.1/labcyte.pod.common.dll                                  | .../job/jobpropertymanager.cs 150 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                   |
| ↳ Potential False Positive (Danaher Corporation): filename is not a direct user input, its constructed dynamically at runtime using time stamps                                                                                                    |                |          |         |                                                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                                           |                                   |
| 557                                                                                                                                                                                                                                                | Likely         | 34       | -       | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.vcode.dll | .../labelplateusingfilejob.cs 124 |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                                    | Location                          |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------------------------|-----------------------------------|
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> This use built in file browser to browse the Barcode file with the hardcoded extension                                                                                                           |                |          |         |                                                                           |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                |                |          |         |                                                                           |                                   |
| 558                                                                                                                                                                                                                                                       | Likely         | 34       | -       | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.vcode.dll | .../labelplateusingfilejob.cs 212 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Part of index file is hardcoded, this file doesn't contain any sensitive data, instead it contains only fields that are used to print labels on a microplate.                                    |                |          |         |                                                                           |                                   |
| ↳ <i>Mitigate by Design (Danaher Corporation):</i> These are Mitigated by OS Environment                                                                                                                                                                  |                |          |         |                                                                           |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                           |                                   |
| 559                                                                                                                                                                                                                                                       | Likely         | 34       | -       | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.vcode.dll | .../labelplateusingfilejob.cs 239 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Part of index file is hardcoded, this file doesn't contain any sensitive data, instead it contains only fields that are used to print labels on a microplate.                                    |                |          |         |                                                                           |                                   |
| ↳ <i>Mitigate by Design (Danaher Corporation):</i> These are Mitigated by OS Environment                                                                                                                                                                  |                |          |         |                                                                           |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                           |                                   |
| 560                                                                                                                                                                                                                                                       | Neutral        | 34       | -       | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.vcode.dll | .../labelplateusingfilejob.cs 267 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                           |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Part of index file is hardcoded, this file doesn't contain any sensitive data, instead it contains only fields that are used to print labels on a microplate.                                    |                |          |         |                                                                           |                                   |
| ↳ <i>Mitigate by Design (Danaher Corporation):</i> These are Mitigated by OS Environment                                                                                                                                                                  |                |          |         |                                                                           |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                           |                                   |
| 561                                                                                                                                                                                                                                                       | Likely         | 34       | -       | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.vcode.dll | .../labelplateusingfilejob.cs 292 |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                                | Location                          |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-----------------------------------------------------------------------|-----------------------------------|
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                       |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Barcode file doesn't contain any sensitive data, instead it contains only fields that are used to print labels on a microplate.                                                                  |                |          |         |                                                                       |                                   |
| ↳ <i>Mitigate by Design (Danaher Corporation):</i> These are Mitigated by OS Environment                                                                                                                                                                  |                |          |         |                                                                       |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                       |                                   |
|                                                                                                                                                                                                                                                           | 933            | Neutral  | 39      | - tempo.exe#3.2.0.1/labbyte.pod.api.dll/labbyte.pod.devices.vcode.dll | .../lookuptableform.cs 221        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                       |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                |                |          |         |                                                                       |                                   |
|                                                                                                                                                                                                                                                           | 932            | Neutral  | 39      | - tempo.exe#3.2.0.1/labbyte.pod.api.dll/labbyte.pod.devices.vcode.dll | .../lookuptableform.cs 280        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                       |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                |                |          |         |                                                                       |                                   |
|                                                                                                                                                                                                                                                           | 672            | Neutral  | 40      | - tempotester.exe#3.2.0.1/labbyte.common.mailer.dll                   | buildagent/.../mailservice.cs 236 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                       |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The parameter to this method is constructed dynamically at run time, and doesn't directly uses any input data from UI.                                                                           |                |          |         |                                                                       |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                       |                                   |
|                                                                                                                                                                                                                                                           | 670            | Neutral  | 40      | - tempo.exe#3.2.0.1/labbyte.common.mailer.dll                         | buildagent/.../mailservice.cs 505 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                       |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> str is constructed dynamically at runtime and doesn't use any direct user input.                                                                                                                 |                |          |         |                                                                       |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                       |                                   |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                           | Exploitability | Module # | Class # | Module                                      | Location                       |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------|--------------------------------|
| 923                                                                                                                                                                                                                                                                                                                                                               | Neutral        | 42       | -       | inhecoclient.exe#3.2.0.1                    | buildagent/.../mainform.cs 260 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                    |                |          |         |                                             |                                |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                        |                |          |         |                                             |                                |
| 922                                                                                                                                                                                                                                                                                                                                                               | Neutral        | 42       | -       | inhecoclient.exe#3.2.0.1                    | buildagent/.../mainform.cs 282 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                    |                |          |         |                                             |                                |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                        |                |          |         |                                             |                                |
| 921                                                                                                                                                                                                                                                                                                                                                               | Neutral        | 42       | -       | inhecoclient.exe#3.2.0.1                    | buildagent/.../mainform.cs 313 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                    |                |          |         |                                             |                                |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                        |                |          |         |                                             |                                |
| 916                                                                                                                                                                                                                                                                                                                                                               | Neutral        | 43       | -       | certusclient.exe#3.2.0.1                    | buildagent/.../mainform.cs 406 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The flagged ReadAllText call is part of a local Windows desktop workflow that loads a user-selected XML file through OpenFileDialog. The application does not construct the path from remote or untrusted server input, and there is no path traversal or arbitrary file access logic in this code path. |                |          |         |                                             |                                |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                        |                |          |         |                                             |                                |
| 915                                                                                                                                                                                                                                                                                                                                                               | Neutral        | 43       | -       | certusclient.exe#3.2.0.1                    | buildagent/.../mainform.cs 481 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The flagged ReadAllText call is part of a local Windows desktop workflow that loads a user-selected XML file through OpenFileDialog. The application does not construct the path from remote or untrusted server input, and there is no path traversal or arbitrary file access logic in this code path. |                |          |         |                                             |                                |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                        |                |          |         |                                             |                                |
| 914                                                                                                                                                                                                                                                                                                                                                               | Neutral        | 43       | -       | certusclient.exe#3.2.0.1                    | buildagent/.../mainform.cs 512 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The flagged ReadAllText call is part of a local Windows desktop workflow that loads a user-selected XML file through OpenFileDialog. The application does not construct the path from remote or untrusted server input, and there is no path traversal or arbitrary file access logic in this code path. |                |          |         |                                             |                                |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                        |                |          |         |                                             |                                |
| 589                                                                                                                                                                                                                                                                                                                                                               | Neutral        | 45       | -       | tempo.exe#3.2.0.1/la<br>byte.pod.devices.mi | .../microservedevice.cs 2881   |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                       | Location                     |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------------------------------------------|------------------------------|
| croserve.dll                                                                                                                                                                                                                                              |                |          |         |                                                              |                              |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                              |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                |          |         |                                                              |                              |
| 590                                                                                                                                                                                                                                                       | Neutral        | 45       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.mi<br>croserve.dll | .../microservedevice.cs 2895 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                              |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                |          |         |                                                              |                              |
| 591                                                                                                                                                                                                                                                       | Neutral        | 45       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.mi<br>croserve.dll | .../microservedevice.cs 2931 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                              |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved for different reason. Paths are program defined. No user input.                                                                                                                               |                |          |         |                                                              |                              |
| 1011                                                                                                                                                                                                                                                      | Neutral        | 49       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br>              | .../plateimportexport.cs 61  |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                              |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                |                |          |         |                                                              |                              |
| 529                                                                                                                                                                                                                                                       | Neutral        | 49       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br>              | .../plateimportexport.cs 325 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                              |                              |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> FileName extension is hardcoded and the value is taken from windows inbuilt SaveFile Dialog                                                                                                      |                |          |         |                                                              |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                |                |          |         |                                                              |                              |
| 530                                                                                                                                                                                                                                                       | Neutral        | 49       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br>              | .../plateimportexport.cs 329 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application.                                                                                                                                          |                |          |         |                                                              |                              |

| Flaw Id                                                                                                                                                                                                                                                                                          | Exploitability | Module # | Class # | Module                                                        | Location                        |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------------|---------------------------------|
| All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                                                         |                |          |         |                                                               |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> FileName extension is hardcoded and the value is taken from windows inbuilt SaveFile Dialog                                                                                                                                             |                |          |         |                                                               |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                       |                |          |         |                                                               |                                 |
| 525                                                                                                                                                                                                                                                                                              | Neutral        | 50       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.platestorag<br>econtrol.dll | .../platestoragegrid.cs 8933    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                        |                |          |         |                                                               |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> A part of root variable is hardcoded and the devicename is originated from DB. Devicename is allowed to modify by the Admin user alone while setting up Tempo initially.                                                                |                |          |         |                                                               |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                       |                |          |         |                                                               |                                 |
| 514                                                                                                                                                                                                                                                                                              | Neutral        | 54       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br>               | buildagent/.../podlogger.cs 123 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                        |                |          |         |                                                               |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> logger.FileName is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                               |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                       |                |          |         |                                                               |                                 |
| 515                                                                                                                                                                                                                                                                                              | Neutral        | 54       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br>               | buildagent/.../podlogger.cs 125 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                        |                |          |         |                                                               |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> logger.FileName is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                               |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                       |                |          |         |                                                               |                                 |
| 516                                                                                                                                                                                                                                                                                              | Neutral        | 54       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br>               | buildagent/.../podlogger.cs 194 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                        |                |          |         |                                                               |                                 |

| Flaw Id                                                                                                                                                                                                                                                                                          | Exploitability | Module # | Class # | Module                                       | Location                           |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|----------------------------------------------|------------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> logger.FileName is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                              |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                       |                |          |         |                                              |                                    |
|  517                                                                                                                                                                                                            | Neutral        | 54       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dll | buildagent/.../podlogger.cs 198    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                        |                |          |         |                                              |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> logger.FileName is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                              |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                       |                |          |         |                                              |                                    |
|  527                                                                                                                                                                                                            | Neutral        | 54       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dll | buildagent/.../podlogger.cs 239    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                        |                |          |         |                                              |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Filename is not a direct user input.                                                                                                                                                                                                    |                |          |         |                                              |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented                                                                                                                                                                                                                        |                |          |         |                                              |                                    |
|  528                                                                                                                                                                                                          | Neutral        | 54       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dll | buildagent/.../podlogger.cs 275    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                        |                |          |         |                                              |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fileName is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.        |                |          |         |                                              |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                       |                |          |         |                                              |                                    |
|  503                                                                                                                                                                                                          | Neutral        | 55       | -       | tempo.exe#3.2.0.1                            | .../symphony/podoptionsform.cs 593 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                        |                |          |         |                                              |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> oldLog is fetched from DB and not a direct user input                                                                                                                                                                                   |                |          |         |                                              |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                       |                |          |         |                                              |                                    |
|  504                                                                                                                                                                                                          | Neutral        | 55       | -       | tempo.exe#3.2.0.1                            | .../symphony/podoptionsform.cs 601 |

| Flaw Id                                                                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                    | Location                           |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------|------------------------------------|
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                           |                                    |
| ↳ Potential False Positive (Danaher Corporation): oldLog is fetched from DB and not a direct user input                                                                                                                                            |                |          |         |                                           |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                           |                                    |
| 1007                                                                                                                                                                                                                                               | Neutral        | 56       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | buildagent/.../podprotocol.cs 2756 |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                           |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                           |                                    |
| 1006                                                                                                                                                                                                                                               | Neutral        | 56       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | buildagent/.../podprotocol.cs 2889 |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                           |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                |                |          |         |                                           |                                    |
| 635                                                                                                                                                                                                                                                | Neutral        | 57       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | buildagent/.../podsystem.cs 152    |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                           |                                    |
| ↳ Potential False Positive (Danaher Corporation): backupFile is constructed dynamically at runtime, logPrefs.FileName is hardcoded. There is not direct user input used here.                                                                      |                |          |         |                                           |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                           |                                    |
| 636                                                                                                                                                                                                                                                | Neutral        | 57       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | buildagent/.../podsystem.cs 152    |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                           |                                    |
| ↳ Potential False Positive (Danaher Corporation): backupFile is constructed dynamically at runtime, logPrefs.FileName is hardcoded. There is not direct user input used here.                                                                      |                |          |         |                                           |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                           |                                    |
| 637                                                                                                                                                                                                                                                | Neutral        | 57       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | buildagent/.../podsystem.cs 161    |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                           |                                    |
| ↳ Potential False Positive (Danaher Corporation): Both errorLog, backupErrorFile are constructed                                                                                                                                                   |                |          |         |                                           |                                    |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                                     | Exploitability | Module # | Class # | Module                                    | Location                           |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------|------------------------------------|
| dynamically at runtime and no direct user input are used here.                                                                                                                                                                                                                                                                                                              |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                                                                  |                |          |         |                                           |                                    |
| 638                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | 57       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | buildagent/.../podsystem.cs 161    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                           |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Both errorLog, backupErrorFile are constructed dynamically at runtime and no direct user input are used here.                                                                                                                                                                                                      |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                                                                  |                |          |         |                                           |                                    |
| 917                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | 58       | -       | robotteachingutility.e<br>xe#3.2.0.1      | .../preciseteachingform.cs 2458    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path. |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                           |                                    |
| 502                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | 59       | -       | tempo.exe#3.2.0.1                         | buildagent/.../program.app.cs 201  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                           |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> LogDirectory is not a value entered from UI                                                                                                                                                                                                                                                                        |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                                                                  |                |          |         |                                           |                                    |
| 673                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | 60       | -       | tempotester.exe#3.2.<br>0.1               | buildagent/.../program.cs 177      |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                           |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> pref.LogDirectory is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.                                                                          |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                                                                  |                |          |         |                                           |                                    |
| 1012                                                                                                                                                                                                                                                                                                                                                                        | Neutral        | 61       | -       | tempo.exe#3.2.0.1                         | buildagent/.../protocolview.cs 473 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                           |                                    |
| 668                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | 62       | -       | tempo.exe#3.2.0.1/la                      | .../jobs/querytimepointjob.cs 243  |

| Flaw Id                                                                                                                                                                                                                                                                                                                       | Exploitability | Module # | Class # | Module                                                                                     | Location                          |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------------------------------------------------------------------------|-----------------------------------|
| bctype.pod.api.dll                                                                                                                                                                                                                                                                                                            |                |          |         |                                                                                            |                                   |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                     |                |          |         |                                                                                            |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> pointsDirectory is program created and a part of it uses saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.         |                |          |         |                                                                                            |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                    |                |          |         |                                                                                            |                                   |
| 669                                                                                                                                                                                                                                                                                                                           | Neutral        | 62       | -       | tempo.exe#3.2.0.1/la<br>bctype.pod.api.dll                                                 | .../jobs/querytimepointjob.cs 326 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                     |                |          |         |                                                                                            |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fullFileName is program created and a part of it uses saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.            |                |          |         |                                                                                            |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                   |                |          |         |                                                                                            |                                   |
| 807                                                                                                                                                                                                                                                                                                                           | Neutral        | 63       | -       | tempo.exe#3.2.0.1/la<br>bctype.pod.api.dll/becl<br>s.pod.devices.volume<br>check.dll       | .../jobs/readvolumesjob.cs 102    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                     |                |          |         |                                                                                            |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fileName and path are constructed dynamically at run time, and it doesn't use the user input directly here.                                                                                                                                                          |                |          |         |                                                                                            |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                   |                |          |         |                                                                                            |                                   |
| 497                                                                                                                                                                                                                                                                                                                           | Neutral        | 67       | -       | tempo.exe#3.2.0.1/la<br>bctype.pod.api.dll/lab<br>cbyte.pod.devices.asm.<br>dll            | .../remoteconnector.cs 315        |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                     |                |          |         |                                                                                            |                                   |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> useDirectory is either hardcoded in the file or fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                                                            |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                    |                |          |         |                                                                                            |                                   |
| 658                                                                                                                                                                                                                                                                                                                           | Neutral        | 65       | -       | tempo.exe#3.2.0.1/la<br>bctype.pod.api.dll/lab<br>cbyte.pod.devices.envi<br>sionreader.dll | .../remoteconnector.cs 325        |

| Flaw Id                                                                                                                                                                                                                                                                                                                | Exploitability | Module # | Class # | Module                                                                  | Location                   |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------------------------|----------------------------|
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                     |                |          |         |                                                                         |                            |
| ↳ Potential False Positive (Danaher Corporation): useDirectory is either hardcoded in the file or fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                                         |                            |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                                                                    |                |          |         |                                                                         |                            |
| 805                                                                                                                                                                                                                                                                                                                    | Neutral        | 66       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.bmg.dll       | .../remoteconnector.cs 331 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                     |                |          |         |                                                                         |                            |
| ↳ Potential False Positive (Danaher Corporation): useDirectory is program created and a part of it uses saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.            |                |          |         |                                                                         |                            |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                                                                    |                |          |         |                                                                         |                            |
| 510                                                                                                                                                                                                                                                                                                                    | Neutral        | 67       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.asm.dll       | .../remoteconnector.cs 355 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                     |                |          |         |                                                                         |                            |
| ↳ Potential False Positive (Danaher Corporation): filePath is constructed dynamically at runtime using timestamp and also the name is not a direct userinput.                                                                                                                                                          |                |          |         |                                                                         |                            |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                                                                    |                |          |         |                                                                         |                            |
| 496                                                                                                                                                                                                                                                                                                                    | Neutral        | 66       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.bmg.dll       | .../remoteconnector.cs 366 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                     |                |          |         |                                                                         |                            |
| ↳ Potential False Positive (Danaher Corporation): filePath is constructed dynamically with appended timespamps during runtime. This is not a direct user input                                                                                                                                                         |                |          |         |                                                                         |                            |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                                                                   |                |          |         |                                                                         |                            |
| 332                                                                                                                                                                                                                                                                                                                    | Neutral        | 65       | -       | envisionreaderclient.exe#9.9.9.0/labcyte.pod.devices.envisionreader.dll | .../remoteconnector.cs 371 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application.                                                                                                                                                                                                              |                |          |         |                                                                         |                            |

| Flaw Id                                                                                                                                                                                                                                                                                                               | Exploitability | Module # | Class # | Module                                                              | Location                           |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------------------|------------------------------------|
| All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                                                                              |                |          |         |                                                                     |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Here filePath is not a direct user input, instead its created dynamically at runtime using timestamp.                                                                                                                                                        |                |          |         |                                                                     |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                            |                |          |         |                                                                     |                                    |
| 666                                                                                                                                                                                                                                                                                                                   | Neutral        | 68       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll                               | .../jobs/removetimepointjob.cs 157 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                             |                |          |         |                                                                     |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> pointsDirectory is program created and a part of it uses saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                                     |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                           |                |          |         |                                                                     |                                    |
| 667                                                                                                                                                                                                                                                                                                                   | Neutral        | 68       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll                               | .../jobs/removetimepointjob.cs 222 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                             |                |          |         |                                                                     |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fullfileName is program created and a part of it uses saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.    |                |          |         |                                                                     |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                            |                |          |         |                                                                     |                                    |
| 540                                                                                                                                                                                                                                                                                                                   | Neutral        | 69       | -       | robotteachingutility.exe#3.2.0.1/labcyte.pod.driver.commonrobot.dll | buildagent/.../robotconfig.cs 443  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                             |                |          |         |                                                                     |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> filename is not a direct user input, its constructed dynamically at runtime using time stamps                                                                                                                                                                |                |          |         |                                                                     |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                           |                |          |         |                                                                     |                                    |
| 541                                                                                                                                                                                                                                                                                                                   | Neutral        | 69       | -       | robotteachingutility.exe#3.2.0.1/labcyte.pod.driver.commonrobot.dll | buildagent/.../robotconfig.cs 443  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                             |                |          |         |                                                                     |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> filename is not a direct user input, its constructed                                                                                                                                                                                                         |                |          |         |                                                                     |                                    |

| Flaw Id                                                                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                                              | Location                             |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------------------|--------------------------------------|
| dynamically at runtime using time stamps                                                                                                                                                                                                           |                |          |         |                                                                     |                                      |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                                     |                                      |
| 678                                                                                                                                                                                                                                                | Neutral        | 69       | -       | robotteachingutility.exe#3.2.0.1/labcyte.pod.driver.commonrobot.dll | buildagent/.../robotconfig.cs 594    |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                     |                                      |
| ↳ Potential False Positive (Danaher Corporation): fName is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                  |                |          |         |                                                                     |                                      |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                                     |                                      |
| 679                                                                                                                                                                                                                                                | Neutral        | 69       | -       | robotteachingutility.exe#3.2.0.1/labcyte.pod.driver.commonrobot.dll | buildagent/.../robotconfig.cs 637    |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                     |                                      |
| ↳ Potential False Positive (Danaher Corporation): fi is derived from fname and fname is derived dynamically at runtime and is not a direct user input from UI.                                                                                     |                |          |         |                                                                     |                                      |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                     |                                      |
| 680                                                                                                                                                                                                                                                | Neutral        | 69       | -       | robotteachingutility.exe#3.2.0.1/labcyte.pod.driver.commonrobot.dll | buildagent/.../robotconfig.cs 649    |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                     |                                      |
| ↳ Potential False Positive (Danaher Corporation): fname is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                  |                |          |         |                                                                     |                                      |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                                                     |                                      |
| 814                                                                                                                                                                                                                                                | Neutral        | 70       | -       | tempoadmin.exe#3.2.0.1                                              | .../autogen/robotconfigloader.cs 68  |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                     |                                      |
| ↳ Potential False Positive (Danaher Corporation): A part of filepath is hardcoded and uses in built folder browser dialog to fetch the folder.                                                                                                     |                |          |         |                                                                     |                                      |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                                                     |                                      |
| 804                                                                                                                                                                                                                                                | Neutral        | 70       | -       | tempoadmin.exe#3.2.0.1                                              | .../autogen/robotconfigloader.cs 143 |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                              | Exploitability | Module # | Class # | Module                                                  | Location                           |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------|------------------------------------|
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                         |                                    |
| ↳ Potential False Positive (Danaher Corporation): filePath is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                                                                                                                                 |                |          |         |                                                         |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                         |                                    |
| 574                                                                                                                                                                                                                                                                                                                                                                  | Neutral        | 71       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.ro<br>bot.dll | buildagent/.../robotdevice.cs 1652 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                         |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved for different reason. Paths are program defined. No user input.                                                                                                                                                                                                                                                 |                |          |         |                                                         |                                    |
| 587                                                                                                                                                                                                                                                                                                                                                                  | Neutral        | 71       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.ro<br>bot.dll | buildagent/.../robotdevice.cs 1674 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                         |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved for different reason. Paths are program defined. No user input.                                                                                                                                                                                                                                                 |                |          |         |                                                         |                                    |
| 588                                                                                                                                                                                                                                                                                                                                                                  | Neutral        | 71       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.devices.ro<br>bot.dll | buildagent/.../robotdevice.cs 1713 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                         |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approved for different reason. Paths are program defined. No user input.                                                                                                                                                                                                                                                 |                |          |         |                                                         |                                    |
| 918                                                                                                                                                                                                                                                                                                                                                                  | Neutral        | 72       | -       | robotteachingutility.e<br>xe#3.2.0.1                    | .../rorzeteachingfrm.cs 1674       |
| ↳ Potential False Positive (Danaher Corporation): The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path. |                |          |         |                                                         |                                    |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                         |                                    |
| 1009                                                                                                                                                                                                                                                                                                                                                                 | Neutral        | 74       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll               | .../protocol/rundefinition.cs 523  |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server                                                                                                                                                  |                |          |         |                                                         |                                    |

| Flaw Id                                                                                                                                                                                                                                                                         | Exploitability | Module # | Class # | Module                                    | Location                          |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------|-----------------------------------|
| path involved here.                                                                                                                                                                                                                                                             |                |          |         |                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                             |                |          |         |                                           |                                   |
| 1008                                                                                                                                                                                                                                                                            | Neutral        | 74       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../protocol/rundefinition.cs 679 |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                         |                |          |         |                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                             |                |          |         |                                           |                                   |
| 623                                                                                                                                                                                                                                                                             | Neutral        | 75       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutormisc.cs 157        |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                              |                |          |         |                                           |                                   |
| ↳ Potential False Positive (Danaher Corporation): fname is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                             |                |          |         |                                           |                                   |
| 599                                                                                                                                                                                                                                                                             | Neutral        | 76       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutorreports.cs 111     |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                              |                |          |         |                                           |                                   |
| ↳ Potential False Positive (Danaher Corporation): ActionLogFilename is derived dynamcially and is not a direct user input from UI                                                                                                                                               |                |          |         |                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                             |                |          |         |                                           |                                   |
| 600                                                                                                                                                                                                                                                                             | Neutral        | 76       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutorreports.cs 222     |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                              |                |          |         |                                           |                                   |
| ↳ Potential False Positive (Danaher Corporation): path is derived dynamically and is not a direct user input from UI                                                                                                                                                            |                |          |         |                                           |                                   |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                            |                |          |         |                                           |                                   |
| 601                                                                                                                                                                                                                                                                             | Neutral        | 76       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutorreports.cs 242     |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                              |                |          |         |                                           |                                   |
| ↳ Potential False Positive (Danaher Corporation): path is derived dynamically and is not a direct user input from UI                                                                                                                                                            |                |          |         |                                           |                                   |

| Flaw Id                                                                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                     | Location                      |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------------------------|-------------------------------|
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                            |                               |
| 618                                                                                                                                                                                                                                                | Neutral        | 76       | -       | tempoadmin.exe#3.2.0.1/labcyte.pod.api.dll | .../runexecutorreports.cs 426 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                            |                               |
| ↳ Potential False Positive (Danaher Corporation): fullFilename is derived dynamically at runtime and is not a direct user input from UI.                                                                                                           |                |          |         |                                            |                               |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                            |                               |
| 619                                                                                                                                                                                                                                                | Neutral        | 76       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll      | .../runexecutorreports.cs 881 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                            |                               |
| ↳ Potential False Positive (Danaher Corporation): fi is originated from path, which is derived dynamically at runtime and is not a direct user input from UI.                                                                                      |                |          |         |                                            |                               |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                            |                               |
| 620                                                                                                                                                                                                                                                | Neutral        | 76       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll      | .../runexecutorreports.cs 886 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                            |                               |
| ↳ Potential False Positive (Danaher Corporation): path, is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                  |                |          |         |                                            |                               |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                               |                |          |         |                                            |                               |
| 621                                                                                                                                                                                                                                                | Neutral        | 76       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll      | .../runexecutorreports.cs 954 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                            |                               |
| ↳ Potential False Positive (Danaher Corporation): fi is originated from path, which is derived dynamically at runtime and is not a direct user input from UI.                                                                                      |                |          |         |                                            |                               |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                |                |          |         |                                            |                               |
| 622                                                                                                                                                                                                                                                | Neutral        | 76       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll      | .../runexecutorreports.cs 959 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                            |                               |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                    | Location                           |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------|------------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> path, is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                  |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                           |                                    |
| 594                                                                                                                                                                                                                                                       | Neutral        | 77       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutorruncreation.cs 621  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                           |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Values are derived dynamically and no direct user input is used here.                                                                                                                            |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                           |                                    |
| 595                                                                                                                                                                                                                                                       | Neutral        | 77       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutorruncreation.cs 626  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                           |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Values are derived dynamically and no direct user input is used here.                                                                                                                            |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                           |                                    |
| 596                                                                                                                                                                                                                                                       | Neutral        | 77       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutorruncreation.cs 725  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                           |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Values are derived dynamically and no direct user input is used here.                                                                                                                            |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                           |                                    |
| 597                                                                                                                                                                                                                                                       | Neutral        | 77       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutorruncreation.cs 3404 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                           |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Values are derived dynamically and no direct user input is used here.                                                                                                                            |                |          |         |                                           |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                           |                                    |
| 598                                                                                                                                                                                                                                                       | Neutral        | 77       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runexecutorruncreation.cs 3415 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has                                  |                |          |         |                                           |                                    |

| Flaw Id                                                                                                                                                                                                                                                          | Exploitability | Module # | Class # | Module                                                                            | Location                     |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-----------------------------------------------------------------------------------|------------------------------|
| direct file access via Explorer.                                                                                                                                                                                                                                 |                |          |         |                                                                                   |                              |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> path is not a direct user input from UI instead derived dynamically using DB value. DB value is RunID which is the database internal key value. This is not supplied by user.                           |                |          |         |                                                                                   |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                      |                |          |         |                                                                                   |                              |
| 954                                                                                                                                                                                                                                                              | Neutral        | 78       | -       | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll/labcyte.pod.storagereprotocol.dll | .../runoptionscontrol.cs 477 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                   |                |          |         |                                                                                   |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                       |                |          |         |                                                                                   |                              |
| 953                                                                                                                                                                                                                                                              | Neutral        | 78       | -       | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll/labcyte.pod.storagereprotocol.dll | .../runoptionscontrol.cs 482 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                   |                |          |         |                                                                                   |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                       |                |          |         |                                                                                   |                              |
| 562                                                                                                                                                                                                                                                              | Neutral        | 79       | -       | tempoadmin.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.tempest.dll        | .../runprotocoljob.cs 361    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.        |                |          |         |                                                                                   |                              |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> ExistingFile is constructed dynamically at run time, using hardcoded search pattern. This is not a direct user input from UI                                                                            |                |          |         |                                                                                   |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                       |                |          |         |                                                                                   |                              |
| 570                                                                                                                                                                                                                                                              | Neutral        | 79       | -       | tempoadmin.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.tempest.dll        | .../runprotocoljob.cs 361    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.        |                |          |         |                                                                                   |                              |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> newFileName is program generated. A part of it is RunID which is fetched from DB. RunID is not passed from UI, instead its primary key autogenerated in DB. There is no direct user input is used here. |                |          |         |                                                                                   |                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                      |                |          |         |                                                                                   |                              |
| 578                                                                                                                                                                                                                                                              | Neutral        | 80       | -       | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll/labcyte.pod.storagereprotocol.dll | .../runreporhandler.cs 26    |

| Flaw Id                                                                                                                                                                                                                                                                                                                         | Exploitability | Module # | Class # | Module                                           | Location                    |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------------------------------|-----------------------------|
| bctype.pod.api.dll                                                                                                                                                                                                                                                                                                              |                |          |         |                                                  |                             |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                       |                |          |         |                                                  |                             |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.                                                                                    |                |          |         |                                                  |                             |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                      |                |          |         |                                                  |                             |
| 579                                                                                                                                                                                                                                                                                                                             | Neutral        | 80       | -       | tempo.exe#3.2.0.1/labctype.pod.api.dll           | .../runreporthandler.cs 94  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                       |                |          |         |                                                  |                             |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> runSummaryPath is either hardcoded in the file or fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                  |                             |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                      |                |          |         |                                                  |                             |
| 580                                                                                                                                                                                                                                                                                                                             | Neutral        | 80       | -       | tempoadmin.exe#3.2.0.1/labctype.pod.api.dll      | .../runreporthandler.cs 135 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                       |                |          |         |                                                  |                             |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> There is no direct user input used here, paths are derived dynamically at runtime, and a part of it is from settings.                                                                                                                                                  |                |          |         |                                                  |                             |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                     |                |          |         |                                                  |                             |
| 581                                                                                                                                                                                                                                                                                                                             | Neutral        | 80       | -       | tempoadmin.exe#3.2.0.1/labctype.pod.api.dll      | .../runreporthandler.cs 135 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                       |                |          |         |                                                  |                             |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> There is no direct user input used here, paths are derived dynamically at runtime, and a part of it is from settings.                                                                                                                                                  |                |          |         |                                                  |                             |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                     |                |          |         |                                                  |                             |
| 509                                                                                                                                                                                                                                                                                                                             | Neutral        | 81       | -       | tempo.exe#3.2.0.1/labctype.pod.userinterface.dll | .../runreportsform.cs 607   |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                       |                |          |         |                                                  |                             |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                  | Location                           |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------|------------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fileName is originated from DB and not a direct user input                                                                                                                                       |                |          |         |                                                         |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                |                |          |         |                                                         |                                    |
| 520                                                                                                                                                                                                                                                       | Neutral        | 81       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.userinterfa<br>ce.dll | .../runreportsform.cs 620          |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                         |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> tempPath- filename is hardcoded, and the pah is program generated. It is not directly user supplied value from UI                                                                                |                |          |         |                                                         |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                |                |          |         |                                                         |                                    |
| 521                                                                                                                                                                                                                                                       | Neutral        | 81       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.userinterfa<br>ce.dll | .../runreportsform.cs 657          |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                         |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> tempPath- filename is hardcoded, and the pah is program generated. It is not directly user supplied value from UI                                                                                |                |          |         |                                                         |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                |                |          |         |                                                         |                                    |
| 522                                                                                                                                                                                                                                                       | Neutral        | 81       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.userinterfa<br>ce.dll | .../runreportsform.cs 662          |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                         |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> tempPath- filename is hardcoded, and the pah is program generated. It is not directly user supplied value from UI                                                                                |                |          |         |                                                         |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                |                |          |         |                                                         |                                    |
| 1010                                                                                                                                                                                                                                                      | Neutral        | 82       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll               | buildagent/.../runscheduler.cs 201 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                         |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                |                |          |         |                                                         |                                    |
| 576                                                                                                                                                                                                                                                       | Neutral        | 83       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll               | .../runscheduler.log.cs 114        |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                         |                                    |

| Flaw Id                                                                                                                                                                                                                                                                                           | Exploitability | Module # | Class # | Module                                    | Location                             |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------|--------------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> dir2 is constructed dynamically at runtime and no direct user input is used here.                                                                                                                                                        |                |          |         |                                           |                                      |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                        |                |          |         |                                           |                                      |
| 577                                                                                                                                                                                                                                                                                               | Neutral        | 83       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runscheduler.log.cs 122          |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                         |                |          |         |                                           |                                      |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fileName is constructed dynamically at runtime and no direct user input is used here.                                                                                                                                                    |                |          |         |                                           |                                      |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                       |                |          |         |                                           |                                      |
| 569                                                                                                                                                                                                                                                                                               | Neutral        | 84       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runscheduler.scheduler.cs 2846   |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                         |                |          |         |                                           |                                      |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> errorLogFileName is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                           |                                      |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                        |                |          |         |                                           |                                      |
| 575                                                                                                                                                                                                                                                                                               | Neutral        | 84       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runscheduler.scheduler.cs 2882   |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                         |                |          |         |                                           |                                      |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> dir2 is constructed dynamically at run time, and no direct user input from UI used here                                                                                                                                                  |                |          |         |                                           |                                      |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                        |                |          |         |                                           |                                      |
| 567                                                                                                                                                                                                                                                                                               | Neutral        | 85       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runscheduler.validationlog.cs 64 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                         |                |          |         |                                           |                                      |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> FileName is constructed dynamically and a part of filepath and name is hardcoded. There is no direct user input from UI used here                                                                                                        |                |          |         |                                           |                                      |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                       |                |          |         |                                           |                                      |
| 568                                                                                                                                                                                                                                                                                               | Neutral        | 85       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../runscheduler.validationlog.cs 78 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application.                                                                                                                                                                                  |                |          |         |                                           |                                      |

| Flaw Id                                                                                                                                                                                                                                                                                                          | Exploitability | Module # | Class # | Module                                                  | Location                         |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------|----------------------------------|
| All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                                                                         |                |          |         |                                                         |                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                                                             |                |          |         |                                                         |                                  |
| 960                                                                                                                                                                                                                                                                                                              | Neutral        | 86       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.userinterfa<br>ce.dll | .../runstatusform.cs 2381        |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                          |                |          |         |                                                         |                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                                                             |                |          |         |                                                         |                                  |
| 566                                                                                                                                                                                                                                                                                                              | Neutral        | 87       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll               | .../runsummarygenerator.cs 129   |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                                         |                                  |
| ↳ Potential False Positive (Danaher Corporation): errorLogFileName is fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature.                       |                |          |         |                                                         |                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                                         |                                  |
| 654                                                                                                                                                                                                                                                                                                              | Neutral        | 88       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll               | .../jobs/savetimepointjob.cs 205 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                                         |                                  |
| ↳ Potential False Positive (Danaher Corporation): a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                         |                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                                         |                                  |
| 655                                                                                                                                                                                                                                                                                                              | Neutral        | 88       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll               | .../jobs/savetimepointjob.cs 250 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                                         |                                  |
| ↳ Potential False Positive (Danaher Corporation): a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                         |                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                                         |                                  |
| 656                                                                                                                                                                                                                                                                                                              | Neutral        | 88       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll               | .../jobs/savetimepointjob.cs 258 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application.                                                                                                                                                                                                        |                |          |         |                                                         |                                  |

| Flaw Id                                                                                                                                                                                                                                                                                                                 | Exploitability | Module # | Class # | Module                                    | Location                               |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------|----------------------------------------|
| All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                                                                                |                |          |         |                                           |                                        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                           |                                        |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                           |                                        |
| 657                                                                                                                                                                                                                                                                                                                     | Neutral        | 88       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../jobs/savetimepointjob.cs 288       |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fullFileName is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                                                                         |                |          |         |                                           |                                        |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                             |                |          |         |                                           |                                        |
| 665                                                                                                                                                                                                                                                                                                                     | Neutral        | 88       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../jobs/savetimepointjob.cs 308       |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> A part of it is program generated and rest are fetched from built in browser.                                                                                                                                                                                  |                |          |         |                                           |                                        |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                           |                                        |
| 624                                                                                                                                                                                                                                                                                                                     | Neutral        | 89       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../protocol/selectechoprotocol.cs 201 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> path is program generated and a part of it uses in built folder browser dialog to fetch the folder name.                                                                                                                                                       |                |          |         |                                           |                                        |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                           |                                        |
| 641                                                                                                                                                                                                                                                                                                                     | Neutral        | 90       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../jobs/staggerplatejob.cs 217        |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                           |                                        |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                           |                                        |

| Flaw Id                                                                                                                                                                                                                                                                                                                 | Exploitability | Module # | Class # | Module                                    | Location                        |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------|---------------------------------|
| 642                                                                                                                                                                                                                                                                                                                     | Neutral        | 90       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../jobs/staggerplatejob.cs 225 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fullFileName is created dynamically and is not a direct user input from UI                                                                                                                                                                                     |                |          |         |                                           |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                             |                |          |         |                                           |                                 |
| 650                                                                                                                                                                                                                                                                                                                     | Neutral        | 90       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../jobs/staggerplatejob.cs 262 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                           |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                           |                                 |
| 651                                                                                                                                                                                                                                                                                                                     | Neutral        | 90       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../jobs/staggerplatejob.cs 343 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                           |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                           |                                 |
| 652                                                                                                                                                                                                                                                                                                                     | Neutral        | 90       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../jobs/staggerplatejob.cs 426 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                           |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                           |                                 |
| 653                                                                                                                                                                                                                                                                                                                     | Neutral        | 90       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../jobs/staggerplatejob.cs 447 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                           |                                 |

| Flaw Id                                                                                                                                                                                                                                                                                                                 | Exploitability | Module # | Class # | Module                                                                                             | Location                        |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|----------------------------------------------------------------------------------------------------|---------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> a part of it is hardcoded in the file and fetched from the saved settings. The directory value is the logdirectory which is set from UI, file name is hardcoded in the code. To get this directory value we use windows in built Browse for directory feature. |                |          |         |                                                                                                    |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                                                                                    |                                 |
| 177                                                                                                                                                                                                                                                                                                                     | V.Likely       | -        | 29      | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/becls.pod.devices.echoms.dll/sciexdevic<br>eprotos.dll | string Write() 22%              |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                                                                                    |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                             |                |          |         |                                                                                                    |                                 |
| 178                                                                                                                                                                                                                                                                                                                     | V.Likely       | -        | 29      | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/becls.pod.devices.echoms.dll/sciexdevic<br>eprotos.dll | string Write() 32%              |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                                                                                    |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Device parameters are passed from Tempo, and these are not untrusted UI user inputs. These values are stored in DB and fetched from there.                                                                                                                     |                |          |         |                                                                                                    |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                              |                |          |         |                                                                                                    |                                 |
| 179                                                                                                                                                                                                                                                                                                                     | V.Likely       | -        | 29      | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/becls.pod.devices.echoms.dll/sciexdevic<br>eprotos.dll | string Write() 55%              |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                                                                                    |                                 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> BatchFilePath is not a direct user input and doesnt contain any untrusted data. Tempo fetches this value from DB and the same is passed to this library.                                                                                                       |                |          |         |                                                                                                    |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                             |                |          |         |                                                                                                    |                                 |
| 180                                                                                                                                                                                                                                                                                                                     | V.Likely       | -        | 29      | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/becls.pod.devices.echoms.dll/sciexdevic<br>eprotos.dll | string Write() 59%              |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                               |                |          |         |                                                                                                    |                                 |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                             |                |          |         |                                                                                                    |                                 |
| 659                                                                                                                                                                                                                                                                                                                     | Likely         | 91       | -       | tempo.exe#3.2.0.1/la                                                                               | buildagent/.../surveyjob.cs 167 |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                                             | Location                                                         |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------------------------------------------------------------------|------------------------------------------------------------------|
|                                                                                                                                                                                                                                                           |                |          |         | bctype.pod.devices.ec<br>ho.dll                                                    |                                                                  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                                    |                                                                  |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> ouputDirectoryName uses the value from the built-in file browser.                                                                                                                                |                |          |         |                                                                                    |                                                                  |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented                                                                                                                                                                                 |                |          |         |                                                                                    |                                                                  |
| 660                                                                                                                                                                                                                                                       | Neutral        | 91       | -       | tempo.exe#3.2.0.1/la<br>bctype.pod.devices.ec<br>ho.dll                            | buildagent/.../surveyjob.cs 249                                  |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                                    |                                                                  |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fullPath uses the value from the built-in file browser.                                                                                                                                          |                |          |         |                                                                                    |                                                                  |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented                                                                                                                                                                                 |                |          |         |                                                                                    |                                                                  |
| 970                                                                                                                                                                                                                                                       | Neutral        | -        | 9       | tempo.exe#3.2.0.1/la<br>bctype.echo.api.dll                                        | System.Collections.Generic.Dictionary<int, string> a(string) 47% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                                    |                                                                  |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                    |                                                                  |
| 997                                                                                                                                                                                                                                                       | Neutral        | -        | 16      | tempo.exe#3.2.0.1/la<br>bctype.echo.api.dll                                        | System.Xml.XmlDocument<br>get_SurveyResults() 22%                |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                                    |                                                                  |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                    |                                                                  |
| 935                                                                                                                                                                                                                                                       | Neutral        | -        | 18      | tempo.exe#3.2.0.1/la<br>bctype.echo.api.dll                                        | System.Xml.XmlDocument<br>get_TransferResults() 26%              |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                                    |                                                                  |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                    |                                                                  |
| 526                                                                                                                                                                                                                                                       | Neutral        | 92       | -       | tempo.exe#3.2.0.1/la<br>bctype.pod.api.dll/lab<br>bctype.pod.driver.storex.<br>dll | .../tcpclientsocketsim.cs 115                                    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                                    |                                                                  |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                                     | Exploitability | Module # | Class # | Module                                                                                         | Location                      |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------------------------------------------------------------------------------|-------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> data is not a direct value from text box, instead its dynamically generated using split function.                                                                                                                                                                                                                  |                |          |         |                                                                                                |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                                                |                               |
| 912                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | 93       | -       | teleshakeclient.exe#3.2.0.1/labcyte.pod.devices.teleshake.dll                                  | .../teleshakedevice.cs 139    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path. |                |          |         |                                                                                                |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                                                |                               |
| 911                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | 93       | -       | teleshakeclient.exe#3.2.0.1/labcyte.pod.devices.teleshake.dll                                  | .../teleshakedevice.cs 214    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                |          |         |                                                                                                |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                                                |                               |
| 523                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | 95       | -       | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll                                                | .../ucechoapplications.cs 868 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                                                                |                               |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> IblEchoProtocol.Text is validated using its extension(.echo) to make sure that only echo files are allowed in the text box.                                                                                                                                                                                        |                |          |         |                                                                                                |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                                                |                               |
| 175                                                                                                                                                                                                                                                                                                                                                                         | V.Likely       | -        | 29      | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/becls.pod.devices.echoms.dll/sciexdeviceprotos.dll | void !ctor(string) 40%        |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                                                                |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                                                |                               |
| 949                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | -        | 11      | tempo.exe#3.2.0.1/labcyte.echo.api.dll                                                         | void a() 95%                  |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                |          |         |                                                                                                |                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                                                |                               |

| Flaw Id                                                                                                                                                                                                                                                         | Exploitability | Module # | Class # | Module                                                                                                                          | Location                                                                                                      |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|
| 903                                                                                                                                                                                                                                                             | Neutral        | -        | 10      | ambistoreclient.exe#3.2.0.1/labcyte.pod.d<br>river.ambistore.dll/la<br>bcyte.pod.common.dl<br>l/labcyte.echo.api.dll            | void<br>a(Labcyte.Echo.Common.Output.OutputD<br>ata,<br>Labcyte.Echo.Common.Output.OutputLog<br>, string) 62% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                  |                |          |         |                                                                                                                                 |                                                                                                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                     |                |          |         |                                                                                                                                 |                                                                                                               |
| 904                                                                                                                                                                                                                                                             | Neutral        | -        | 10      | ambistoreclient.exe#3.2.0.1/labcyte.pod.d<br>river.ambistore.dll/la<br>bcyte.pod.common.dl<br>l/labcyte.echo.api.dll            | void<br>a(Labcyte.Echo.Common.Output.OutputD<br>ata,<br>Labcyte.Echo.Common.Output.OutputLog<br>, string) 90% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The directory path is based on the local machine hostname in a desktop application and is used only for local app data/logging storage. Server access or remote path manipulation is not possible here |                |          |         |                                                                                                                                 |                                                                                                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                     |                |          |         |                                                                                                                                 |                                                                                                               |
| 908                                                                                                                                                                                                                                                             | Neutral        | -        | 9       | ambistoreclient.exe#3.2.0.1/labcyte.pod.d<br>river.ambistore.dll/la<br>bcyte.pod.common.dl<br>l/labcyte.echo.api.dll            | void a(object, System.EventArgs) 40%                                                                          |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                  |                |          |         |                                                                                                                                 |                                                                                                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                     |                |          |         |                                                                                                                                 |                                                                                                               |
| 948                                                                                                                                                                                                                                                             | Neutral        | -        | 11      | tempo.exe#3.2.0.1/la<br>bcyte.echo.api.dll                                                                                      | void Add(PexLogItem) 48%                                                                                      |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                  |                |          |         |                                                                                                                                 |                                                                                                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                     |                |          |         |                                                                                                                                 |                                                                                                               |
| 910                                                                                                                                                                                                                                                             | Neutral        | -        | 9       | ambistoreclient.exe#3.2.0.1/labcyte.pod.d<br>river.ambistore.dll/la<br>bcyte.pod.common.dl<br>l/labcyte.echo.api.dll            | void b(object, System.EventArgs) 81%                                                                          |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                  |                |          |         |                                                                                                                                 |                                                                                                               |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                     |                |          |         |                                                                                                                                 |                                                                                                               |
| 902                                                                                                                                                                                                                                                             | Neutral        | -        | 10      | temposcheduleviewew<br>r.exe#3.2.0.1/labcyte<br>.pod.scheduleviewer.<br>dll/labcyte.pod.comm<br>on.dll/labcyte.echo.a<br>pi.dll | void CreateOutputDirectories() 39%                                                                            |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                                                                     | Location                           |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------------------------------------------------------------------------------------------|------------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                                                            |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                                            |                                    |
| 901                                                                                                                                                                                                                                                       | Neutral        | -        | 10      | temposcheduleviewer.exe#3.2.0.1/labcyte.pod.scheduleviewer.dll/labcyte.pod.common.dll/labcyte.echo.api.dll | void CreateOutputDirectories() 57% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                                                            |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                                            |                                    |
| 900                                                                                                                                                                                                                                                       | Neutral        | -        | 10      | temposcheduleviewer.exe#3.2.0.1/labcyte.pod.scheduleviewer.dll/labcyte.pod.common.dll/labcyte.echo.api.dll | void CreateOutputDirectories() 89% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                                                            |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                                            |                                    |
| 162                                                                                                                                                                                                                                                       | V.Likely       | -        | 23      | tempo.exe#3.2.0.1/labcyte.pod.userinterface.dll/labcyte.echo.labware.dll                                   | void Load(string) 42%              |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                                                            |                                    |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> fname is the value passed to the external library, and this value is constructed dynamically in code, and it doesn't contain any untrusted data. fname is not a complete UI value.               |                |          |         |                                                                                                            |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                                            |                                    |
| 952                                                                                                                                                                                                                                                       | Neutral        | -        | 17      | tempoadmin.exe#3.2.0.1/labcyte.pod.common.dll/labcyte.echo.api.dll                                         | void Load(string) 47%              |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.            |                |          |         |                                                                                                            |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                               |                |          |         |                                                                                                            |                                    |
| 947                                                                                                                                                                                                                                                       | Neutral        | -        | 15      | tempoadmin.exe#3.2.0.1/labcyte.pod.com                                                                     | void Load(string) 71%              |

| Flaw Id | Exploitability | Module # | Class # | Module                                                                                                                                                                                                                                                                                                                                                                      | Location                 |
|---------|----------------|----------|---------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------|
|         |                |          |         | mon.dll/labcyte.echo.<br>api.dll                                                                                                                                                                                                                                                                                                                                            |                          |
|         |                |          |         | ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                          |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                          |
| 938     | Neutral        | -        | 14      | tempo.exe#3.2.0.1/la<br>bcyte.echo.api.dll                                                                                                                                                                                                                                                                                                                                  | void LoadCSV(string) 91% |
|         |                |          |         | ↳ <i>Potential False Positive (Danaher Corporation):</i> The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path. |                          |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                          |
| 134     | Neutral        | -        | 28      | tempo.exe#3.2.0.1/la<br>bcyte.pod.api.dll/becl<br>s.pod.devices.echoms<br>.dll/sciexdeviceclient.<br>dll                                                                                                                                                                                                                                                                    | void MoveNext() 34%      |
|         |                |          |         | ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                          |
|         |                |          |         | ↳ <i>Potential False Positive (Danaher Corporation):</i> The value for DeviceName originated from DB, and doesn't contain any untrusted user input.                                                                                                                                                                                                                         |                          |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                                                                  |                          |
| 135     | Neutral        | -        | 28      | tempoadmin.exe#3.2.<br>0.1/labcyte.pod.api.dl<br>l/becls.pod.devices.e<br>choms.dll/sciexdevic<br>eclient.dll                                                                                                                                                                                                                                                               | void MoveNext() 36%      |
|         |                |          |         | ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                          |
|         |                |          |         | ↳ <i>Potential False Positive (Danaher Corporation):</i> The value for DeviceName originated from DB and doesn't contain any untrusted user input.                                                                                                                                                                                                                          |                          |
|         |                |          |         | ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                                                                                                                                                  |                          |
| 136     | Neutral        | -        | 28      | tempoadmin.exe#3.2.<br>0.1/labcyte.pod.api.dl<br>l/becls.pod.devices.e<br>choms.dll/sciexdevic<br>eclient.dll                                                                                                                                                                                                                                                               | void MoveNext() 37%      |
|         |                |          |         | ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                          |
|         |                |          |         | ↳ <i>Potential False Positive (Danaher Corporation):</i> The value for DeviceName originated from DB, and                                                                                                                                                                                                                                                                   |                          |

| Flaw Id                                                                                                                                                                                                                                 | Exploitability | Module # | Class # | Module                                                                           | Location              |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|----------------------------------------------------------------------------------|-----------------------|
| doesn't contain any untrusted user input.                                                                                                                                                                                               |                |          |         |                                                                                  |                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                     |                |          |         |                                                                                  |                       |
| 943                                                                                                                                                                                                                                     | Neutral        | -        | 13      | tempoadmin.exe#3.2.0.1/labcyte.pod.com<br>mon.dll/labcyte.echo.<br>api.dll       | void Save(string) 25% |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                  |                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                  |                       |
| 972                                                                                                                                                                                                                                     | Neutral        | -        | 23      | tempo.exe#3.2.0.1/labcyte.pod.userinterfa<br>ce.dll/labcyte.echo.la<br>bware.dll | void Save(string) 30% |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                  |                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                  |                       |
| 946                                                                                                                                                                                                                                     | Neutral        | -        | 15      | tempoadmin.exe#3.2.0.1/labcyte.pod.com<br>mon.dll/labcyte.echo.<br>api.dll       | void Save(string) 42% |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                  |                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                  |                       |
| 944                                                                                                                                                                                                                                     | Neutral        | -        | 13      | tempoadmin.exe#3.2.0.1/labcyte.pod.com<br>mon.dll/labcyte.echo.<br>api.dll       | void Save(string) 49% |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                  |                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                  |                       |
| 971                                                                                                                                                                                                                                     | Neutral        | -        | 23      | tempo.exe#3.2.0.1/labcyte.pod.userinterfa<br>ce.dll/labcyte.echo.la<br>bware.dll | void Save(string) 50% |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                  |                       |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                    |                |          |         |                                                                                  |                       |
| 951                                                                                                                                                                                                                                     | Neutral        | -        | 17      | tempoadmin.exe#3.2.0.1/labcyte.pod.com<br>mon.dll/labcyte.echo.<br>api.dll       | void Save(string) 64% |

| Flaw Id                                                                                                                                                                                                                                        | Exploitability | Module # | Class # | Module                                                                                                                                    | Location                                     |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |                |          |         |                                                                                                                                           |                                              |
| 950                                                                                                                                                                                                                                            | Neutral        | -        | 17      | tempoadmin.exe#3.2.0.1/labcyte.pod.common.dll/labcyte.echo.api.dll                                                                        | void Save(string) 87%                        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |                |          |         |                                                                                                                                           |                                              |
| 945                                                                                                                                                                                                                                            | Neutral        | -        | 15      | tempoadmin.exe#3.2.0.1/labcyte.pod.common.dll/labcyte.echo.api.dll                                                                        | void Save(string) 87%                        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |                |          |         |                                                                                                                                           |                                              |
| 895                                                                                                                                                                                                                                            | Neutral        | -        | 10      | bioshakeclient.exe#3.2.0.1/labcyte.pod.driver.bioshake.dll/labcyte.pod.drivers.interfaces.dll/labcyte.pod.common.dll/labcyte.echo.api.dll | void SaveAppData(string, string, string) 21% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |                |          |         |                                                                                                                                           |                                              |
| 896                                                                                                                                                                                                                                            | Neutral        | -        | 10      | ambistoreclient.exe#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.echo.api.dll                                  | void SaveAppData(string, string, string) 40% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here. |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                    |                |          |         |                                                                                                                                           |                                              |
| 894                                                                                                                                                                                                                                            | Neutral        | -        | 10      | bioshakeclient.exe#3.2.0.1/labcyte.pod.driver.bioshake.dll/labcyte.pod.drivers.interfaces.dll/labcyte.pod.common.dll/labcyte.echo.api.dll | void SaveAppData(string, string, string) 57% |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                                     | Exploitability | Module # | Class # | Module                                                                                                                                    | Location                                     |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------|
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                                                                                           |                                              |
| 899                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | -        | 10      | bioshakeclient.exe#3.2.0.1/labcyte.pod.driver.bioshake.dll/labcyte.pod.drivers.interfaces.dll/labcyte.pod.common.dll/labcyte.echo.api.dll | void SaveAppData(string, string, string) 76% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                                                                                           |                                              |
| 898                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | -        | 10      | bioshakeclient.exe#3.2.0.1/labcyte.pod.driver.bioshake.dll/labcyte.pod.drivers.interfaces.dll/labcyte.pod.common.dll/labcyte.echo.api.dll | void SaveAppData(string, string, string) 78% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                                                                                           |                                              |
| 897                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | -        | 10      | ambistoreclient.exe#3.2.0.1/labcyte.pod.driver.ambistore.dll/labcyte.pod.common.dll/labcyte.echo.api.dll                                  | void SaveAppData(string, string, string) 80% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                                                                                           |                                              |
| 936                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | -        | 14      | tempo.exe#3.2.0.1/labcyte.echo.api.dll                                                                                                    | void SaveCSV(string, int) 18%                |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path. |                |          |         |                                                                                                                                           |                                              |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                                                                                           |                                              |
| 937                                                                                                                                                                                                                                                                                                                                                                         | Neutral        | -        | 14      | tempo.exe#3.2.0.1/labcyte.echo.api.dll                                                                                                    | void SaveCSV(string, int) 56%                |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server                                                                                                                                                  |                |          |         |                                                                                                                                           |                                              |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                              | Exploitability | Module # | Class # | Module                                                                  | Location                     |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------------------------|------------------------------|
| path involved here.                                                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                         |                              |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                         |                              |
| ✓ 806                                                                                                                                                                                                                                                                                                                                                                | Neutral        | 96       | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/becls.pod.devices.volumecheck.dll | .../volumecheckdevice.cs 163 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                                                                                                   |                |          |         |                                                                         |                              |
| ↳ Potential False Positive (Danaher Corporation): outputFilePath is derived dynamically at runtime and is not a direct user input from UI.                                                                                                                                                                                                                           |                |          |         |                                                                         |                              |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                                                                                                                 |                |          |         |                                                                         |                              |
| ✓ 920                                                                                                                                                                                                                                                                                                                                                                | Neutral        | 97       | -       | robotteachingutility.exe#3.2.0.1                                        | .../xcfgfileviewer.cs 43     |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                                                                                                              |                |          |         |                                                                         |                              |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                         |                              |
| ✓ 919                                                                                                                                                                                                                                                                                                                                                                | Neutral        | 97       | -       | robotteachingutility.exe#3.2.0.1                                        | .../xcfgfileviewer.cs 61     |
| ↳ Potential False Positive (Danaher Corporation): The value passed to StreamWriter comes directly from a local SaveFileDialog selection. It is a user-chosen save path in a Windows desktop application, not a path constructed from server data, remote input, or any other attacker-controlled source. The code does not read or fetch information from that path. |                |          |         |                                                                         |                              |
| ↳ Approve Mitigation (Danaher Corporation): Approve                                                                                                                                                                                                                                                                                                                  |                |          |         |                                                                         |                              |

## → Untrusted Search Path(1 flaw)

Associated Flaws by CWE ID:

## → Uncontrolled Search Path Element (CWE ID 427)(1 flaw)

### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                           | Exploitability | Module # | Class # | Module                   | Location                   |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------|----------------------------|
| ✓ 382                                                                                                                                                                                                                                             | Neutral        | 44       | -       | koyoclickplc.exe#3.2.0.1 | .../mainwindow.xaml.cs 603 |
| ↳ Mitigate by OS Environment (Danaher Corporation): This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer |                |          |         |                          |                            |
| ↳ Approve Mitigation (Danaher Corporation): Approved but for a different reason. Implementation does                                                                                                                                              |                |          |         |                          |                            |

| Flaw Id                  | Exploitability | Module # | Class # | Module | Location |
|--------------------------|----------------|----------|---------|--------|----------|
| not use untrusted input. |                |          |         |        |          |

→ Information Leakage(5 flaws)

Associated Flaws by CWE ID:

→ Improper Restriction of XML External Entity Reference (CWE ID 611)(5 flaws)

Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                                        | Exploitability | Module # | Class # | Module                                  | Location                                         |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-----------------------------------------|--------------------------------------------------|
| ✓ 811                                                                                                                                                                                                                                                                          | Neutral        | 14       | -       | tempoadmin.exe#3.2.0.1                  | .../deviceinfoviewmodel.cs 446                   |
| ↳ Mitigate by Design (Danaher Corporation): Code uses XmlReader.Create() with XmlReaderSettings before loading into XmlDocument, Input is from local file dialog selection only (trusted source) in internal desktop application.                                              |                |          |         |                                         |                                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                            |                |          |         |                                         |                                                  |
| ✓ 925                                                                                                                                                                                                                                                                          | Neutral        | 20       | -       | cytomatclient.exe#3.2.0.1               | buildagent/.../frmdevcaps.cs 103                 |
| ↳ Mitigate by Design (Danaher Corporation): The input file is selected locally through OpenFileDialog, not supplied through a remote server-facing interface. Given parser defaults and local-only input path, this issue is mitigated by design for this application context. |                |          |         |                                         |                                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                           |                |          |         |                                         |                                                  |
| ✓ 435                                                                                                                                                                                                                                                                          | Neutral        | 23       | -       | cytomatclient.exe#3.2.0.1               | buildagent/.../frmmain.cs 1202                   |
| ↳ Mitigate by Design (Danaher Corporation): Code uses XmlReader.Create() with XmlReaderSettings before loading into XmlDocument, Input is from local file dialog selection only (trusted source) in internal desktop application.                                              |                |          |         |                                         |                                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented.                                                                                                                                                                                                            |                |          |         |                                         |                                                  |
| ✓ 996                                                                                                                                                                                                                                                                          | Neutral        | -        | 16      | tempo.exe#3.2.0.1/la bcyte.echo.api.dll | System.Xml.XmlDocument get_SurveyResults() 22%   |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                        |                |          |         |                                         |                                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                           |                |          |         |                                         |                                                  |
| ✓ 934                                                                                                                                                                                                                                                                          | Neutral        | -        | 18      | tempo.exe#3.2.0.1/la bcyte.echo.api.dll | System.Xml.XmlDocument get_TransferResults() 26% |
| ↳ Potential False Positive (Danaher Corporation): Unauthorized access to files on server is not applicable in this context, as the application is deployed on a local machine, not on a webserver. No remote server path involved here.                                        |                |          |         |                                         |                                                  |
| ↳ Approve Mitigation (Danaher Corporation): Approved                                                                                                                                                                                                                           |                |          |         |                                         |                                                  |

| Flaw Id | Exploitability | Module # | Class # | Module | Location |
|---------|----------------|----------|---------|--------|----------|
|---------|----------------|----------|---------|--------|----------|

### → Cryptographic Issues(2 flaws)

Associated Flaws by CWE ID:

### → Insufficient Entropy (CWE ID 331)(2 flaws)

#### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                      | Exploitability | Module # | Class # | Module                                                                                                                             | Location                                                                                                             |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------|
| 114                                                                                                                                                                                                                                          | Unlikely       | -        | 4       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll/becl<br>s.pod.devices.volume<br>check.dll/beckmanco<br>ulter.devices.volumec<br>heck.dll | double[] GetPlateVolumes(string, int, int,<br>System.Threading.CancellationToken,<br>System.Action<string, int>) 81% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> System.Random is used for non-security purposes. These are NOT used for cryptographic operations, session keys, or authentication                                                   |                |          |         |                                                                                                                                    |                                                                                                                      |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> This internally uses another dll, which has random used. But this is not used for any security purposes, and NOT used for cryptographic operations, session keys, or authentication |                |          |         |                                                                                                                                    |                                                                                                                      |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved. Verified in 3rd party code                                                                                                                                                      |                |          |         |                                                                                                                                    |                                                                                                                      |
| 551                                                                                                                                                                                                                                          | Unlikely       | 30       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.common.dl<br>                                                                                    | .../jobdevicepoolingcalculations.cs 72                                                                               |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> System.Random is used for non-security purposes. These are NOT used for cryptographic operations, session keys, or authentication                                                   |                |          |         |                                                                                                                                    |                                                                                                                      |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved as documented.                                                                                                                                                                   |                |          |         |                                                                                                                                    |                                                                                                                      |

### → Time and State(1 flaw)

Associated Flaws by CWE ID:

### → Insecure Temporary File (CWE ID 377)(1 flaw)

#### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                     | Location                              |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|----------------------------|---------------------------------------|
| 810                                                                                                                                                                                                | Neutral        | 36       | -       | tempoadmin.exe#3.2.<br>0.1 | .../viewmodels/licenseviewmodel.cs 67 |
| ↳ <i>Mitigate by Design (Danaher Corporation):</i> File contains license information already displayed in the UI (non-secret to the logged-in user). Access restricted to user's profile directory |                |          |         |                            |                                       |

| Flaw Id                                                            | Exploitability | Module # | Class # | Module | Location |
|--------------------------------------------------------------------|----------------|----------|---------|--------|----------|
| ↳ Approve Mitigation (Danaher Corporation): Approved as documented |                |          |         |        |          |

→ Server Configuration(3 flaws)

Associated Flaws by CWE ID:

→ Unintended Proxy or Intermediary ('Confused Deputy') (CWE ID 441)(3 flaws)

Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    | Exploitability | Module # | Class # | Module                                                                 | Location                    |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------------------------------------------------------|-----------------------------|
|  735                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | Neutral        | 2        | -       | ambistoreclient.exe#3.2.0.1/labcyte.pod.driver.ambistore.dll           | .../ambistoredriver.cs 526  |
| ↳ Mitigate by Design (Danaher Corporation): All files on same local machine. This establishes point-to-point connections with owned laboratory equipment. No remediation required.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         |                |          |         |                                                                        |                             |
| ↳ Approve Mitigation (Danaher Corporation): Approved for functionality                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |                |          |         |                                                                        |                             |
|  555                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      | Neutral        | 3        | -       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.driver.ambistore.dll | .../ambistoredriver2.cs 294 |
| ↳ Mitigate by Design (Danaher Corporation): The Socket.Connect operation is used for communication with devices directly connected to the local machine, or within a highly restricted local network segment. This communication is not exposed to external networks (e.g., the internet). The "client" providing the data is an internal component of the application, or a specific, controlled hardware device. This implies a lower likelihood of malicious data injection from an untrusted source. For a malicious payload to be introduced, an attacker would likely require physical access to the local machine or the connected device, adding a significant barrier to exploitation. It is not propagating malicious data to other servers or external systems. |                |          |         |                                                                        |                             |
| ↳ Approve Mitigation (Danaher Corporation): Approved for functionality                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |                |          |         |                                                                        |                             |
|  676                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    | Neutral        | 46       | -       | microserveclient.exe#3.2.0.1/labcyte.pod.driver.microserve.dll         | .../microservedriver.cs 442 |
| ↳ Mitigate by Design (Danaher Corporation): All files on same local machine. This establishes point-to-point connections with owned laboratory equipment. No remediation required.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         |                |          |         |                                                                        |                             |
| ↳ Approve Mitigation (Danaher Corporation): Approved for functionality                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |                |          |         |                                                                        |                             |

## Low (5 flaws)

### Information Leakage(5 flaws)

#### Associated Flaws by CWE ID:

### Generation of Error Message Containing Sensitive Information (CWE ID 209)(3 flaws)

#### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                                                      | Exploitability | Module # | Class # | Module                                            | Location                           |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------|------------------------------------|
| 501                                                                                                                                                                                                                                                                                          | Neutral        | 40       | -       | tempotester.exe#3.2.0.1/labcyte.common.mailer.dll | buildagent/.../mailservice.cs 552  |
| ↳ <i>Mitigate by Design (Danaher Corporation):</i> Tempo has internal email notification system only. All recipients are internal who require detailed information for troubleshooting. Email addresses configured administratively, not user-supplied.                                      |                |          |         |                                                   |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved. Information remains within trusted envelope.                                                                                                                                                                                    |                |          |         |                                                   |                                    |
| 518                                                                                                                                                                                                                                                                                          | Neutral        | 40       | -       | tempotester.exe#3.2.0.1/labcyte.common.mailer.dll | buildagent/.../mailservice.cs 552  |
| ↳ <i>Mitigate by Design (Danaher Corporation):</i> Tempo has internal email notification system only. All recipients are internal who require detailed information for troubleshooting. Email addresses configured administratively, not user-supplied.                                      |                |          |         |                                                   |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved. Information remains within trusted envelope.                                                                                                                                                                                    |                |          |         |                                                   |                                    |
| 500                                                                                                                                                                                                                                                                                          | Neutral        | 61       | -       | tempo.exe#3.2.0.1                                 | buildagent/.../protocolview.cs 480 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> Detailed error messages are necessary for operational support and do not create security vulnerabilities as this is a desktop application. All files on same local machine. User physically present and authenticated by Windows. |                |          |         |                                                   |                                    |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved. Information remains within trusted envelope.                                                                                                                                                                                    |                |          |         |                                                   |                                    |

### Insertion of Sensitive Information Into Sent Data (CWE ID 201)(2 flaws)

#### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                 | Exploitability | Module # | Class # | Module                                            | Location                          |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|---------------------------------------------------|-----------------------------------|
| 519                                                                                                                                                                                                                                                     | Unlikely       | 40       | -       | tempotester.exe#3.2.0.1/labcyte.common.mailer.dll | buildagent/.../mailservice.cs 552 |
| ↳ <i>Mitigate by Design (Danaher Corporation):</i> Tempo has internal email notification system only. All recipients are internal who require detailed information for troubleshooting. Email addresses configured administratively, not user-supplied. |                |          |         |                                                   |                                   |
| ↳ <i>Approve Mitigation (Danaher Corporation):</i> Approved. Information remains within trusted envelope.                                                                                                                                               |                |          |         |                                                   |                                   |
| 534                                                                                                                                                                                                                                                     | Unlikely       | 40       | -       | tempotester.exe#3.2.0.1/labcyte.common            | buildagent/.../mailservice.cs 552 |

| Flaw Id                                                                                                                                                                                                                                          | Exploitability | Module # | Class # | Module | Location |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------|----------|
| mailer.dll                                                                                                                                                                                                                                       |                |          |         |        |          |
| ↳ Mitigate by Design (Danaher Corporation): Tempo has internal email notification system only. All recipients are internal who require detailed information for troubleshooting. Email addresses configured administratively, not user-supplied. |                |          |         |        |          |
| ↳ Approve Mitigation (Danaher Corporation): Approved. Information remains within trusted envelope.                                                                                                                                               |                |          |         |        |          |

Appendix C: Proposed Mitigated Flaws (by Danaher Corporation)

NOTE: Except in circumstances where the Customer has purchased Veracode's Mitigation Proposal Review Solution, Veracode does not review the mitigation strategy described below and is not responsible for its contents or the accuracy of any statements provided.

Very High (4 flaws) Fix Required by Policy

→ Command or Argument Injection(4 flaws)

Associated Flaws by CWE ID:

→ Improper Neutralization of Special Elements used in an OS Command ('OS Command Injection') (CWE ID 78)(4 flaws)

Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                                                                                                          | Exploitability | Module # | Class # | Module                                                | Location                               |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------|----------------------------------------|
|  809                                                                                                                                                                                                                                                            | Neutral        | 14       | -       | tempoadmin.exe#3.2.0.1                                | .../deviceinfoviewmodel.cs 342         |
| <p>↳ Accept the risk (Danaher Corporation): We are accepting that they're vulnerability, but they're part of the design of the product. There's not a good mitigation without removing the functionality that's desired.</p>                                                                                                                     |                |          |         |                                                       |                                        |
| <p>↳ Accept the risk (Danaher Corporation): Added all possible validations for both executable path and command line parameters. We cannot whitelist it fully. hence accepting the risk and raised a Jira ticket- <a href="https://lsjira.beckman.com/browse/ISON-1338">https://lsjira.beckman.com/browse/ISON-1338</a></p>                      |                |          |         |                                                       |                                        |
|  582                                                                                                                                                                                                                                                           | Neutral        | 35       | -       | tempo.exe#3.2.0.1                                     | .../launchprogramcommand.cs 42         |
| <p>↳ Accept the risk (Danaher Corporation): We are accepting that they're vulnerability, but they're part of the design of the product. There's not a good mitigation without removing the functionality that's desired.</p>                                                                                                                     |                |          |         |                                                       |                                        |
| <p>↳ Accept the risk (Danaher Corporation): We are accepting the risk as these implementations are part of the design.</p> <p>Below Jira ticket has been raised - <a href="https://lsjira.beckman.com/browse/ISON-1338">https://lsjira.beckman.com/browse/ISON-1338</a></p>                                                                      |                |          |         |                                                       |                                        |
|  1005                                                                                                                                                                                                                                                         | Neutral        | 53       | -       | koyopl.exe#3.2.0.1/labcyte.pod.common.dll             | buildagent/.../device/poddevice.cs 772 |
| <p>↳ Accept the risk (Danaher Corporation): Added all possible validations for both executable path and command line parameters. We cannot whitelist it fully. hence accepting the risk and raised a Jira ticket- <a href="https://lsjira.beckman.com/browse/ISON-1338">https://lsjira.beckman.com/browse/ISON-1338</a></p>                      |                |          |         |                                                       |                                        |
|  823                                                                                                                                                                                                                                                          | Neutral        | 81       | -       | tempotester.exe#3.2.0.1/labcyte.pod.userinterface.dll | .../runreportsform.cs 750              |
| <p>↳ Accept the risk (Danaher Corporation): Command line parameters are constructed dynamically in code and executable path is fetched from registry. We cannot whitelist it fully, hence accepting the risk and raised a Jira ticket- <a href="https://lsjira.beckman.com/browse/ISON-1338">https://lsjira.beckman.com/browse/ISON-1338</a></p> |                |          |         |                                                       |                                        |

# Medium (18 flaws)

## Credentials Management(3 flaws)

Associated Flaws by CWE ID:

### Use of Hard-coded Password (CWE ID 259)(3 flaws)

#### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                     | Exploitability | Module # | Class # | Module                           | Location                             |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|----------------------------------|--------------------------------------|
| 384                                                                                                                                                                                                         | Likely         | 41       | -       | koyopl.exe#3.2.0.1               | buildagent/.../koyopl/mainform.cs 36 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The hardcoded password is used solely as a UI safety lock to prevent accidental UI values modifications, not for authentication or authorization.  |                |          |         |                                  |                                      |
| ↳ <i>Accept the risk (Danaher Corporation):</i> Accepting the risk and opened Jira ticket to track it <a href="https://lsjira.beckman.com/browse/ISON-1335">https://lsjira.beckman.com/browse/ISON-1335</a> |                |          |         |                                  |                                      |
| 682                                                                                                                                                                                                         | Likely         | 58       | -       | robotteachingutility.exe#3.2.0.1 | .../preciseteachingform.cs 27        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The hardcoded password is used solely as a UI safety lock to prevent accidental UI values modifications, not for authentication or authorization.  |                |          |         |                                  |                                      |
| ↳ <i>Accept the risk (Danaher Corporation):</i> Accepting the risk and opened Jira ticket to track it <a href="https://lsjira.beckman.com/browse/ISON-1335">https://lsjira.beckman.com/browse/ISON-1335</a> |                |          |         |                                  |                                      |
| 681                                                                                                                                                                                                         | Likely         | 72       | -       | robotteachingutility.exe#3.2.0.1 | .../rorzeteachingfrm.cs 18           |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The hardcoded password is used solely as a UI safety lock to prevent accidental UI values modifications, not for authentication or authorization.  |                |          |         |                                  |                                      |
| ↳ <i>Accept the risk (Danaher Corporation):</i> Accepting the risk and opened Jira ticket to track it <a href="https://lsjira.beckman.com/browse/ISON-1335">https://lsjira.beckman.com/browse/ISON-1335</a> |                |          |         |                                  |                                      |

## Cryptographic Issues(3 flaws)

Associated Flaws by CWE ID:

### Use of a Broken or Risky Cryptographic Algorithm (CWE ID 327)(3 flaws)

#### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                         | Exploitability | Module # | Class # | Module                                    | Location                               |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------|----------------------------------------|
| 625                                                                                                                                                                                                                             | Likely         | 89       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll | .../protocol/selectechoprotocol.cs 203 |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> MD5 is being used for file integrity checking (checksums), NOT for password hashing or cryptographic security.                                                         |                |          |         |                                           |                                        |
| ↳ <i>Accept the risk (Danaher Corporation):</i> Accepting the risk for now. Opened below Jira ticket for tracking the fix <a href="https://lsjira.beckman.com/browse/ISON-1337">https://lsjira.beckman.com/browse/ISON-1337</a> |                |          |         |                                           |                                        |

| Flaw Id                                                                                                                                                                                                                                                                                         | Exploitability | Module # | Class # | Module                                                         | Location                             |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|----------------------------------------------------------------|--------------------------------------|
| 524                                                                                                                                                                                                                                                                                             | Likely         | 95       | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.userinterfa<br>ce.dll        | .../ucechoapplications.cs 870        |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> MD5 is being used for file integrity checking (checksums), NOT for password hashing or cryptographic security.                                                                                                                         |                |          |         |                                                                |                                      |
| ↳ <i>Accept the risk (Danaher Corporation):</i> Accepting the risk for now. Opened below Jira ticket for tracking the fix <a href="https://lsjira.beckman.com/browse/ISON-1337">https://lsjira.beckman.com/browse/ISON-1337</a>                                                                 |                |          |         |                                                                |                                      |
| 160                                                                                                                                                                                                                                                                                             | Likely         | -        | 8       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll/btilh<br>crunner.dll | void a(ref string /*1*/, string) 45% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> MD5 is being used for file integrity checking (checksums), NOT for password hashing or cryptographic security.                                                                                                                         |                |          |         |                                                                |                                      |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> Reported the issue to Library maintainer and we are accepting the risk. Below Jira ticket has been created to track it.<br><a href="https://lsjira.beckman.com/browse/ISON-1333">https://lsjira.beckman.com/browse/ISON-1333</a> |                |          |         |                                                                |                                      |

## → Directory Traversal(5 flaws)

### Associated Flaws by CWE ID:

## → External Control of File Name or Path (CWE ID 73)(5 flaws)

### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                                                         | Exploitability | Module # | Class # | Module                                                         | Location                           |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|----------------------------------------------------------------|------------------------------------|
| 145                                                                                                                                                                                                                                                                                             | Neutral        | -        | 5       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll/btilh<br>crunner.dll | void a(object) 99%                 |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                       |                |          |         |                                                                |                                    |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> Reported the issue to Library maintainer and we are accepting the risk. Below Jira ticket has been created to track it.<br><a href="https://lsjira.beckman.com/browse/ISON-1333">https://lsjira.beckman.com/browse/ISON-1333</a> |                |          |         |                                                                |                                    |
| 141                                                                                                                                                                                                                                                                                             | Neutral        | -        | 5       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll/btilh<br>crunner.dll | void a(string, object, object) 98% |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer.                                       |                |          |         |                                                                |                                    |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> Reported the issue to Library maintainer and we are accepting the risk. Below Jira ticket has been created to track it.                                                                                                          |                |          |         |                                                                |                                    |

| Flaw Id                                                                                                                                                                                                                                                   | Exploitability | Module # | Class # | Module                                                         | Location                                              |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|----------------------------------------------------------------|-------------------------------------------------------|
| <a href="https://lsjira.beckman.com/browse/ISON-1333">https://lsjira.beckman.com/browse/ISON-1333</a>                                                                                                                                                     |                |          |         |                                                                |                                                       |
| 143                                                                                                                                                                                                                                                       | Neutral        | -        | 5       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll/btilh<br>crunner.dll | void a(string, object, object, object, object)<br>99% |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                |                                                       |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> Reported the issue to Library maintainer and we are accepting the risk. Below Jira ticket has been created to track it.                                                                    |                |          |         |                                                                |                                                       |
| <a href="https://lsjira.beckman.com/browse/ISON-1333">https://lsjira.beckman.com/browse/ISON-1333</a>                                                                                                                                                     |                |          |         |                                                                |                                                       |
| 140                                                                                                                                                                                                                                                       | Neutral        | -        | 5       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll/btilh<br>crunner.dll | void ao() 96%                                         |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                |                                                       |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> Reported the issue to Library maintainer and we are accepting the risk. Below Jira ticket has been created to track it.                                                                    |                |          |         |                                                                |                                                       |
| <a href="https://lsjira.beckman.com/browse/ISON-1333">https://lsjira.beckman.com/browse/ISON-1333</a>                                                                                                                                                     |                |          |         |                                                                |                                                       |
| 142                                                                                                                                                                                                                                                       | Neutral        | -        | 5       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll/btilh<br>crunner.dll | void q(string, object) 98%                            |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer. |                |          |         |                                                                |                                                       |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> Reported the issue to Library maintainer and we are accepting the risk. Below Jira ticket has been created to track it.                                                                    |                |          |         |                                                                |                                                       |
| <a href="https://lsjira.beckman.com/browse/ISON-1333">https://lsjira.beckman.com/browse/ISON-1333</a>                                                                                                                                                     |                |          |         |                                                                |                                                       |

## → Information Leakage(5 flaws)

Associated Flaws by CWE ID:

## → Improper Restriction of XML External Entity Reference (CWE ID 611)(2 flaws)

### Instances found via Static Scan

| Flaw Id                                                                                                         | Exploitability | Module # | Class # | Module                                                                              | Location                    |
|-----------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-------------------------------------------------------------------------------------|-----------------------------|
| 832                                                                                                             | Neutral        | 7        | -       | tempo.exe#3.2.0.1/la<br>bbyte.pod.api.dll/lab<br>byte.pod.devices.cer<br>tuflex.dll | .../certusflexdevice.cs 520 |
| ↳ <i>Accept the risk (Danaher Corporation):</i> We are accepting the risk for now and have raised a Jira ticket |                |          |         |                                                                                     |                             |

| Flaw Id                                                                                                                                                                                                                                                        | Exploitability | Module # | Class # | Module                 | Location                       |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------|--------------------------------|
| to track this issue. <a href="https://beckmancoulter.atlassian.net/browse/ISON-1388">https://beckmancoulter.atlassian.net/browse/ISON-1388</a>                                                                                                                 |                |          |         |                        |                                |
| 828                                                                                                                                                                                                                                                            | Neutral        | 14       | -       | tempoadmin.exe#3.2.0.1 | .../deviceinfoviewmodel.cs 432 |
| ↳ <i>Accept the risk (Danaher Corporation):</i> We are accepting the risk for now and have raised a Jira ticket to track this issue. <a href="https://beckmancoulter.atlassian.net/browse/ISON-1388">https://beckmancoulter.atlassian.net/browse/ISON-1388</a> |                |          |         |                        |                                |

## Server-Side Request Forgery (SSRF) (CWE ID 918)(3 flaws)

### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                            | Exploitability | Module # | Class # | Module                                                                                                           | Location            |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------------------------------------------------------------------------------------------------|---------------------|
| 211                                                                                                                                                                                                | Neutral        | -        | 3       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.samplestore.dll/beckmancoulter.devices.samplestore.dll | void MoveNext() 23% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The "MoveNext" reference is a compiler-generated async state machine method name, not a different vulnerability. No remediation required. |                |          |         |                                                                                                                  |                     |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> <a href="https://lsjira.beckman.com/browse/ISON-1334">https://lsjira.beckman.com/browse/ISON-1334</a>                               |                |          |         |                                                                                                                  |                     |
| 210                                                                                                                                                                                                | Neutral        | -        | 2       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.samplestore.dll/beckmancoulter.devices.samplestore.dll | void MoveNext() 27% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The "MoveNext" reference is a compiler-generated async state machine method name, not a different vulnerability. No remediation required. |                |          |         |                                                                                                                  |                     |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> <a href="https://lsjira.beckman.com/browse/ISON-1334">https://lsjira.beckman.com/browse/ISON-1334</a>                               |                |          |         |                                                                                                                  |                     |
| 209                                                                                                                                                                                                | Neutral        | -        | 1       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.samplestore.dll/beckmancoulter.devices.samplestore.dll | void MoveNext() 60% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The "MoveNext" reference is a compiler-generated async state machine method name, not a different vulnerability. No remediation required. |                |          |         |                                                                                                                  |                     |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> <a href="https://lsjira.beckman.com/browse/ISON-1334">https://lsjira.beckman.com/browse/ISON-1334</a>                               |                |          |         |                                                                                                                  |                     |

## Untrusted Search Path(2 flaws)

### Associated Flaws by CWE ID:

## Uncontrolled Search Path Element (CWE ID 427)(2 flaws)

### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                                                         | Exploitability | Module # | Class # | Module                                                 | Location                             |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|--------------------------------------------------------|--------------------------------------|
| 139                                                                                                                                                                                                                                                                                             | Neutral        | -        | 5       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/btilhcrunner.dll | System.Type a(string, string) 29%    |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer                                        |                |          |         |                                                        |                                      |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> Reported the issue to Library maintainer and we are accepting the risk. Below Jira ticket has been created to track it.<br><a href="https://lsjira.beckman.com/browse/ISON-1333">https://lsjira.beckman.com/browse/ISON-1333</a> |                |          |         |                                                        |                                      |
| 158                                                                                                                                                                                                                                                                                             | Neutral        | -        | 7       | tempo.exe#3.2.0.1/labcyte.pod.api.dll/btilhcrunner.dll | void b(object, System.EventArgs) 94% |
| ↳ <i>Mitigate by OS Environment (Danaher Corporation):</i> This is a desktop application, not a web application. All files on same local machine. User physically present and authenticated by Windows. User already has direct file access via Explorer                                        |                |          |         |                                                        |                                      |
| ↳ <i>Reported to Library Maintainer (Danaher Corporation):</i> Reported the issue to Library maintainer and we are accepting the risk. Below Jira ticket has been created to track it.<br><a href="https://lsjira.beckman.com/browse/ISON-1333">https://lsjira.beckman.com/browse/ISON-1333</a> |                |          |         |                                                        |                                      |

## Low (3 flaws)

### Information Leakage(3 flaws)

#### Associated Flaws by CWE ID:

### Insertion of Sensitive Information Into Sent Data (CWE ID 201)(3 flaws)

### Instances found via Static Scan

| Flaw Id                                                                                                                                                                                                                                                                                                                                        | Exploitability | Module # | Class # | Module                                                                                                                 | Location            |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|------------------------------------------------------------------------------------------------------------------------|---------------------|
| 499                                                                                                                                                                                                                                                                                                                                            | Unlikely       | -        | 2       | tempotester.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.samplestore.dll/beckmancoulter.devices.samplestore.dll | void MoveNext() 27% |
| ↳ <i>Potential False Positive (Danaher Corporation):</i> The data (serial number) is NOT sensitive. The fact that it originates from the Registry is irrelevant, as it contains non-sensitive data.                                                                                                                                            |                |          |         |                                                                                                                        |                     |
| ↳ <i>Accept the risk (Danaher Corporation):</i> It passes username and password to start the session and to communicate with the device. This is intended use and is expected. Accepting the risk and Jira ticket is raised for the same <a href="https://lsjira.beckman.com/browse/ISON-1334">https://lsjira.beckman.com/browse/ISON-1334</a> |                |          |         |                                                                                                                        |                     |

| Flaw Id                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                          | Exploitability | Module # | Class # | Module                                                                                                                | Location            |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|----------|---------|-----------------------------------------------------------------------------------------------------------------------|---------------------|
| 308                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              | Unlikely       | -        | 1       | tempoadmin.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.samplestore.dll/beckmancoulter.devices.samplestore.dll | void MoveNext() 60% |
| <p>↳ <i>Potential False Positive (Danaher Corporation):</i> The data (serial number) is NOT sensitive. The fact that it originates from the Registry is irrelevant, as it contains non-sensitive data.</p> <p>↳ <i>Accept the risk (Danaher Corporation):</i> It passes username and password to start the session and to communicate with the device. This is intended use and is expected. Accepting the risk and Jira ticket is raised for the same <a href="https://lsjira.beckman.com/browse/ISON-1334">https://lsjira.beckman.com/browse/ISON-1334</a></p> |                |          |         |                                                                                                                       |                     |
| 309                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              | Unlikely       | -        | 1       | tempoadmin.exe#3.2.0.1/labcyte.pod.api.dll/labcyte.pod.devices.samplestore.dll/beckmancoulter.devices.samplestore.dll | void MoveNext() 60% |
| <p>↳ <i>Potential False Positive (Danaher Corporation):</i> The data (serial number) is NOT sensitive. The fact that it originates from the Registry is irrelevant, as it contains non-sensitive data.</p> <p>↳ <i>Accept the risk (Danaher Corporation):</i> It passes username and password to start the session and to communicate with the device. This is intended use and is expected. Accepting the risk and Jira ticket is raised for the same <a href="https://lsjira.beckman.com/browse/ISON-1334">https://lsjira.beckman.com/browse/ISON-1334</a></p> |                |          |         |                                                                                                                       |                     |

## Appendix D: Referenced Source Files

| Id | Filename                        | Path                                                                                |
|----|---------------------------------|-------------------------------------------------------------------------------------|
| 1  | ambistoredevice.cs              | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.ambistore/          |
| 2  | ambistoredriver.cs              | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.driver.ambistore/           |
| 3  | ambistoredriver2.cs             | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.driver.ambistore/           |
| 4  | barcodeform.cs                  | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.userinterface/              |
| 5  | bioshakedevice.cs               | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.bioshake/           |
| 6  | bulkfilljob.cs                  | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.tempest/            |
| 7  | certusflexdevice.cs             | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.certusflex/         |
| 8  | clsmousekeyboardevents.cs       | buildagent/work/f81ed2ccc88a1e9a/thirdparty/agcsn/                                  |
| 9  | combinldevice.cs                | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.combinl/            |
| 10 | configfileencryptor.cs          | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.driver.commonrobot/config/  |
| 11 | configfileencryptor.cs          | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/podadmin/autogen/                     |
| 12 | cytomatdevice.cs                | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.cytomat/            |
| 13 | delayplatejob.cs                | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/multijobs/              |
| 14 | deviceinfoviewmodel.cs          | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/podadmin/viewmodels/                  |
| 15 | devicereportmanager.cs          | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.common/device/              |
| 16 | devicerundefinition.cs          | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.samplestore/rundef/ |
| 17 | devicetestform.cs               | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.userinterface/              |
| 18 | echodevice.cs                   | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.echo/               |
| 19 | filestreamlock.cs               | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.common/                     |
| 20 | frmdevcaps.cs                   | buildagent/work/f81ed2ccc88a1e9a/pod/common/cytomatclient/                          |
| 21 | frmmain.cs                      | buildagent/work/f81ed2ccc88a1e9a/pod/common/microserveclient/                       |
| 22 | frmmain.cs                      | buildagent/work/f81ed2ccc88a1e9a/pod/common/ambistoreclient/                        |
| 23 | frmmain.cs                      | buildagent/work/f81ed2ccc88a1e9a/pod/common/cytomatclient/                          |
| 24 | genrobotxcfg.cs                 | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/podadmin/                             |
| 25 | gxdevice.cs                     | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.gx/                 |
| 26 | inhecocpacdevice.cs             | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.inhecocpac/         |
| 27 | inhecomethodfile.cs             | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.inhecoodtc/         |
| 28 | inhecoodtcdevice.cs             | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.inhecoodtc/         |
| 29 | intellixcapdevice.cs            | github/tempo/pod/common/labcyte.pod.devices.intellixcap/                            |
| 30 | jobdevicepoolingcalculations.cs | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.common/job/                 |
| 31 | jobmanager.cs                   | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.common/job/                 |
| 32 | jobmanagerbase.cs               | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.common/job/                 |
| 33 | jobpropertymanager.cs           | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.common/job/                 |

| Id | Filename                  | Path                                                                                |
|----|---------------------------|-------------------------------------------------------------------------------------|
| 34 | labelplateusingfilejob.cs | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.vcode/              |
| 35 | launchprogramcommand.cs   | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/symphony/                             |
| 36 | licenseviewmodel.cs       | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/podadmin/viewmodels/                  |
| 37 | lickey.cs                 | github/tempo/pod/common/labcyte.pod.base/                                           |
| 38 | lightcyclerdevice.cs      | github/tempo/pod/common/labcyte.pod.devices.lightcycler/                            |
| 39 | lookuptableform.cs        | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.vcode/              |
| 40 | mailservice.cs            | buildagent/work/f81ed2ccc88a1e9a/common/labcyte.common.mailer/                      |
| 41 | mainform.cs               | buildagent/work/f81ed2ccc88a1e9a/pod/common/koyoplrc/koyoplrc/                      |
| 42 | mainform.cs               | buildagent/work/f81ed2ccc88a1e9a/pod/common/inhecoclient/                           |
| 43 | mainform.cs               | buildagent/work/f81ed2ccc88a1e9a/pod/common/certusclient/                           |
| 44 | mainwindow.xaml.cs        | buildagent/work/f81ed2ccc88a1e9a/pod/common/koyoclickplc/koyoclickplc/              |
| 45 | microservedriver.cs       | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.microserv<br>e/     |
| 46 | microservedriver.cs       | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.driver.microserve/          |
| 47 | mobilerobotdoordevice.cs  | buildagent/work/f81ed2ccc88a1e9a/pod/common/becls.pod.devices.mobilerobot<br>door/  |
| 48 | modbusserialtransport.cs  | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.driver.koyo/modb<br>us/io/  |
| 49 | plateimportexport.cs      | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.common/                     |
| 50 | platestoragegrid.cs       | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.platestoragecontr<br>ol/    |
| 51 | platestoreengine.cs       | github/tempo/pod/common/labcyte.pod.common/platestorage/                            |
| 52 | podddb.designer.cs        | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.database/                   |
| 53 | poddevice.cs              | github/tempo/pod/common/labcyte.pod.common/device/                                  |
| 54 | podlogger.cs              | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.common/                     |
| 55 | podoptionsform.cs         | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/symphony/                             |
| 56 | podprotocol.cs            | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/protocol/               |
| 57 | podsystem.cs              | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                        |
| 58 | preciseteachingform.cs    | buildagent/work/f81ed2ccc88a1e9a/pod/common/robotteachingutility/                   |
| 59 | program.app.cs            | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/symphony/                             |
| 60 | program.cs                | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/podtester/                            |
| 61 | protocolview.cs           | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/symphony/protocols/                   |
| 62 | querytimepointjob.cs      | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/jobs/                   |
| 63 | readvolumesjob.cs         | buildagent/work/f81ed2ccc88a1e9a/pod/common/becls.pod.devices.volumechec<br>k/jobs/ |
| 64 | register.aspx.cs          | github/tempo/pod/common/robotteachclient/account/                                   |
| 65 | remoteconnector.cs        | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.envision/           |
| 66 | remoteconnector.cs        | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.bmg/                |
| 67 | remoteconnector.cs        | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.asm/                |
| 68 | removetimepointjob.cs     | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/jobs/                   |

| Id | Filename                      | Path                                                                                    |
|----|-------------------------------|-----------------------------------------------------------------------------------------|
| 69 | robotconfig.cs                | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.driver.commonrobot/config/      |
| 70 | robotconfigloader.cs          | buildagent/work/f81ed2ccc88a1e9a/pod/symphony/podadmin/autogen/                         |
| 71 | robotdevice.cs                | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.robot/                  |
| 72 | rorzeteachingfrm.cs           | buildagent/work/f81ed2ccc88a1e9a/pod/common/robotteachingutility/                       |
| 73 | run1536experimentjob.cs       | github/tempo/pod/common/labcyte.pod.devices.lightcycler1536/                            |
| 74 | rundefinition.cs              | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/protocol/                   |
| 75 | runexecutormisc.cs            | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 76 | runexecutorreports.cs         | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 77 | runexecutorruncreation.cs     | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 78 | runoptionscontrol.cs          | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.storageprotocol/                |
| 79 | runprotocoljob.cs             | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.tempest/                |
| 80 | runreporhandler.cs            | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 81 | runreportsform.cs             | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.userinterface/                  |
| 82 | runscheduler.cs               | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 83 | runscheduler.log.cs           | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 84 | runscheduler.scheduler.cs     | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 85 | runscheduler.validationlog.cs | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 86 | runstatusform.cs              | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.userinterface/                  |
| 87 | runsummarygenerator.cs        | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/                            |
| 88 | savetimepointjob.cs           | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/jobs/                       |
| 89 | selectechoprotocol.cs         | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/protocol/                   |
| 90 | staggerplatejob.cs            | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.api/jobs/                       |
| 91 | surveyjob.cs                  | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.echo/                   |
| 92 | tcpclientsocketsim.cs         | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.driver.storex/                  |
| 93 | teleshakedevice.cs            | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.devices.teleshake/              |
| 94 | uccapcheckimages.cs           | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.userinterface/usercontrols/     |
| 95 | ucechoapplications.cs         | buildagent/work/f81ed2ccc88a1e9a/pod/common/labcyte.pod.userinterface/protocolcontrols/ |
| 96 | volumecheckdevice.cs          | buildagent/work/f81ed2ccc88a1e9a/pod/common/becls.pod.devices.volumecheck/              |
| 97 | xcfgfileviewer.cs             | buildagent/work/f81ed2ccc88a1e9a/pod/common/robotteachingutility/                       |

### Appendix E: Referenced Classpaths

| Id | Path                                                                                                                       |
|----|----------------------------------------------------------------------------------------------------------------------------|
| 1  | beckmancoulter_devices_samplestore_dll.BeckmanCoulter.Devices.SampleStore.OAuth.OAuthRequestHandler._3CSendRequest_3Ed__9  |
| 2  | beckmancoulter_devices_samplestore_dll.BeckmanCoulter.Devices.SampleStore.OAuth.OAuthTokenGenerator._3CGetToken_3Ed__1     |
| 3  | beckmancoulter_devices_samplestore_dll.BeckmanCoulter.Devices.SampleStore.OAuth.OAuthTokenGenerator._3CRefreshToken_3Ed__2 |
| 4  | beckmancoulter_devices_volumecheck_dll.BeckmanCoulter.Devices.VolumeCheck.SimVolumeCheckDevice                             |
| 5  | btilhcrunner_dll.g                                                                                                         |
| 6  | btilhcrunner_dll.j                                                                                                         |
| 7  | btilhcrunner_dll.m                                                                                                         |
| 8  | btilhcrunner_dll.n                                                                                                         |
| 9  | labcyte_echo_api_dll.Labcyte.Echo.API.ImportPickListForm                                                                   |
| 10 | labcyte_echo_api_dll.Labcyte.Echo.API.OutputGenerator                                                                      |
| 11 | labcyte_echo_api_dll.Labcyte.Echo.API.PexLogger                                                                            |
| 12 | labcyte_echo_api_dll.Labcyte.Echo.API.PickList.EchoPickList                                                                |
| 13 | labcyte_echo_api_dll.Labcyte.Echo.API.PlateInfo                                                                            |
| 14 | labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.BulkFillOperation                                                           |
| 15 | labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.EchoBase                                                                    |
| 16 | labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.EchoPlateSurvey                                                             |
| 17 | labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.EchoProtocol                                                                |
| 18 | labcyte_echo_api_dll.Labcyte.Echo.API.Protocol.EchoTransferMap                                                             |
| 19 | labcyte_echo_api_dll.Labcyte.Echo.API.SurveyLog.SurveyData                                                                 |
| 20 | labcyte_echo_api_dll.Labcyte.Echo.API.SurveyLog.SurveyHistoryManager                                                       |
| 21 | labcyte_echo_api_dll.k                                                                                                     |
| 22 | labcyte_echo_common_dll.Labcyte.Echo.Common.Output.OutputData                                                              |
| 23 | labcyte_echo_labware_dll.Labcyte.Echo.Labware.PlateInfoManager                                                             |
| 24 | roche_htc1536_aim_api_dll.Roche.HTC1536.AIM.API.HTC1536Aim                                                                 |
| 25 | roche_htc1536_aim_api_dll.Roche.HTC1536.AIM.API.HTC1536AimClient                                                           |
| 26 | roche_htc1536_aim_api_dll.Roche.HTC1536.AIM.API.HTC1536AimSimulator                                                        |
| 27 | roche_htc1536_aim_api_dll.Roche.HTC1536.AIM.API.Simulator.ExperimentFileHelper                                             |
| 28 | sciexdeviceclient_dll.BeckmanCoulter.Devices.EchoMS.SciexDeviceClientWrapper._3CRunSamples_3Ed__19                         |
| 29 | sciexdeviceprotos_dll.BeckmanCoulter.Devices.EchoMS.Protos.BatchFileParser                                                 |